

Arun Sanyasi

report

 Quick Submit Quick Submit College of Science and Technology, Bhutan

Document Details

Submission ID

trn:oid:::1:3583124267

Submission Date

May 31, 2026, 4:26 PM GMT+5:30

Download Date

May 31, 2026, 4:36 PM GMT+5:30

File Name

The_Final.pdf

File Size

7.5 MB

154 Pages

24,811 Words

139,659 Characters

*% detected as AI

AI detection includes the possibility of false positives. Although some text in this submission is likely AI generated, scores below the 20% threshold are not surfaced because they have a higher likelihood of false positives.


Caution: Review required.

It is essential to understand the limitations of AI detection before making decisions about a student's work. We encourage you to learn more about Turnitin's AI detection capabilities before using the tool.

Disclaimer

Our AI writing assessment is designed to help educators identify text that might be prepared by a generative AI tool. Our AI writing assessment may not always be accurate (i.e., our AI models may produce either false positive results or false negative results), so it should not be used as the sole basis for adverse actions against a student. It takes further scrutiny and human judgment in conjunction with an organization's application of its specific academic policies to determine whether any academic misconduct has occurred.

DESIGN AND DEVELOPMENT OF BASIC HEALTH MONITORING STATION

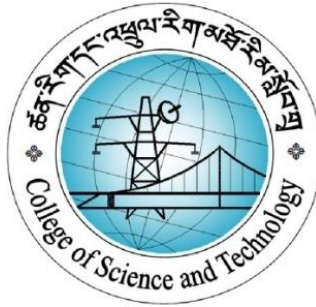

**Project report submitted in partial fulfillment of the requirements for the award of
the degree of Bachelor Engineering in Instrumentation and Control Engineering**

By

Mr. Arun Sanyasi

Ms. Sherab Zangmo

Ms. Tashi Chozom

Under the guidance of:

Mrs. Karma Kelzang Eudon (Guide)

Mr. Kinley Penjor (External Guide)

ELECTRICAL AND ELECTRONIC ENGINEERING DEPARTMENT

COLLEGE OF SCIENCE AND TECNOLOGY

RINCHENDING: PHUENTSHONLING

CHHUKHA DZONGKHAG: BHUTAN

June 2026

ROYAL UNIVERSITY OF BHUTAN
COLLEGE OF SCIENCE AND TECHNOLOGY
DEPARTMENT OF INSTRUMENTATION AND CONTROL ENGINEERING



CERTIFICATE

This is to certify that the B.E. project titled **“Basic Health Monitoring System”**, which is being submitted by **Mr. Arun Sanyasi(02220230)**, **Ms. Tashi Chozom(02220244)**, and **Ms. Sherab Zangmo (02220238)** in partial fulfillment of the requirement for the award of the degree of Bachelor of Engineering in Instrumentation and Control

Mr. Kinley Penjor
External Guide

Ms. Karma Kelzang Eudon
Project Guide

ACKNOWLEDGEMENT

We acknowledge guidance and support of each other at every stage of the development of our project “Basic Health Monitoring Station”. First, we would like to thank the College of Science and Technology, and Instrumentation and Control Engineering program which provided the facilities, resources and conducive working environment needed to successfully conclude this project.

We are grateful to our guide, Ms. Karma Kelzang Eudon, as she has been with us during project development as well. In this regard, she supported us with wisdom and encouragement with no delay. Her suggestions and technical knowledge have made us able to finish this task successfully.

Finally, we would like to extend our heartfelt appreciation to our external guide Mr. Kinley Penjor for taking the time required to peruse our work and providing us guidance, with positive comments, suggestions and continued assistance throughout the work process.

We would also like to thank Ms. Deepa Sharma who’s currently pursuing Bachelor of Nursing at Royal Academy of Thimphu who helped us understand the linkages between each parameter and help us understand the core concept of providing diagnoses and health tips for our applications.

We would like to express our gratitude to Chief Medical Officer and the staff of Phuentsholing General Hospital for providing the working space and access to all the clinical devices we needed to calibrate and test our system. They were generous enough to let us calibrate our device when necessary. We appreciate the team at Fab-Lab CST and their support to us in using their 3D printing capabilities while doing our project. Timber Lab staff, Mr.Vijay, helped us during the fabrication of the physical body of our system, therefore, we would like to thank him for his support. Finally, we would like to thank our fellow classmates, friends and family members for the encouragement, and support throughout this process.

Arun Sanyasi

Tashi Chozom

Sherab Zangmo

ABSTRACT

As the importance of preventative healthcare and early diagnosis of medical conditions continues to grow, there is an urgent need for effective methods of testing for health issues. In most establishments, pulse rate, respiratory rate, SpO₂, body temperature, body height, and body weight must be measured individually through different pieces of equipment, making the process longer and necessitating manual recording. The purpose of this research project is to develop a prototype of a Basic Health Monitoring Station. The system utilizes the ESP32 microcontroller for live acquisition of the sensor data and Raspberry Pi 4B for data processing, graphical interface management, and report generation. The system incorporates sensors like MAX30102 for measuring pulse rate and SpO₂, MLX90614 for non-contact measurement of body temperature, VL53L0X for calculating body height, and a load cell sensor with the HX711 amplifier for body weight calculation. The BMI is automatically calculated using height and weight measurements. The output is categorized as low, normal, or high. The output is presented on the touch screen display. Finally, the QR code is generated to produce the health report. Thus, the total cost of building this device is Nu. 20,813, making it economical and feasible for colleges, rural hospitals, community centers, and offices.

List of Abbreviations

Abbreviation	Full Form
BMI	Body Mass Index
SpO ₂	Blood Oxygen Saturation
BPM	Beats Per Minute
°C	Degree Celsius
ESP32	Espressif 32-bit Microcontroller
GPIO	General Purpose Input/Output
LCD	Liquid Crystal Display
GUI	Graphical User Interface
QR	Quick Response
IoT	Internet of Things
ToF	Time of Flight
ADC	Analog-to-Digital Converter
I ² C	Inter-Integrated Circuit
IDE	Integrated Development Environment
USB	Universal Serial Bus

List of Figures

<i>Figure 1 : Methodology</i>	<i>17</i>
<i>Figure 2 : System Overview</i>	<i>18</i>
<i>Figure 3 : System Flow Chart</i>	<i>19</i>
<i>Figure 4 : Schematic Diagram</i>	<i>24</i>
<i>Figure 5 : Simulation Circuit</i>	<i>25</i>
<i>Figure 6 : Communication Process Flow of the System</i>	<i>27</i>
<i>Figure 7 : User Interface Home Screen</i>	<i>32</i>
<i>Figure 8 : User Interface Measurement Screen</i>	<i>32</i>
<i>Figure 9 : User Interface Result Display Screen</i>	<i>32</i>
<i>Figure 10 : QR Code Process</i>	<i>33</i>
<i>Figure 11 : System Model</i>	<i>34</i>
<i>Figure 12 : Views of 3D Design</i>	<i>35</i>
<i>Figure 13 : Woodwork</i>	<i>39</i>
<i>Figure 14 : Integration of Load Cells and HX711</i>	<i>41</i>
<i>Figure 15 : Integration of Sensors in Sensor Box</i>	<i>41</i>
<i>Figure 16 : Mounting the temperature sensor</i>	<i>41</i>
<i>Figure 17 : Final Assembly of Basic Health Monitoring Station</i>	<i>42</i>
<i>Figure 18 : Testing of ESP32</i>	<i>46</i>
<i>Figure 19 : Testing of LCD and Raspberry pi 4 with Raspberry pi</i>	<i>46</i>
<i>Figure 20 : Testing of HX711 with load cell</i>	<i>46</i>
<i>Figure 21 : Testing of MLX90614ESF-DCI</i>	<i>47</i>
<i>Figure 22 : Testing of VLX5L0X</i>	<i>47</i>
<i>Figure 23 : Testing of Pulse rate sensor</i>	<i>47</i>
<i>Figure 24 : Testing of MAX301202</i>	<i>47</i>
<i>Figure 25 : Integrated System Testing</i>	<i>54</i>
<i>Figure 25 .1: Uploading of APPLICATIONS</i>	<i>55</i>

TABLE OF CONTENT

ACKNOWLEDGEMENT	1
ABSTRACT.....	1
List of Abbreviations	2
List of Figures	3
CHAPTER 1: INTRODUCTION	9
1.1 Background	9
1.2 Problem statement.....	10
1.3 Motivation	11
1.4 Expected Outcomes	12
1.5 Aim and Objectives.....	12
1.6 Literature Review	12
CHAPTER 2: SYSTEM DESIGN	15
2.1 Methodology.....	16
2.2 System Overview.....	17
2.3 System Flow Chart.....	18
<i>Figure 3: System Flow Chart</i>	19
2.4 Voice-Assisted User Flow	19
CHAPTER 3: HARDWARE AND SOFTWARE DESIGN	21
3.1 Component Requirement	21
3.1.1 Hardware Components	21
3.1.2 Software Components	22
3.2 Schematic Diagram.....	23
3.3 Simulation Circuit.....	25
3.4 Software Architecture and System Implementation	25
3.5 Data Acquisition Using ESP32.....	26
3.6 Communication Between ESP32 and Raspberry Pi	26
3.7 Feeding Data Using Sensors.....	28

3.7.3 SpO ₂ and Pulse Rate Data	28
3.7.5 Temperature Data.....	29
3.8 Back-end Processing in Raspberry Pi.....	29
3.9 BMI Calculation	30
3.10 Classification of Health Status Parameters	31
3.11 User Interface Design.....	31
3.12 QR Code and Health Report Creation.....	32
CHAPTER 4: MECHANICAL, DESIGN AND FABRICATION	33
4.1 Mechanical Design Overview	33
4.2 3D Model of the System Body	34
4.3 Dimensions of Components of the System Body	36
4.5 Modification of Body Frame	40
4.7 Assembly of the System.....	40
4.7 Material Selection.....	43
4.8 Planning for the College Infirmary and User Guidance	43
CHAPTER 5: TESTING AND CALIBRATION.....	45
5.1 Individual component testing	46
5.2 Integrated testing	52
5.4 Uploading of APP to the Capacitive Touch Display Screen.	54
6.1 Introduction	56
6.2 ESP32 and Raspberry Pi Programming	56
6.3 GUI Design.....	57
6.5.1 Saved Reports Screen.....	76
6.5.2 Storage Maintenance Screen.....	77
6.5.3 <i>Settings Screen</i>	78
6.5.4 Admin Panel.....	79
6.5.5 Admin Login Screen.....	79

6.5.6 Instruction and Mode Selection Screen	80
6.5.7 Welcome Screen.....	81
6.5.8 Measurement Result Screen	81
6.5.9 Health Diagnosis Screen	82
6.5.10 Consultation Guidance Screen.....	83
6.5.11 Temperature Detail Screen.....	84
6.5.12 Blood Oxygen Detail Screen	84
6.5.13 Pulse Rate Detail Screen	85
6.5.14 BMI Detail Screen	86
6.6 Voice Assistance Module	Error! Bookmark not defined.
6.7 Serial Communication.....	Error! Bookmark not defined.
6.8 Data Processing Algorithm.....	Error! Bookmark not defined.
6.9 BMI Calculation	Error! Bookmark not defined.
6.10 Health Classification Logic	90
6.11 QR Report Generation.....	Error! Bookmark not defined.
CHAPTER 7: TESTING WITH NUMBER OF SAMPLES	93
7.1 Introduction	93
7.2 Individual Component Testing	Error! Bookmark not defined.
7.3 Calibration Procedures.....	Error! Bookmark not defined.
7.3.1 MAX30102 Calibration.....	Error! Bookmark not defined.
7.3.2 MLX90614 Calibration	Error! Bookmark not defined.
7.3.3 Pulse Sensor Calibration.....	Error! Bookmark not defined.
7.3.4 VL53L0X Calibration	Error! Bookmark not defined.
7.3.5 Load Cell Calibration	96
CHAPTER 8: RESULTS AND DISCUSSION	98
8.1 Introduction	Error! Bookmark not defined.
8.2 Integrated System Testing	Error! Bookmark not defined.

8.3 Data Collection	Error! Bookmark not defined.
8.4 Comparison with Medical Devices.....	Error! Bookmark not defined.
8.5 Error Analysis	104
8.6 Time Efficiency Analysis	Error! Bookmark not defined.
8.7 Performance Evaluation	Error! Bookmark not defined.
8.8 Trend Analysis	Error! Bookmark not defined.
8.9 Feasibility Analysis	Error! Bookmark not defined.
8.10 Cost Analysis	115
CHAPTER 9: CHALLENGES AND LIMITATIONS	118
9.1 Challenges Faced.....	Error! Bookmark not defined.
9.2 System Limitations	Error! Bookmark not defined.
CHAPTER 10: CONCLUSION AND FUTURE SCOPE.....	120
10.1 Conclusion.....	120
10.2 Contributions of the Project	120
10.3 Future Scope.....	121
REFERENCES.....	122
APPENDICES	124
Appendix A: Complete BHMS Data Analysis Graphs	124
Appendix B: Additional Notes on Voice Assistance...	Error! Bookmark not defined.
Appendix C: Detailed Gui Screen Functional Documentation.	Error! Bookmark not defined.
Appendix D: Consolidated Data Interpretation And Workflow Notes	147
Appendix F: Source Code Summary And Key Program Excerpts..	Error! Bookmark not defined.

CHAPTER 1: INTRODUCTION

1.1 Background

Countries around the world are now prioritizing preventive care over reactive management. Detection of any unusual health condition at an early stage is very important to decrease the disease burden as well as reduce medical expenses and enhance patient's health. Periodic observation on vital parameters (pulse rate, respiratory rate, blood oxygen saturation (SpO₂), body temperature, height, weight and Body Mass Index (BMI)) is of significant interest for understanding the health conditions of a patient. These measurements are important for determining cardiovascular health, respiratory efficiency, metabolic performance and well-being. Vital signs measurements based on independent devices are also carried out in most institutions such as colleges, offices, and smaller nursing clinics. For instance, temperature can be obtained by thermometer, SpO₂ by a pulse oximeter, weight by digital weighing-scale, height by a hand measuring scale. Then BMI is calculated by manually calculating height and weight values. This stepwise approach results in unnecessary screening time, cumbersome workflow, manual records of data, and human error. Moreover, it gets tedious to maintain organized digital health records when measurements are taken separately. The rapid progress of embedded systems and microcontrollers, coupled with biomedical sensors, allows the integration of multiple functions like health monitoring to create a single automated platform. Microcontrollers like ESP32 enable real-time acquisition of data, while single-board computers such as Raspberry Pi enable advanced processing, graphical interface development, and data management capabilities. These modern biomedical sensors offer precise, small-volume, and budget-friendly methods for physiological measurement. This Basic Health Monitoring Station aims to replace the drawbacks of mainstream health screening techniques by establishing an integrated, automated health kiosk integrated together with the capabilities to measure several crucial parameters in one system. It incorporates hardware integration, embedded programming, and graphical user interface development to design an intuitive and compact health monitoring solution. Pulse rate and SpO₂ are measured with MAX30102, body temperature with MLX90614 non-contact sensor, height with VL53L0X time-of-flight sensor, and weight with load cell with HX711 amplifier. BMI is calculated automatically from measurements of height and weight. The collected sensor data are processed and divided into health categories defined as low, normal, or high. The results are displayed on screen through a touchscreen, while a QR code function allows users to

access a detailed digital self-report. This seamless strategy decreases screening time, increases operational effectiveness and facilitates systematic documentation of digital health management. The system is a portable, affordable, and deployable tool suitable in not only college infirmaries but also rural health centers, community centers, corporate offices and public health camps. At a design cost of Nu. 24,143, this project is affordability and scalability based and can be used in Bhutanese setting. In solving this project, one can use their practical knowledge with embedded systems, biomedical devices, circuit design, software development, system implementation, and testing methods to solve a real-world problem. The Basic Health Monitoring Station showcases an undergraduate education for engineering in technology-based solution that integrates creativity, pragmatism, and social impact to help in preventive healthcare, and this is the main objective to be met with the project.

1.2 Problem statement

Monitoring vital signs is one of the most basic components of primary care. Pulse rate, respiratory rate, blood oxygen saturation (SpO_2), body temperature, height, weight, and Body Mass Index (BMI) can also provide useful data about a person's physiological condition. These parameters are measured on separate and independent devices that do not correlate with one another in many institutional settings, like college infirmaries, community health centers, and small medical units. The piecemeal approach to health screening poses several operational as well as technological hurdles. Firstly, using multiple standalone instruments makes one health check-up take longer. Every parameter needs to be measured individually, recorded manually, and interpreted. Due to that, the process becomes inefficient, particularly during health camps, peak hours, or emergency situations where multiple people need screening. Screening time is also negatively impacted in the sense that it becomes time-consuming and creates extra work burden among healthcare professionals. Second, human error is involved in noting the manual data of recording and computation. For example, BMI is a value that must be calculated manually using measured height and weight. Any flaw in recording or calculation can lead to misrepresentation of health status. The same goes for paper-based data, which can be misplaced, damaged, and are not consistent. Lack of a digital recording system eliminates the ability to keep an organized and reliable system for health data collection.

The other problem is that measurement and reporting systems lack integration. Conventional systems offer individual measurements and do not provide a summative summary of health status. There is no comprehensive system that combines different parameters and makes them into clear indicators, such as low, normal, or high state. As a result, preliminary screening of this sort may well be less organized and less informative. Additionally, many existing integrated health monitoring systems are quite pricey and are largely for use at hospital levels. These systems may not be cost-effective for Bhutanese micro-institutions, health systems in rural areas, or schools and communities. This lack of a cheap, locally designed approach creates a disconnect between technology and practical availability.

Also, given the increasing importance of preventive healthcare and digital transformation, there is demand for automated, technology-enabled health solutions. Institutions cannot promote self-assessment, digital record keeping, or structured health monitoring without a well-integrated system. Thus, the fundamental issue identified in this project is the absence of an inexpensive, automated, and integrated health monitoring system able to acquire numerous key parameters that are tracked in a singular compact device, providing real-time digital results and reporting in a defined way. This problem is addressed in the introduction of the Basic Health Monitoring Station, which seeks to provide a unique system for comprehensive biomedical sensors, embedded systems, automated data processing, and digital reporting as a low-cost and comprehensive solution for institutional and community distribution.

1.3 Motivation

The inspiration behind this project is to design an inexpensive and comprehensive health screening device for basic health screening purposes. Many universities have their own clinics, while many others use community health clinics. All of them employ individual equipment to screen vital parameters like heart rate, SpO₂, body temperature, height, weight, and BMI.

The motivation behind the implementation of the project lies in the requirement for an effective yet simple and economical health monitoring system used in basic health screenings. For instance, in college infirmaries, community health centers, and small medical facilities, vital signs such as pulse rate, SpO₂, body temperature, body height, weight, and BMI are monitored using individual devices. As a result, the whole process becomes quite cumbersome, which may lead to errors when documenting the data manually.

To address this issue in Bhutan, there is a need for a system capable of calculating various health parameters, generating real-time results, determining BMI automatically, and providing digital documentation capabilities.

This project will also allow me to utilize the knowledge I gained about embedded systems, biomedical sensors, electronic circuitry, programming skills, and other techniques in solving a real-world problem in the field of health care.

1.4 Expected Outcomes

The BMI, pulse rate, SpO₂, breathing rate, and body temperature will be measured using the proposed system. The values will be provided through an intuitive display system, along with a health index that determines whether the individual is healthy or requires any medical help.

1.5 Aim and Objectives

1.5.1 Aim: To design a health monitoring system that measures four vital signs and BMI for basic check up.

1.5.2 Objectives:

1. To design an integrated system for the measurement of pulse rate, respiratory rate, SpO₂, temperature, and BMI.
2. To publish vital signs and health insights using ESP32 and Raspberry Pi.
3. To install the system in the infirmary for quick and easy health check-ups.

Together, these objectives ensure that the project is technically feasible, practically applicable, and academically aligned with undergraduate engineering learning outcomes. The defined objectives also provide a structured road map for successful implementation and evaluation of the Basic Health Monitoring Station.

1.6 Literature Review

The literature review includes previous research on health monitoring systems, biomedical sensors, healthcare kiosks, and Internet-of-Things (IoT) based health care. The emphasis will be made on measuring pulse rate, respiration rate, SpO₂, body temperature, and Body Mass Index (BMI) in an automated manner.

1.6.1 Real Time Vital Sign Monitoring System

In this research paper, Dhingra & Sharma (2019), designed a real-time vital sign monitoring system with the help of microcontroller-based sensors. This paper was based on the development of an embedded system for real-time measurements of pulse rate and temperature, which are important health parameters. Thus, it can be said that embedded systems can be effectively used to perform real time monitoring of fundamental human health conditions.

This paper had a limitation in that their work was only based on certain health parameters; there was no measure of BMI and reporting was absent from the system. Thus, this research paper is relevant to the proposed project as this research paper provides the basis for the proposed Basic Health Monitoring Station. The major difference between the two projects is that the proposed Basic Health Monitoring Station will include all the mentioned parameters of this project along with BMI, height, and weight.

1.6.2 Automated Body Mass Index System

Patrick O. and Kazeem A. (2020) developed an Automated Body Mass Index system to reduce manual errors in BMI calculation. Their system used load cell and height measurement technology to calculate BMI automatically. The main advantage of this system was that it reduced the need for manual measurement and manual calculation. It also improved accuracy and made BMI assessment faster and more convenient.

Although the system was useful for BMI measurement, it was limited only to one health parameter. It did not measure pulse rate, respiratory rate, oxygen saturation, or body temperature. Therefore, it could not provide a broader health condition assessment. This work is important for the present project because BMI is one of the main parameters measured by the proposed station. The present project improves this by integrating BMI measurement with other vital signs in a single health monitoring station.

1.6.3 Guidelines for Recording Vital Signs

The Cleveland Clinic (2025) was the source that supplied information related to the documentation of vital signs including the pulse rate, respiration rate, temperature, and

oxygen saturation levels. The source highlights the normal ranges for vital signs and aids in distinguishing between various readings to know whether they fall within the acceptable health range. Standard values of health parameters are significant when designing the health monitoring system because the system will compare readings with health parameters.

This source cannot serve as a hardware implementation study because it does not contain details regarding the development of the monitoring system or the design process. The source contains information concerning standard reference values for different health parameters. For the current project, the guideline will help in specifying the classification values for parameters like the pulse rate, respiratory rate, temperature, and SpO₂.

1.6.4 Healthcare Kiosk Systems Utilizing Vital Signs

The Bhutani et al. study from 2025 analyzed 36 papers on healthcare kiosks for vital sign measurements such as heart rate, oxygen saturation, temperature, and body mass index (BMI). Healthcare kiosks have many benefits but also some problems like technology and cost constraints. The relevance of the study lies in the fact that the suggested Basic Health Monitoring Station kiosk is of the same kind.

eMedicineHealth provides information about determining the normal SpO₂ level according to age in order to distinguish between normal and critical values of oxygen saturation. Although there are no technical features described, it can be useful for recognizing user oxygen saturation levels.

1.6.5 Patient Evaluation and Measurement of Vital Signs

According to the AAST Technical Guideline (2022), patient evaluation and measurement of vital signs were mentioned in the document. The document stresses that it is important to ensure accuracy in measuring, documenting, and recording the measured vital signs properly. Reliable patient evaluation is not only connected with measurement but also needs correct documentation for the future evaluation process.

The limitation in this case is the lack of attention to the use of hardware equipment in order to automate and digitize the measurement process. The technical guidelines provided clinical protocols in the procedure of vital signs measurement. The current research is related to the project as it includes development of hardware for structured and digitized presentation of data collected.

1.6.6 Integrated Biomedical Sensor Systems

In their research on integrated biomedical sensor systems based on microcontrollers, Sharma and Singh (2021) revealed that the combination of various sensors into one system increases the effectiveness of health monitoring. Multisensor systems are valuable since they make it possible to measure a variety of health indicators simultaneously, instead of requiring numerous standalone tools. Thus, the time and convenience of users are improved. Unfortunately, the above-mentioned paper did not pay sufficient attention to the cost of development, portability, and ease of installation in organizations. Such aspects are critical when creating a device for a college infirmary or health screening zone. In contrast, the current project involves the application of this idea of multisensor integration in a real-life health monitoring station powered by ESP32, Raspberry Pi, and a touchscreen display.

1.6.7 Preventive Health Monitoring System Using IoT

According to Chen et al. (2019), IoT can be used remotely to monitor people's health, which offers many advantages, including data sharing and accessibility and prolonged life. However, using IoT may increase costs and complexity, as well as pose serious problems related to privacy. Institutional systems require reliability and easy use. This paper contributes to the current research by recommending specific aspects of IoT that should be improved, although IoT is not used at present.

1.6.8 Summary of Literature Review

In most cases, current systems deal with either vital signs monitoring or body mass index computation but not both at the same time. While some previous studies have made measurements such as pulse rate and temperature using microcontrollers, and other studies solely concentrated on BMI computation, there is yet to be seen a combination of all these features along with medically accepted limits within the same hardware system. Thus, an existing gap within the body of knowledge is that there is no device which acts as a health monitoring station and can measure body mass index, pulse rate, breathing rate, oxygen saturation, and body temperature together.

CHAPTER 2: SYSTEM DESIGN

This chapter presents the overall design logic of the Basic Health Monitoring Station. The system design was prepared before fabrication and programming so that the hardware arrangement, software workflow, user interaction, voice guidance, report generation and data flow could work as one complete system. The design stage was important because the project combines biomedical sensors, embedded control, Raspberry Pi processing, touchscreen interaction and QR-based report access in a single station.

The system was designed with three main requirements. First, the station should measure multiple health parameters in one guided sequence. Second, the interface should be simple enough for users without technical knowledge. Third, the measured data should be stored and presented as a clear health report. These requirements guided the selection of sensors, controller, display, mechanical structure and software architecture.

2.1 Methodology

Methodology followed the step-by-step engineering design process. It began with the identification of the problem in basic health screening, selecting measurable parameters, reviewing literature, choosing appropriate sensors, designing circuits and mechanics, developing software, testing, calibration, and finally installation. Each of the steps was performed successfully before proceeding to the next step so that problems would be identified and rectified.

The suggested methodology also takes into account usability. As the proposed health station is meant to be used in an infirmary setting, the user needs to be capable of conducting measurements almost without any help. That is why such aspects as GUI development, button functionality, voice assistance, report creation and accessing QR code became an integral part of the development rather than add-on features.

Testing and calibration took place prior to the installation of the system in the college infirmary. That allowed ensuring that all components including sensors, ESP32 communication, Raspberry Pi application, touchscreen, report creation, and voice assistance have been properly tested.

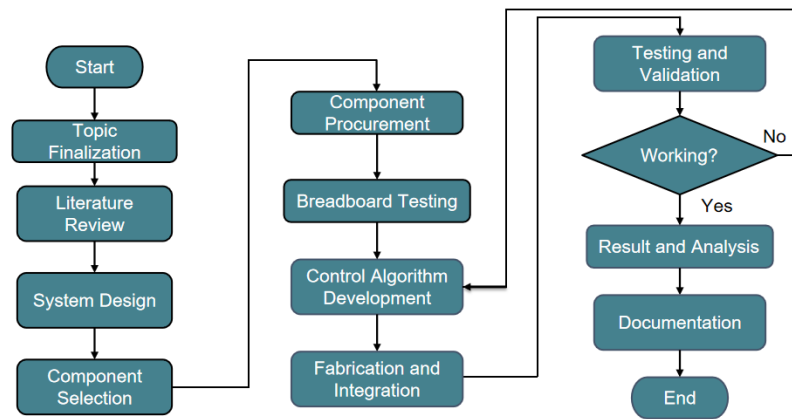


Figure 1: Methodology

2.2 System Overview

The Basic Health Monitoring Station consists of four layers: sensing layer, controller layer, processing and interface layer, and reporting layer. In the sensing layer, there is the MAX30102, MLX90614, VL53L0X, pulse sensor, and load-cell assembly. The controller layer uses ESP32 to read sensor values and prepare labeled packets of data. In the processing and interface layer, the Raspberry Pi 4B operates the PyQt software, calculates BMI, makes readings classification, and provides control of the touchscreen. In the reporting layer, there is the generation of a PDF health report and a QR code.

The station is able to monitor height, weight, BMI, body temperature, SpO2 level, pulse rate, and respiratory rate. Height and weight parameters are used for automatic calculation of the BMI. Other parameters are displayed as health indicators and used for the diagnosis logic. The results will appear on the results screen and can be opened as a PDF report via the QR code viewer.

The voice assistant is part of the system overview since it is a component within the measurement process. Through the DFPlayer Mini and speaker, verbal instructions will be delivered telling the user when to turn the temperature screw, put the left finger, put the right finger, remain steady during weight measurement, stand up for height, and see the results.

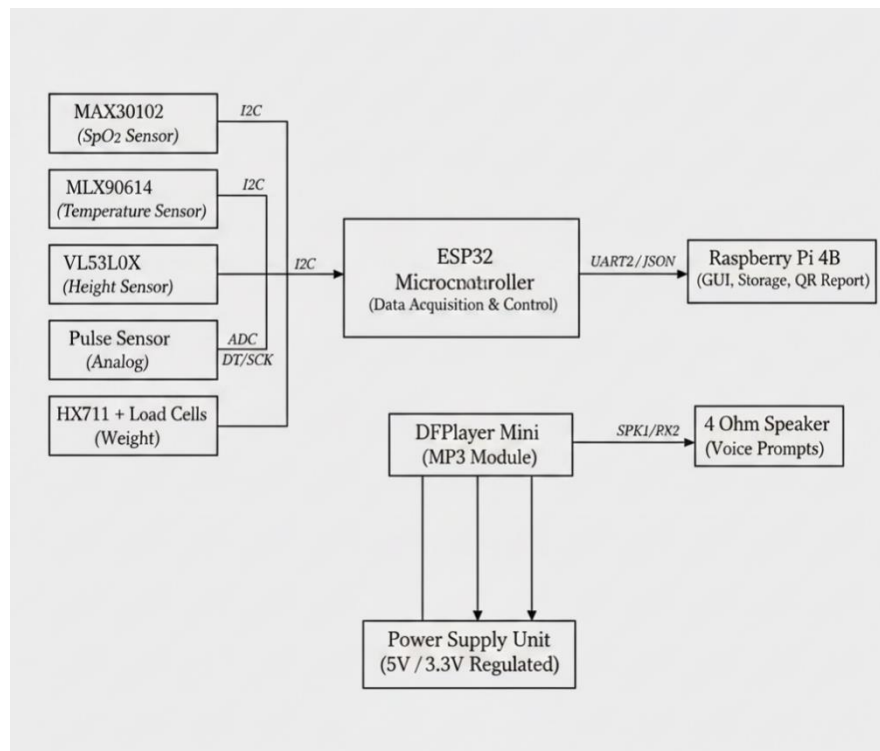


Figure 2: System Overview

2.3 System Flow Chart

The system flow starts when the user launches the application from the welcome screen. After the instruction screen, the user starts the check-up in hardware mode. The station then guides the user through each measurement stage in a fixed order. The chosen order is temperature measurement, SpO2 and respiratory-rate measurement, pulse measurement, weight measurement, height measurement, BMI calculation, result display, diagnosis and report generation.

The flow is arranged in this order to improve usability and reduce movement during measurement. The upper-body sensors are used first, followed by weight and height measurement. After all readings are captured, the Raspberry Pi verifies and displays the final values. The user can then select Report, Open QR, Diagnosis, Consult, Retake or Refresh depending on the need.

The final flow of the system is therefore not only a sensor-reading sequence but a complete user-service sequence. It includes measurement, processing, classification, report

generation, QR access and storage. This makes the BHMS suitable for a self-service screening environment.

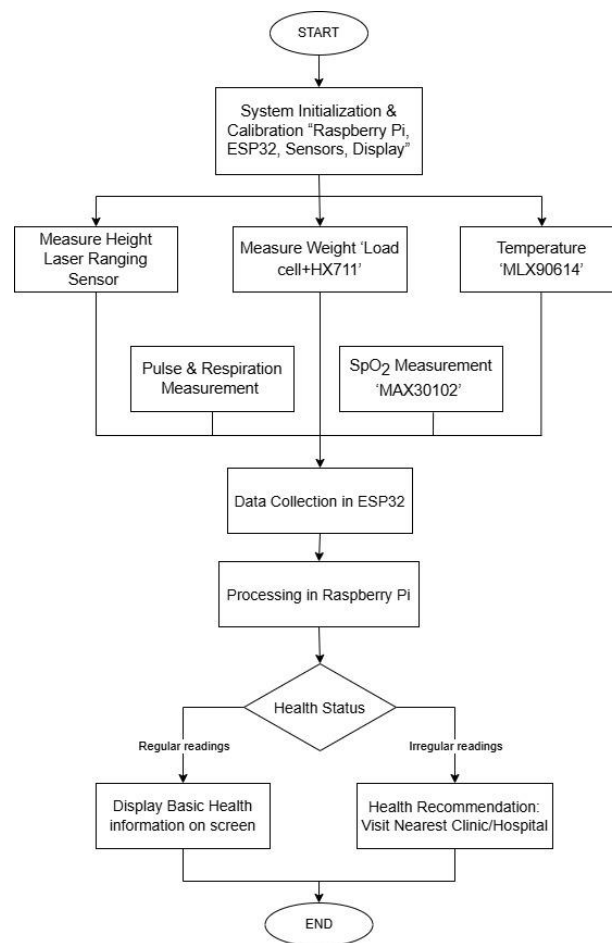


Figure 3: System Flow Chart

2.4 Voice-Assisted User Flow

The voice-assisted flow was added to reduce user confusion and to make the station easier to operate without continuous staff supervision. The voice prompts are synchronized with the measurement stages. The prompts do not replace the GUI; instead, they support the GUI by giving spoken instructions at the exact time when the user needs to move to the next action.

Table 2.1: Voice-assisted measurement flow used in the BHMS.

Stage	User Action	System Response
Welcome	User enters the kiosk interface	Plays welcome prompt and opens the instruction screen
Temperature	User adjusts temperature screw over the head	Records non-contact temperature using MLX90614
SpO2 and RR	User places left index finger in the sensor box	Records oxygen saturation and respiratory rate
Pulse	User places right index finger in the pulse sensor hole	Records BPM and confirms reading
Weight	User stands still on the platform	Records stabilized load-cell value
Height	User stands straight under the height sensor	Records height and calculates BMI
Result and Report	User views result, opens diagnosis or scans QR	Generates PDF report and QR viewer page

CHAPTER 3: HARDWARE AND SOFTWARE DESIGN

3.1 Component Requirement

3.1.1 Hardware Components

Sl.No	Name	Voltage (V)	Current (mA)	Image
1.	Raspberry Pi 4	5	3000	
2.	Capacitive Touchscreen LCD Display	5	1500-2000	
3.	ESP32 WROOM-32	3	500	
4.	4PCS 50KG Load Cell	2.6-5.5	1.5	
5.	VL53L0X- Laser ranging sensor	2.6-5.5	5	
6.	MAX30102	3.3-5	0.6	

7.	MLX90614ESF-DCI	3.3	1.5	
8.	Mini UPS Power Bank	5	3000	
9.	Jumper wire			
10.	15m electric wire			

3.1.2 Software Components

Sl.No	Software	Description	Logo
1.	Arduino IDE	An open-source software used to write, compile, and upload code to Arduino boards and other compatible microcontrollers. It provides a simple interface for programming	

		and debugging embedded systems.	
2.	AutodeskFusion 360	Cloud-based CAD, CAM, and CAE software used for 3D modeling, product design, simulation, and manufacturing.	
3.	Proteus 9	Electronic design automation (EDA) software used for circuit design, simulation, and PCB layout. It allows users to simulate microcontroller-based circuits before hardware implementation.	
4.	Fritzing	Open-source electronics design tool that helps users create circuit diagrams, breadboard layouts, and PCB designs, making it easier to document and prototype electronic projects.	
5.	Easy EDA	Web-based electronic design automation tool used for schematic capture, circuit simulation, and PCB design. It allows users to design and order PCBs directly from the platform.	
6.	VS code studio	A powerful integrated development environment (IDE) used for developing software applications, debugging code, and building programs in various programming languages.	
7.	Cursor IDE	A IDE that enables the usage of the PyQt6 with the ability to correct and remove the delays in full functionality of the APP,	

3.2 Schematic Diagram

The ESP32 acts as the main acquisition controller and communicates with the Raspberry Pi through UART to transfer processed sensor readings. The MAX30102, MLX90614, and

VL53L0X sensors communicate over the I²C bus, while the pulse sensor uses an analog input and the HX711 amplifier receives the load cell signals.

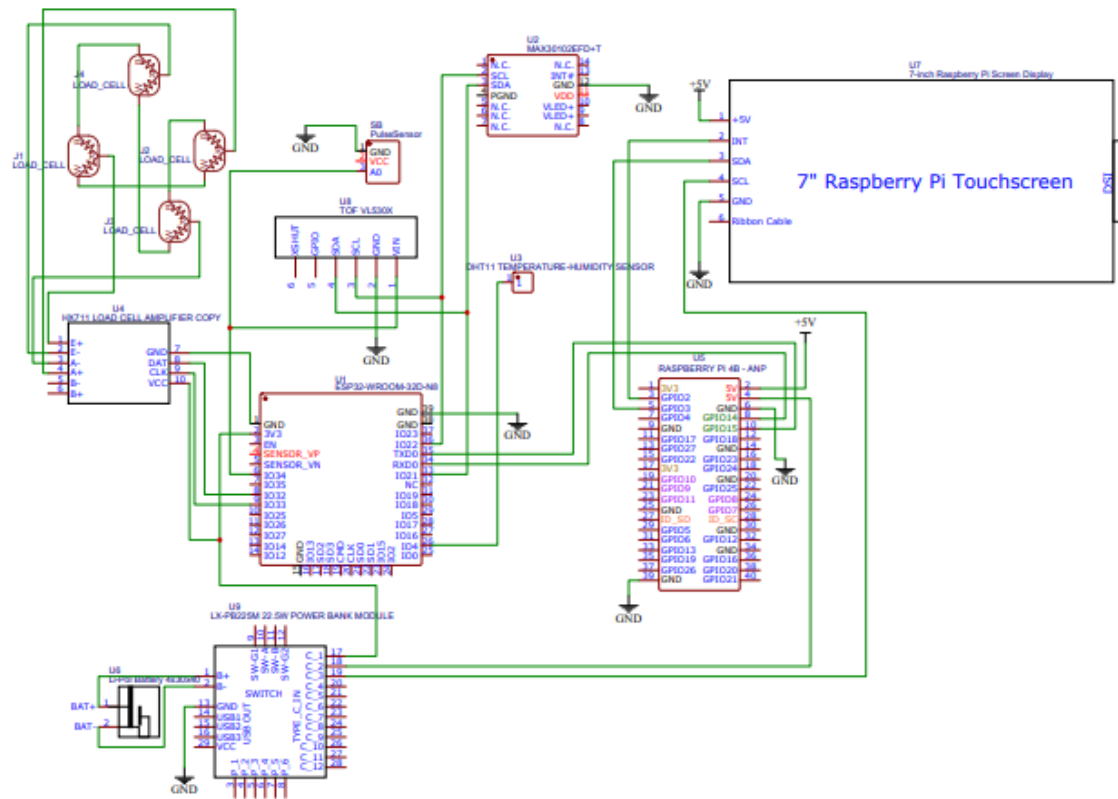


Figure 4: Schematic Diagram

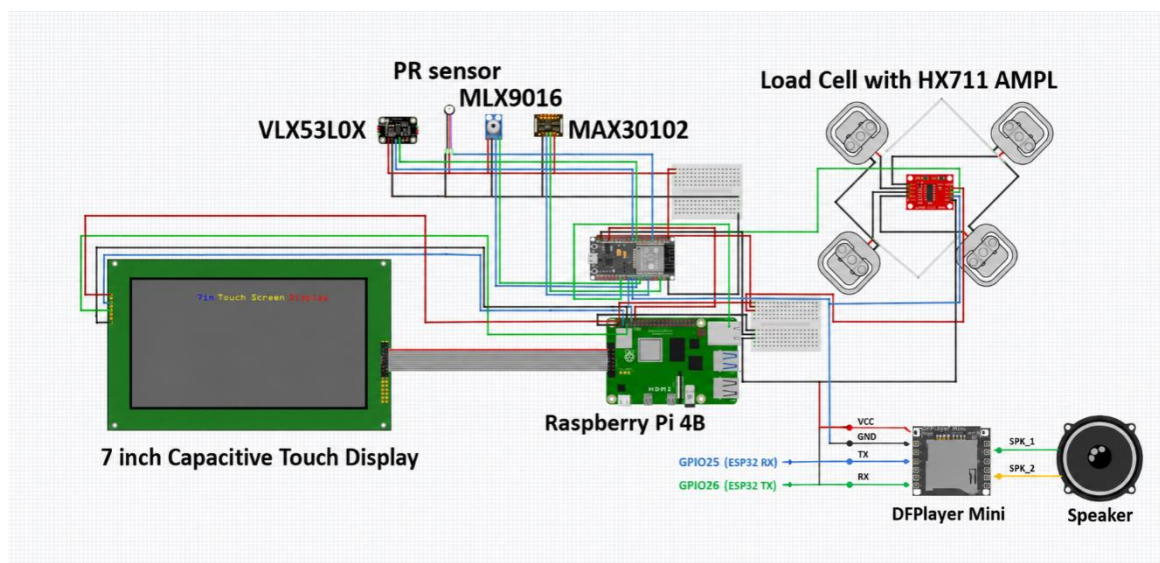


Figure 4.1: Circuit diagram showing capacitive touch display, Raspberry Pi 4B, ESP32, sensors, and load cell with HX711 amplifier connections.

3.3 Simulation Circuit

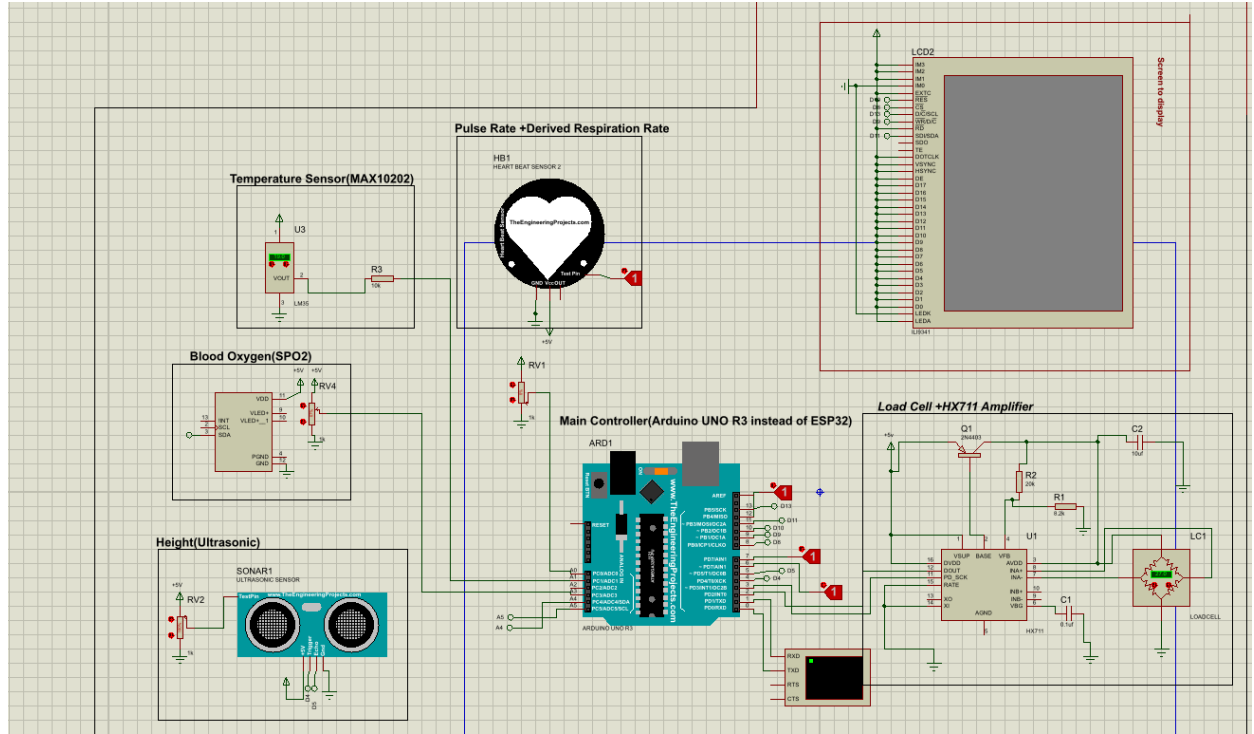


Figure 5: Simulation Circuit

3.4 Software Architecture and System Implementation

The software architecture of the Basic Health Monitoring Station includes all hardware elements in one system including sensors, ESP32 board, Raspberry Pi 4B, LCD screen, and reporting.

The functions were assigned to two control units, with the ESP32 assigned to collect data from such sensors as the MAX30102, MLX90614, VL53L0X, and the load cell via HX711. The Raspberry Pi 4B was assigned to process and output data which include data collection, correction, BMI calculation, interface management, and reporting by QR codes.

It made the development process much simpler since the sensor data collection was differentiated from other operations.

3.5 Data Acquisition Using ESP32

The ESP32 microcontroller was selected for real-time data acquisition from all sensors due to its low cost, compact size, and adequate processing capability. It collects raw readings from each sensor, applies calibration, and prepares the data for transmission to the Raspberry Pi. Specifically, the ESP32 performs the following functions:

- collecting readings from biomedical and measuring sensors;
- receiving information about weight using the HX711 load cell amplifier;
- receiving distance information related to height using the VL53L0X sensor;
- collecting pulse rate and SpO₂ values from the MAX30102 sensor;
- collecting temperature information using the MLX90614 sensor;
- preparing collected values for transfer to the Raspberry Pi;
- managing real-time data collection of the whole system.

Thus, the ESP32 acts as an intermediary between the hardware sensors and the main computing unit.

3.6 Communication Between ESP32 and Raspberry Pi

There will be an essential need for communication between ESP32 and Raspberry Pi in order for the Basic Health Monitoring Station to operate effectively. This is because ESP32 records the values of sensors, while the Raspberry Pi analyzes and shows the results; hence, there will be a need for communication. By communication, ESP32 sends values of sensors to the Raspberry Pi in the form of code, which recognizes them as individual values.

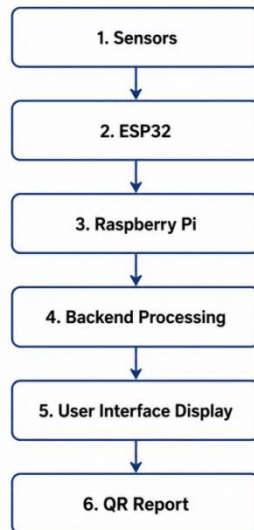
Communication Process Flow

Figure 6: Communication Process Flow of the System

That is why it is necessary to note that the sensors transmit data not to the display but to the ESP32 first, next to the Raspberry Pi, and from there to the touchscreen interface. This allows for a more structured and dependable system.

The sensors first send their information to the ESP32, which then sends it to the Raspberry Pi, and finally, the touchscreen interface. This provides a much better organized system compared to previous versions. This system uses named parameters for sending data, ensuring clarity and eliminating errors. For instance:

- Height: 164.0 cm
- Weight: 57.9 kg
- BMI: 21.5 kg/m²
- SpO₂: 98%
- Pulse: 78 BPM
- Respiration: 16 breaths/min
- Temperature: 36.3°C

In this way, the Raspberry Pi will know where to place each value. Such an approach is necessary to make sure the system can link the hardware measurement part to the software display/reporting part.

3.7 Feeding Data Using Sensors

Feeding data refers to collecting sensor readings and using the collected data for the system's operation. Feeding of data in case of the Basic Health Monitoring Station starts when the health check starts from the user end. Each of the sensors gathers only one kind of data and feeds it to ESP32.

Feeding Data Process:

3.7.1 Height Measurement Data

Height is determined by the help of VL53L0X laser range finder sensor. The sensor determines the distance between the sensor and the person's head. Considering the fact that the sensor is fitted at a certain height, the person's height can be calculated through deducting the distance detected by the sensor from the sensor's height.

For example:

$$\text{User Height} = \text{Fixed Sensor Height} - \text{Measured Distance}$$

This calculated height is then sent for BMI calculation. Proper positioning of the sensor is very important because height measurement depends on alignment between the sensor and the user's head.

3.7.2 Weight Measurement Data

The measurement of weight is done using four load cells attached at the bottom of the base platform. The instant the user stands on the base platform, the weight exerted will be sensed through the load cells. Since the output signal of load cells is quite small, the HX711 amplification module is employed to amplify the signals and obtain the digital value.

The ESP32 will then read the output signals of the HX711 module. Correction factors are introduced during the calibration process in order to ensure that the weight value measured is very close to the weight of the standard weight. After obtaining the weight, it is combined with height to determine the BMI.

The use of load cells in the project is very crucial since any errors in weighing will lead to erroneous BMI results.

3.7.3 SpO₂ and Pulse Rate Data

The MAX30102 sensor helps in measuring the SpO₂ and pulse rate levels. For measuring these, the user places their finger on the sensor, and the sensor then makes use of optical sensing to calculate changes in blood volume. It provides the pulse rate and SpO₂ values depending upon the measured values.

SpO₂ measures the blood oxygen saturation value in the body, and the pulse rate measures the heartbeat count per minute. These parameters play an important role in basic health check-ups as they give information related to the breathing and heart conditions.

The data measured by the MAX30102 sensor is then sent to the ESP32 device and finally processed using the Raspberry Pi. Calibration correction is done wherever required.

3.7.4 Respiratory Rate Data

The respiratory rate describes the breathing frequency in breaths per minute. This parameter has been incorporated into this study as one of the vital signs that can be assessed by the system. This respiratory rate is derived from the sensing and processing technique adopted in the system.

Once the value has been derived, it will be processed together with the other vital signs for display. Once processed, the data will be compared with the normal respiratory rate range to establish whether the patient's respiratory rate is low, high, or within normal range.

Respiratory rate has been included because it provides important information about a person's breathing condition.

3.7.5 Temperature Data

For measuring body temperature, an MLX90614 infrared sensor is used. By using this sensor, there is no need for physical contact with the device; therefore, it becomes more hygienic for multiple users.

In this sensor, body temperature can be calculated by measuring infrared radiation without any physical contact of the object. First, raw body temperature is calculated, and then the result is modified by the offset to get the final output temperature.

It was found during the calibration process that raw body temperature was less than that of reference temperature. So, to make the temperature measurement more accurate, an offset correction technique was introduced.

3.8 Back-end Processing in Raspberry Pi

The backend tasks are carried out primarily through the Raspberry Pi 4B. The Raspberry Pi receives the data sent by the ESP32, and thereafter the processor checks the data, corrects the data, arranges the data, and finally displays the data on the display screen.

The backend must perform the following functions:

- receive sensor readings from ESP32;
- divide the parameter from the received readings;
- check the correctness of the received values;
- calibrate the values of reading;
- calculate the BMI index based on the parameters of height and weight;
- compare the results with the norm;
- determine the overall health condition;
- send the processed values to the user interface;
- produce the health condition report QR-code.

The back end is important since it plays a vital role in converting the readings from the sensors into health readings. While it is correct that the sensor will simply provide a reading, it is also worth noting that the back end provides meaning to the readings. As an example, the back end does not only receive readings of height and weight but rather the one that helps figure out the BMI of the patient. It is also not simply about receiving the reading of SpO2 levels or body temperature but those that can tell whether the readings are normal or not.

3.9 BMI Calculation

The Body Mass Index will be automatically calculated using the height and body weights obtained. The Body Mass Index is a very important indicator to help determine whether one's body weight is in proportion to his/her height.

The formula used in computing BMI is:

$$\text{BMI} = \text{Weight (kg)} / [\text{Height (m)}]^2$$

In this case, the user's height would be determined by the VL53L0X sensor, while the user's weight would be detected by the load cells. The RPi uses these measurements to calculate the BMI of the individual according to the formula discussed above.

As an example, suppose that the weight of the individual is 58 kg and height is 1.64 meters,

$$\text{BMI} = 58 / (1.64 * 1.64)$$

$$\text{BMI} = 21.6 \text{ kg/m}^2$$

This information can then be displayed on the touch screen panel. BMI is part of the health report of the individual.

Calculating the BMI automatically reduces the burden of doing such calculations manually. Any possible error in the process is therefore avoided by automatic calculations.

3.10 Classification of Health Status Parameters

Once the readings have been obtained, the system categorizes them using easy to understand classifications. The main aim of this feature is to make it easier for the user to know whether the values recorded fall within acceptable parameters.

The classifications include:

- Low
- Normal
- High

The normal readings of oxygen levels would show normal readings while increased readings of body temperatures would indicate high temperature. The readings of body mass index would be high, normal, or low based on actual readings. It must be noted here that

the analysis of the health status does not have anything to do with the identification of any illness or disease. It serves as an indication for users to go to the doctor or hospital if there is something wrong with them.

3.11 User Interface Design

Interface development was done in a manner that makes it easy for non-experts to interact with at an academic health center and within the community. The health metrics recorded will be shown on the screen monitor, which is attached to the Raspberry Pi board. The user interface design has been done with the aim of achieving the following objectives:

- To display all key health data systematically
- To minimize chances of confusion during use
- To prevent manual recording of results
- To enable BMI calculation automatically
- To indicate health status easily
- To enable access to a digital health report through QR code

This interface will act as the primary interface through which the user interacts with the monitoring station. After completing the measurements, all the information will be displayed on the screen and therefore the user will have access to all his health information on one platform.

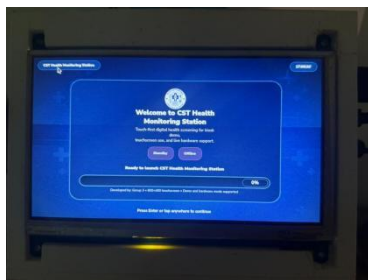


Figure 7: User Interface
Home Screen

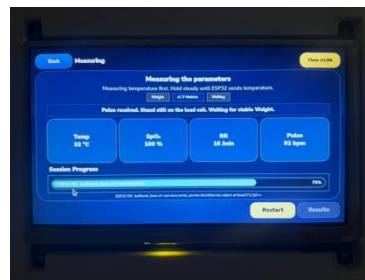


Figure 8: User Interface
Measurement Screen

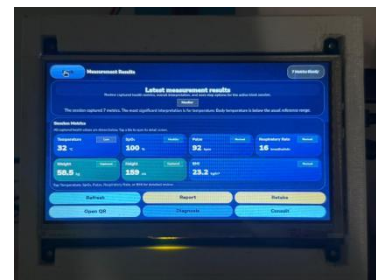


Figure 9: User Interface
Result Display Screen

3.12 QR Code and Health Report Creation

QR code based health report generation is yet another important feature of this system. After taking the measurement, the health report with the readings and the health status of the patient is generated by using Raspberry Pi.

The information included in the report is:

- Height;
- Weight;
- Body mass index;
- SpO₂;
- Pulse rate;
- Breathing rate;
- Body temperature;
- Health condition.

The flow chart of the QR code generation is presented as below:

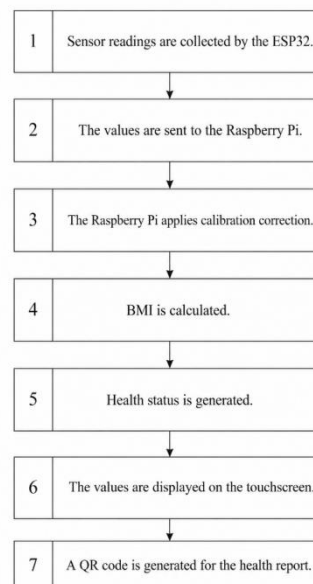


Figure 10: QR Code Process

This feature enhances the system's practicality and modernity through electronic documentation. Furthermore, it increases convenience, as results become available immediately upon completion of the measurement.

CHAPTER 4: MECHANICAL, DESIGN AND FABRICATION

4.1 Mechanical Design Overview

The mechanical body was designed in a way that ensured the strength and stability would be high enough when making different health measurements. This is because the design of

the mechanical body had to make sure the person remained upright during all the measurements, especially those concerning measuring his or her height and weight.

It was important to build the mechanical body in such a way that it would house the load cell platform, the height measurement rail, the sensor box, touch screen display, electrical components, and wiring system. It was important to make the base part of the mechanical body strong enough to hold load cells with room for height measurement.

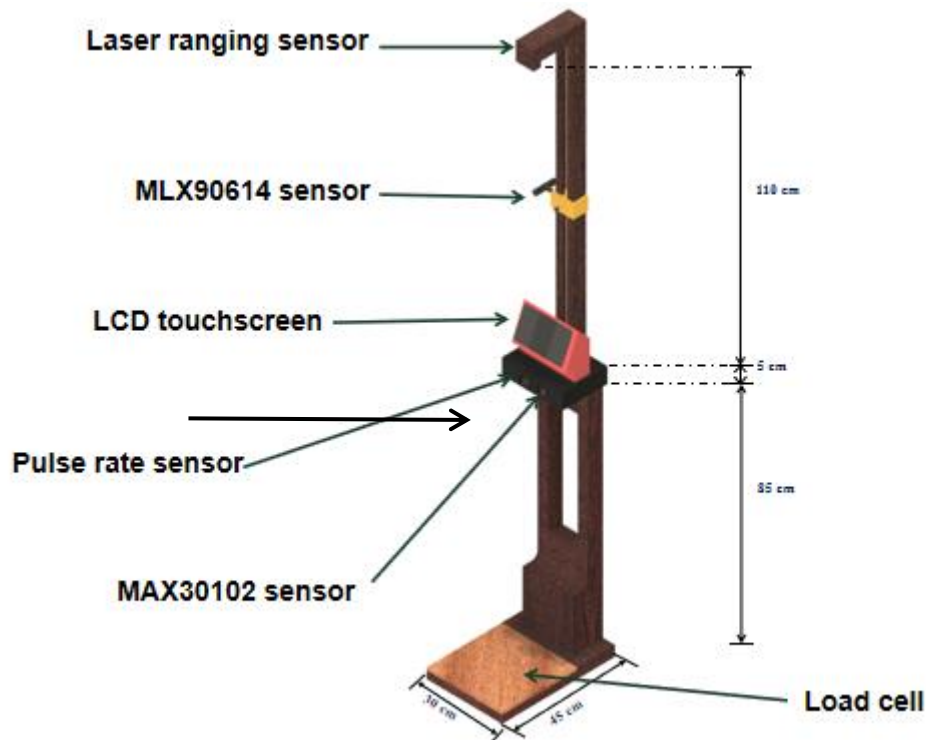


Figure 11: System Model

4.2 3D Model of the System Body

The motive for building the 3D model was to get an overall view of the entire set up before beginning the manufacturing process. Some of the view included the front view, side view, back view, top view, and bottom view. The 3D model helped get an understanding of the

base, vertical, sensor box, display stand, and the height sensor rail. The 3D model also gave an idea of the dimensions that would be required during the manufacturing process. The base was utilized to fix the load cell.

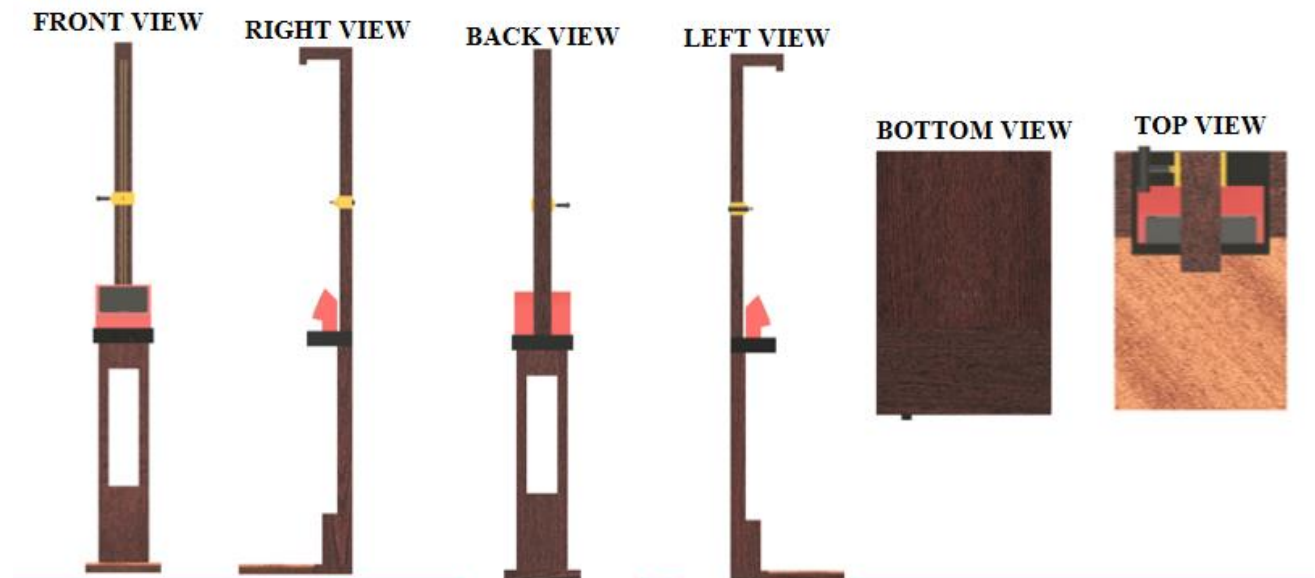


Figure 12: Views of 3D Design

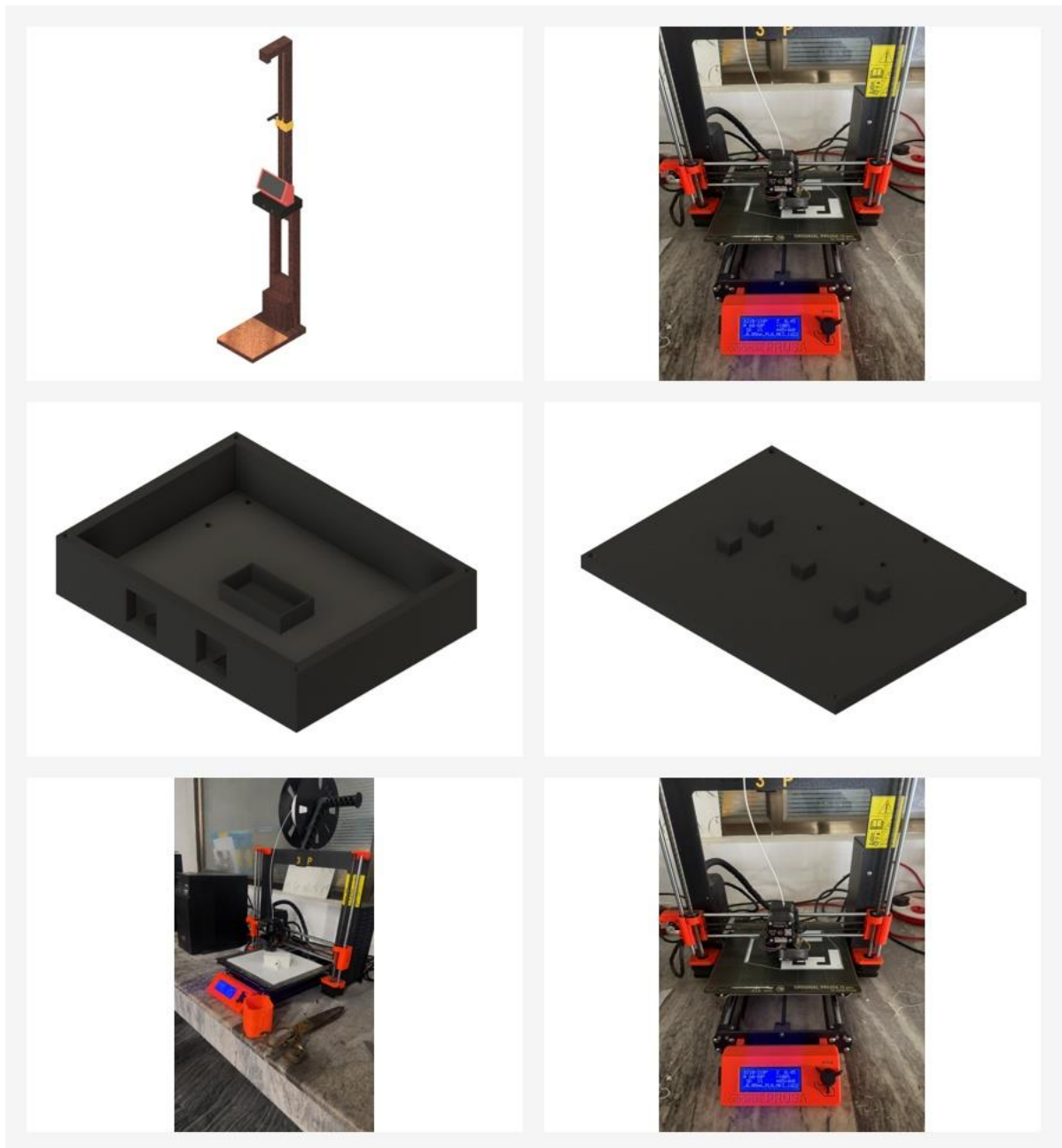
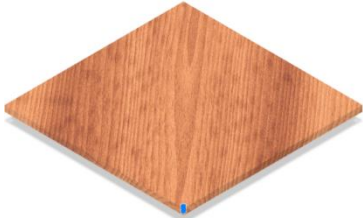
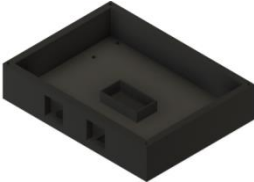
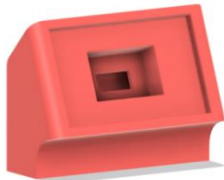
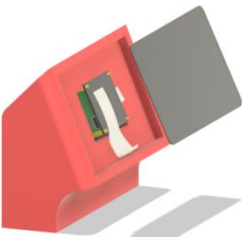
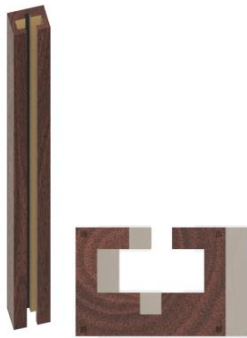
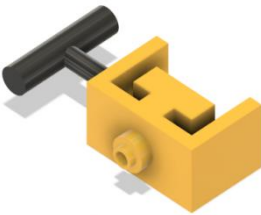
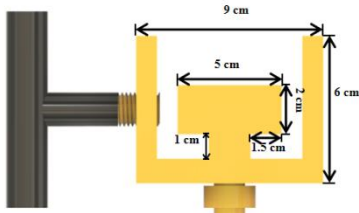
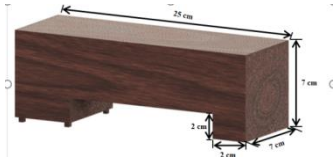


Figure 12(a): 3D-printed sensor box, protective covers, and prototype components fabricated for the BHMS body and sensor housing.

4.3 Dimensions of Components of the System Body

The body system is made up of the base plate having its cover, extended plate, mid box, sensor box having its cover, screen stand, the rail, the clamp, and the height box. All these were manufactured based on the positioning of sensors and the comfort of the user. The base plate is used to secure the load cell and provide support to the weight of the user. The extended plate enhances the stability of the system and prevents the user from touching the top portion of the system.

Design	Name	Dimension
	Base plate	30cm x 30cmx4cm
	Base plate cover	30cm x 30cm
	Extended plate	30cm x 15cm x 4cm
	Extended plate cover	30cm x 15cm
	Mid box	85 cm height
	Sensor box	24cm x 18cm x 6cm
	Sensor box cover	24cmx 18cm

 	Screen stand	22cm x 6cm base
	Rail	110cm height
   	Clamp	
	Height box	

The dimensions were chosen having in mind the positioning of the sensors, screen, load cell, and the standing position of the person being measured. This made it possible to choose the correct rail height, which would enable measuring the height of the person.

4.4 Carpentry Work

The Prototype of the Basic Health Monitoring Station was mostly constructed from wooden materials due to their cost-effectiveness, availability, easy manufacturing process, and durability. Wood was used to construct the base, the main body of the machine, the load cell body, and the screen body. The manufacturing process involves cutting, making, assembling, and finally finishing to get the desired sizes. A lot of consideration was taken into account in designing the base body to ensure that the weight of the individual is supported for proper measurements. The load cell cover was designed to protect the load cells.



Figure 13: Woodwork

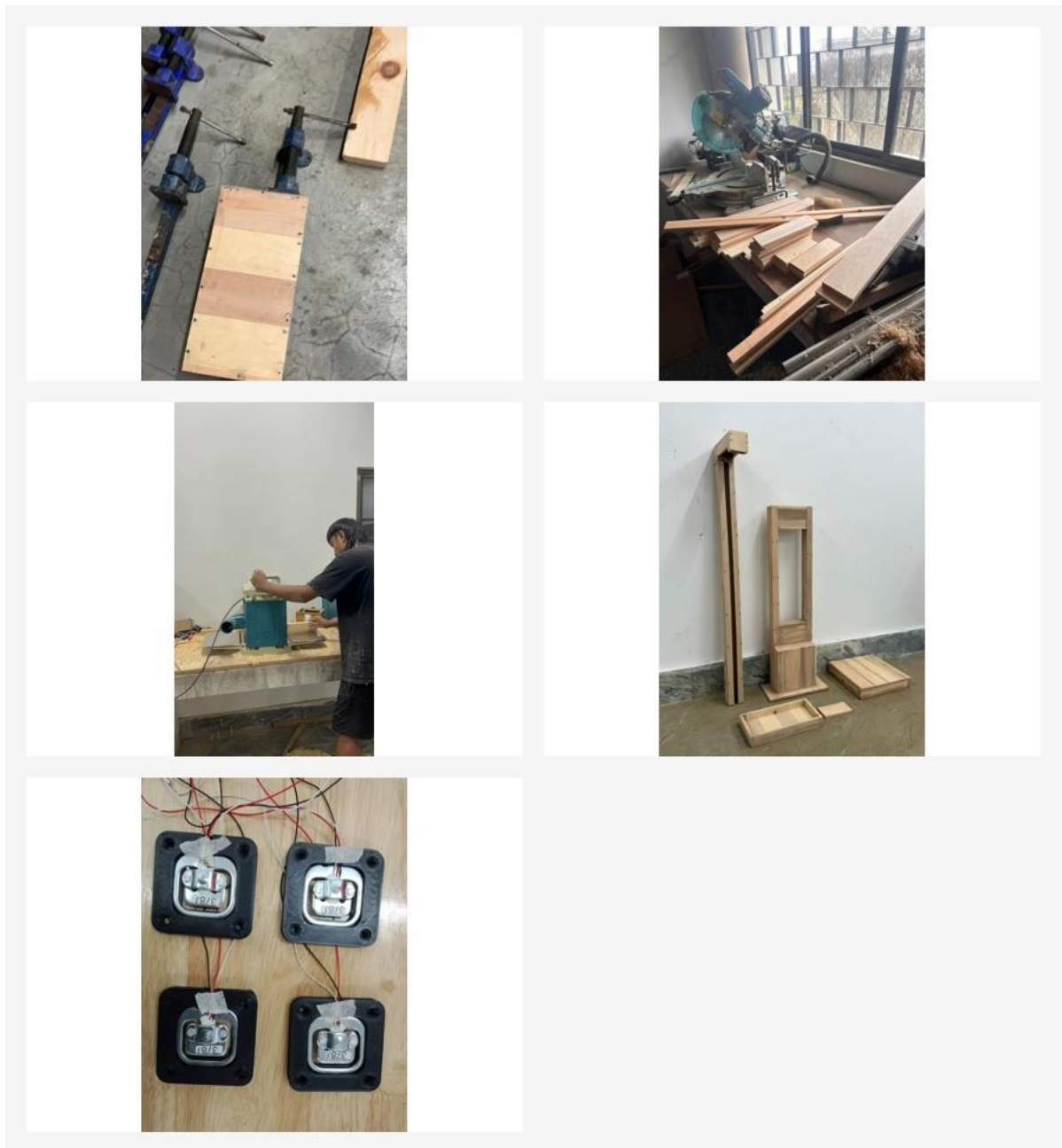


Figure 13(a): Carpentry work, wooden-frame preparation, and fabrication of the main structural parts of the BHMS prototype.

4.5 Modification of Body Frame

Once it had been installed, there was an issue that arose regarding the placement of the sensor box and the LCD, whereby they were placed too close to the body of the individual using the product. This could lead to a situation whereby someone would accidentally touch the top when standing on the base, especially individuals with larger builds. This led to discomfort, as well as instability, as touching the top meant that the sensors would be affected, particularly the height sensors.

4.7 Assembly of the System

Assembly of the wooden frame and 3D components occurred after the completion of all manufacturing processes. Firstly, four load cells were mounted on the base plate. Then, they were attached to the HX711 module. Secondly, the laser rangefinder was mounted at the height measurement module position. Thirdly, MAX30102, pulse rate sensor, and temperature sensor were put into the sensor box.

Connection of ESP32 with sensors enabled data monitoring online. The Raspberry Pi 4B handled computation tasks, GUI creation, mounting of the touch screen LCD, and printing of health reports. Mounting of the touch screen onto the holder was done. All wires were hidden in the body. Power wires powered all components of the project. Testing was done to ensure that all hardware and software worked properly.

Consequently, the final model could be used to calculate a patient's height, weight, BMI, heart rate, respiration rate, SpO₂ level, and body temperature. The information appeared on the display screen alongside QR codes.

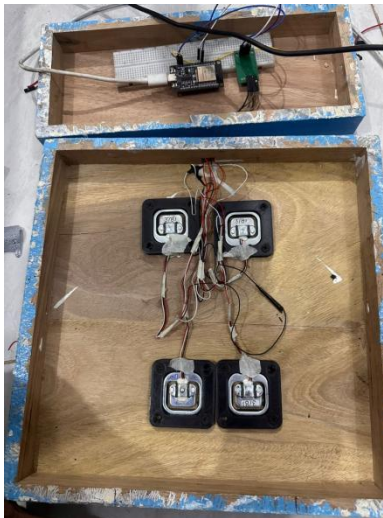


Figure 14: Integration of Load Cells and HX711

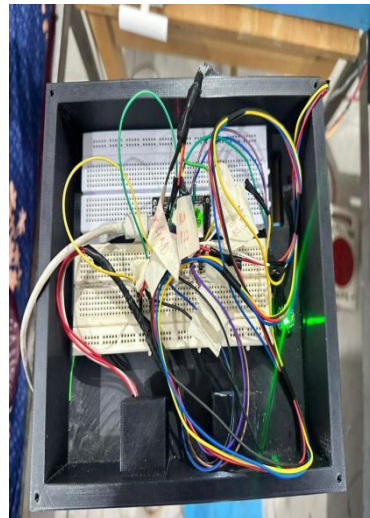


Figure 15: Integration of Sensors in Sensor Box

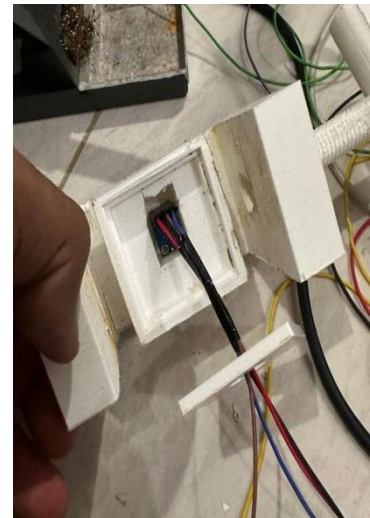


Figure 16: Mounting the temperature sensor



Figure 17: Final Assembly of Basic Health Monitoring Station

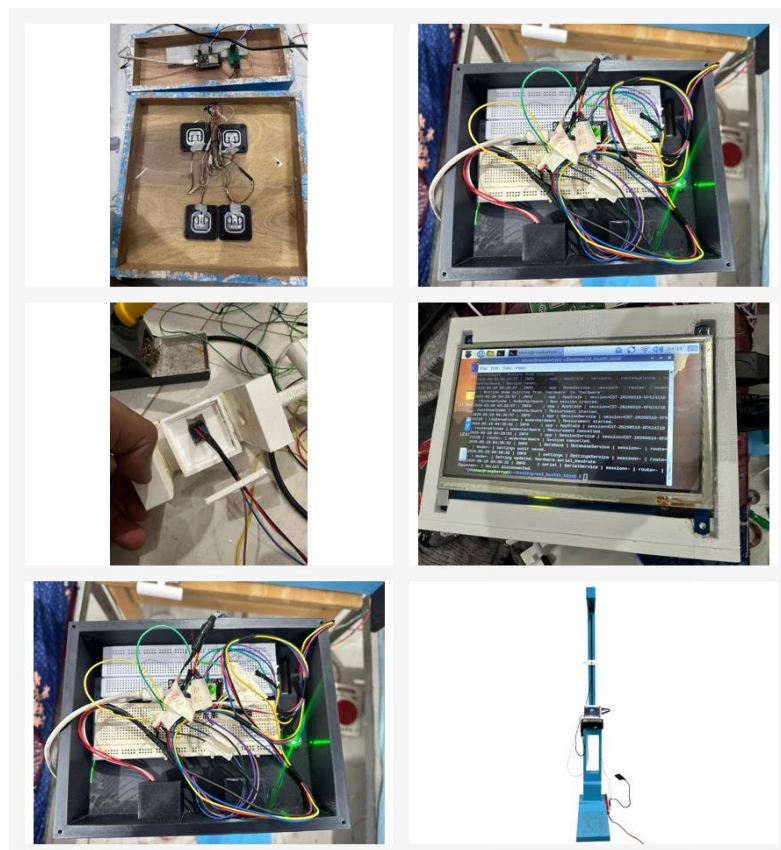


Figure 17(a): Assembly activities showing load-cell mounting, internal wiring, touchscreen integration, and the completed upright prototype.

4.7 Material Selection

The mechanical body was primarily made using wood while the sensors were located inside structures that were produced through 3D printing technology. Wood provided a solid structure for the body and made it stronger, while 3D printing allowed creating a perfectly fitting housing for the sensors. Using these materials made it possible to reduce costs without sacrificing functionality.

Materials used in the prototype were chosen taking into account the following aspects:

- strength to bear the weight of the user
- simplicity in manufacturing
- local availability
- cost effectiveness
- simplicity in modification
- appropriateness for prototype building

The selection of materials was appropriate for making the prototype. But when it comes to future improvements, the mechanical body could be made stronger with the use of metals or any other enclosure material.

4.8 Planning for the College Infirmary and User Guidance

The section describes the installation setting, space layout, and strategy for guidance while installing the Basic Health Monitoring Station within the college infirmary. The installation and usage of the prototype took place after it was successfully tested and calibrated through the tests and calibrations discussed in Chapter 5.

The chosen site for installation was the college infirmary since this setting provided an appropriate environment for health checkups of the individuals. The intended users of the system were adult individuals within the age group of 19 to 35 years.

Guidance labels and on-screen instructions for correct posture, finger positioning, and order of measurements were created as part of the preparation process to enable systematic introduction of the kiosk in the infirmary once the validation phase was complete. Preparation activities allowed for a better understanding of how the kiosk was to be placed and controlled in the infirmary once the validation phase was complete.

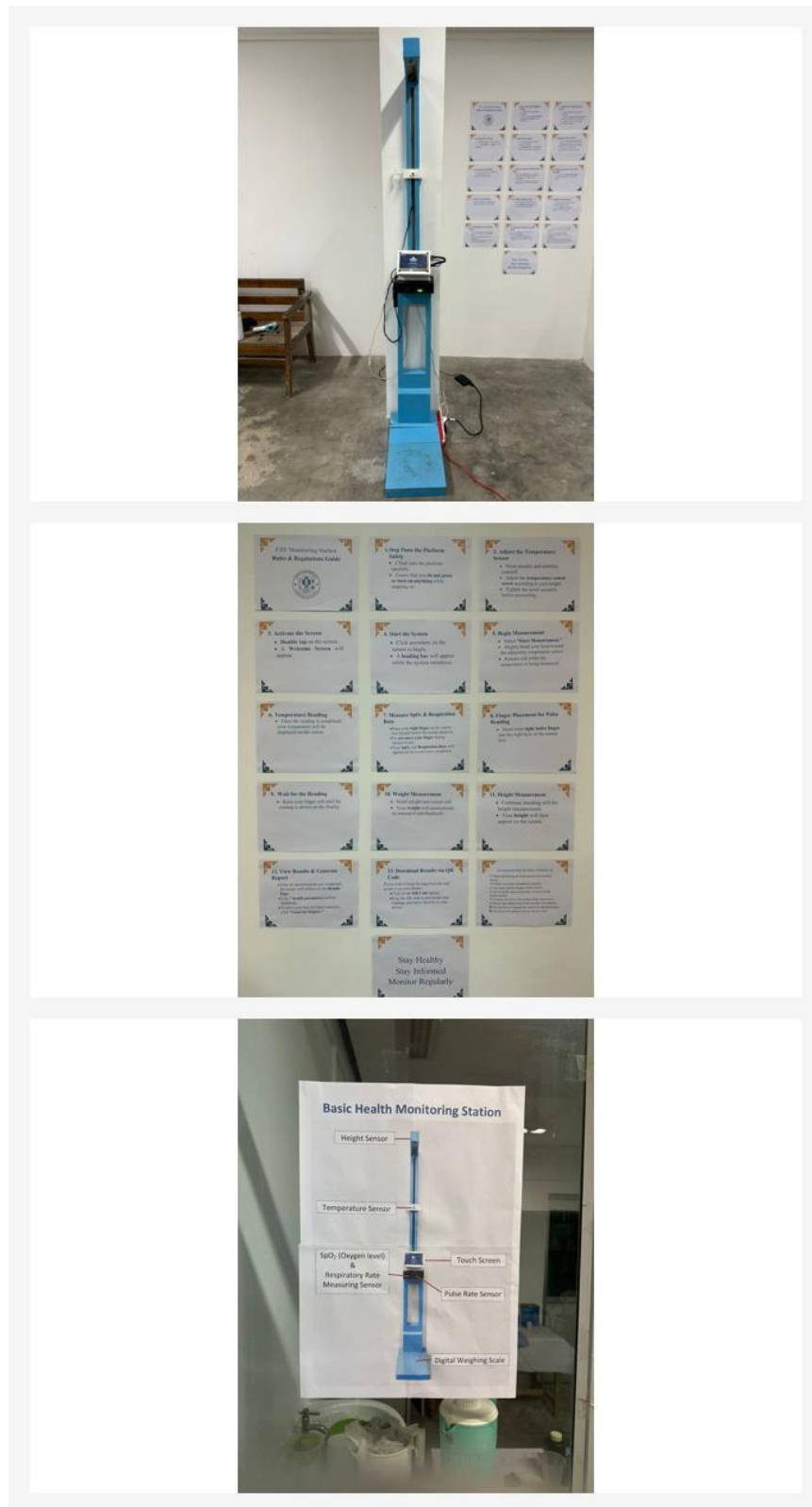


Figure 17(b): Prototype placement in the infirmary area together with the user-guidance poster and installation environment.

Figure 17.1 User Guideline

CHAPTER 5: TESTING AND CALIBRATION

The process of testing and calibration became critical when developing the Basic Health Monitoring Station, as it is necessary to measure quite a few health indicators. Therefore, before assembling, every single sensor was tested and calibrated. Testing ensures that everything works fine – sensors, microcontroller, Raspberry Pi, monitor, and interface included.

Then, following the testing process, each sensor was calibrated to improve the accuracy of the readings taken. The offset, or the difference between the true reading and the measured value, could be defined.

5.1 Individual component testing

Before joining the entire system together, tests for individual hardware were done to find any faults in connection, sensors, and the software interaction. The important things tested included:



Figure 18: Testing of ESP32



Figure 19: Testing of LCD and Raspberry pi 4 with Raspberry pi

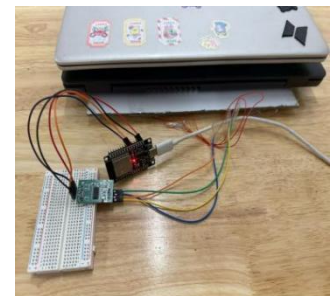


Figure 20: Testing of HX711 with load cell

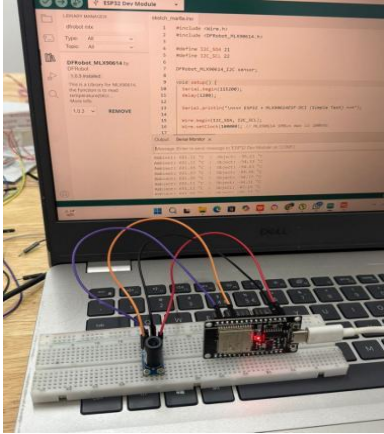


Figure 21: Testing of
MLX90614ESF-DCI

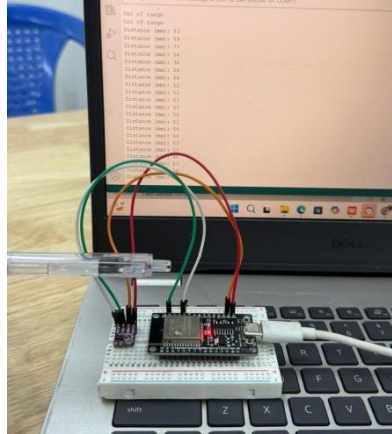


Figure 22: Testing of
VLX5L0X

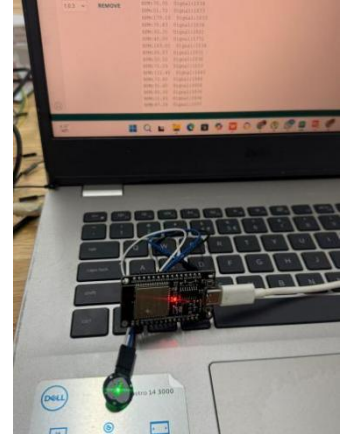


Figure 23: Testing of
Pulse rate sensor

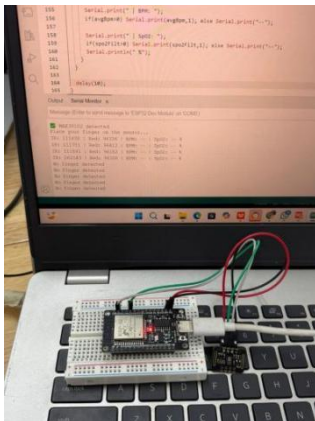


Figure 24: Testing of
MAX30120

The test was carried out on ESP32 to find out whether the microcontroller can be able to obtain data from the sensors. Testing of the Raspberry Pi was performed by installation of the operating system and the usage of the Graphical User Interface. Touchscreen display test was performed to ascertain proper display of the output data. HX711 and load cell testing was conducted to verify weight measurement functionality. On the other hand, testing of oxygen sat

The test was carried out on ESP32 to find out whether the microcontroller can be able to obtain data from the sensors. Testing of the Raspberry Pi was performed by installation of the operating system and the usage of the Graphical User Interface. Touchscreen display test was performed to ascertain proper display of the output

data. HX711 and load cell testing was conducted to verify weight measurement functionality. On the other hand, testing of oxygen saturation (SpO₂), heart rate, and body temperature measurements was performed using MAX30102 and MLX90614 sensors respectively.uration (SpO₂), heart rate, and body temperature measurements was performed using MAX30102 and MLX90614 sensors respectively.

Calibration of sensors

1. MAX30102

Calibrating:



```

Message (Enter to send message to 'ESP32 Dev Module on 'COM4')
Sensor SpO2 = 88.8
Enter Medical Device SpO2 value:
Medical Entered = 0.00
Sample 12/15

Sensor SpO2 = 84.8
Enter Medical Device SpO2 value:
Medical Entered = 89.00
Sample 13/15

Sensor SpO2 = 84.8
Enter Medical Device SpO2 value:
Medical Entered = 0.00
Sample 14/15

Sensor SpO2 = 84.8
Enter Medical Device SpO2 value:
Medical Entered = 91.00
Sample 15/15

===== CALIBRATION RESULT =====
Slope (A): -1.6955
Offset (B): 197.2204

Use these in your main code:
const float SpO2_A = -1.6955;
const float SpO2_B = 197.2204;
  
```

After Calibration:



```

No finger detected
No finger detected
No finger detected
No finger detected
No finger detected
No finger detected
MAX30102 Pulse Oximeter
Place finger on sensor
Finger detected
BPM: 63.4 | SpO2: 97.0 %
BPM: 63.4 | SpO2: 96.7 %
BPM: 63.4 | SpO2: 97.5 %
BPM: 84.5 | SpO2: 97.5 %
BPM: 99.1 | SpO2: 97.2 %
BPM: 97.5 | SpO2: 97.1 %
BPM: 87.8 | SpO2: 96.9 %
No finger detected
  
```

Blood oxygen saturation level was determined using the MAX30102 sensor. In the process of calibration, the readings of SpO₂ from the MAX30102 sensor were compared with the readings of the reference device. It was noted that the sensor readings were a little bit smaller compared to those of the reference device. Thus, an offset was determined.

The offset was calculated using the following formula:

$$\text{Offset} = \text{Raw Sensor Reading} - \text{Medical Device Reading}$$

Table 5.1: MAX30102 SpO₂ Calibration Before Correction

Reading	Medical Device SpO ₂ (%)	MAX30102 Raw SpO ₂ (%)	Offset (%)
1	98	96.31	-1.69
2	97	95.42	-1.58
3	99	97.26	-1.74
4	96	94.28	-1.72
5	98	96.37	-1.63
6	97	95.29	-1.71
7	99	97.41	-1.59
8	98	96.25	-1.75
9	97	95.36	-1.64
10	96	94.35	-1.65

Offset for the readings provided by the MAX30102 SpO₂ sensor was around -1.67%. The raw data being smaller than that of the medical device, the correction required is +1.67%.

Table 5.2: MAX30102 SpO₂ Calibration After Correction

Reading	Medical Device SpO ₂ (%)	MAX30102 Raw (%)	Correction Applied (%)	Corrected MAX30102 (%)
1	98	96.46	+1.67	98.13
2	97	95.21	+1.67	96.88
3	99	97.34	+1.67	99.01
4	96	94.41	+1.67	96.08
5	98	96.18	+1.67	97.85
6	97	95.47	+1.67	97.14
7	99	97.12	+1.67	98.79
8	98	96.52	+1.67	98.19
9	97	95.33	+1.67	97.00
10	96	94.26	+1.67	95.93

When the correction factor was applied, the MAX30102 measurements became more in line with the reference medical device. It made the measurement of the SpO₂ more accurate.

2. MLX90614ESF-DCI

Calibration of the MLX90614 temperature sensor was done through the analysis of the relationship between the raw values recorded by the sensor and the actual temperature. It was observed that during the five trials carried out, the readings taken were slightly lower than the actual temperature.

The correction factor was +2.30°C.

Corrected temperature = Temperature raw value + 2.30°C

Following correction, the readings taken were close to the actual temperature.

Table 5.3: MLX90614 Temperature Calibration

Reaadings	Actual Temperature (°C)	Sensor Raw Temperature (°C)	Offset Applied (°C)	Corrected Tempera ture (°C)
1	36.30	34.02	+2.30	36.32
2	36.10	33.79	+2.30	36.09
3	35.90	33.61	+2.30	35.91
4	36.40	34.08	+2.30	36.38
5	36.20	33.95	+2.30	36.25

3. Pulse Rate Sensor

Calibration of the pulse rate sensor involved checking the BPM information from the sensor against the pulse rate readings. There were five readings of the information in the sensor in the same conditions as compared to find the difference between the two readings.

The results revealed that there was high BPM reading from the sensor as opposed to actual readings. An offset correction of -0.30 BPM was used to balance the difference between the two sets of readings.

Corrected Pulse Rate = Sensor Raw Pulse Rate - 0.30 BPM

As a result of offset correction, the pulse rates were near to the actual readings.

Table 5.4: Pulse Rate Calibration

Readings	Actual Pulse Rate (BPM)	Sensor Raw Pulse Rate (BPM)	Offset Applied (BPM)	Corrected Pulse Rate (BPM)
1	77.0	78.4	-0.30	78.1
2	78.0	76.9	-0.30	76.6
3	76.0	76.8	-0.30	76.5
4	79.0	79.6	-0.30	79.3
5	75.0	75.2	-0.30	74.9

4. VLX53L0X

The calibration of VL53L0X range finder sensor was done by taking raw values and then comparing them with the actual values of distance, using six different data sets. There were some variations observed in the raw values, and thus the offset value, which is nearly equal to -40 mm, was computed as follows:

Calibrated Distance = Sensor Raw Distance - 40 mm

Using the calibrated values, the height in cm was found.

Table 5.5: VL53L0X Height Sensor Calibration

Readings	Sensor Raw Distance (mm)	Calibrated Distance (mm)	Offset Applied (mm)	Calculated Person Height (cm)
1	404.00	363.55	-40.45	163.6
2	405.00	364.45	-40.55	163.6
3	405.60	364.99	-40.61	163.5
4	392.60	353.30	-39.30	164.7
5	395.80	356.18	-39.62	164.4
6	400.40	360.31	-40.09	164.0

5. HX711 with Load Cells

Load cell calibration was carried out at Druk Wang by applying loads starting from 20kg up to 200kg using increments of 20kg. At the start, before loading anything on the load cell, the load cell was first set to zero in order to show the load reading as 0kg. There was an offset observed in the reading of load cell as compared to the actual loading. The offset is given by the subtraction of the actual load from the load cell reading. The offset came out as +2.01kg after analyzing ten samples with 1.83% error.

Table 5.6: Load Cell Calibration Using Standard Weights

Reading No.	Standard Weight Applied, kg	Raw Load Cell Reading, kg	Offset / Correction, kg	Corrected Reading, kg
0	0	0.00	0.00	0.00
1	20	19.62	+0.38	20.00
2	40	39.28	+0.72	40.00
3	60	58.93	+1.07	60.00
4	80	78.55	+1.45	80.00
5	100	98.21	+1.79	100.00
6	120	117.80	+2.20	120.00
7	140	137.42	+2.58	140.00
8	160	157.10	+2.90	160.00
9	180	176.65	+3.35	180.00
10	200	196.30	+3.70	200.00

5.2 Integrated testing

Once all sensors were calibrated, they were connected together to complete the Basic Health Monitoring Station. Test runs for the complete system revealed that all sensors, along with the ESP32, Raspberry Pi, touchscreen, and interface, functioned properly.

When the components of the unit were assembled and tested, a person was positioned on the platform in order to determine the body weight through load cells and body height with the use of the VL53L0X sensor. SpO2 and pulse rate readings were determined through MAX30102 sensor and pulse sensor, respectively.

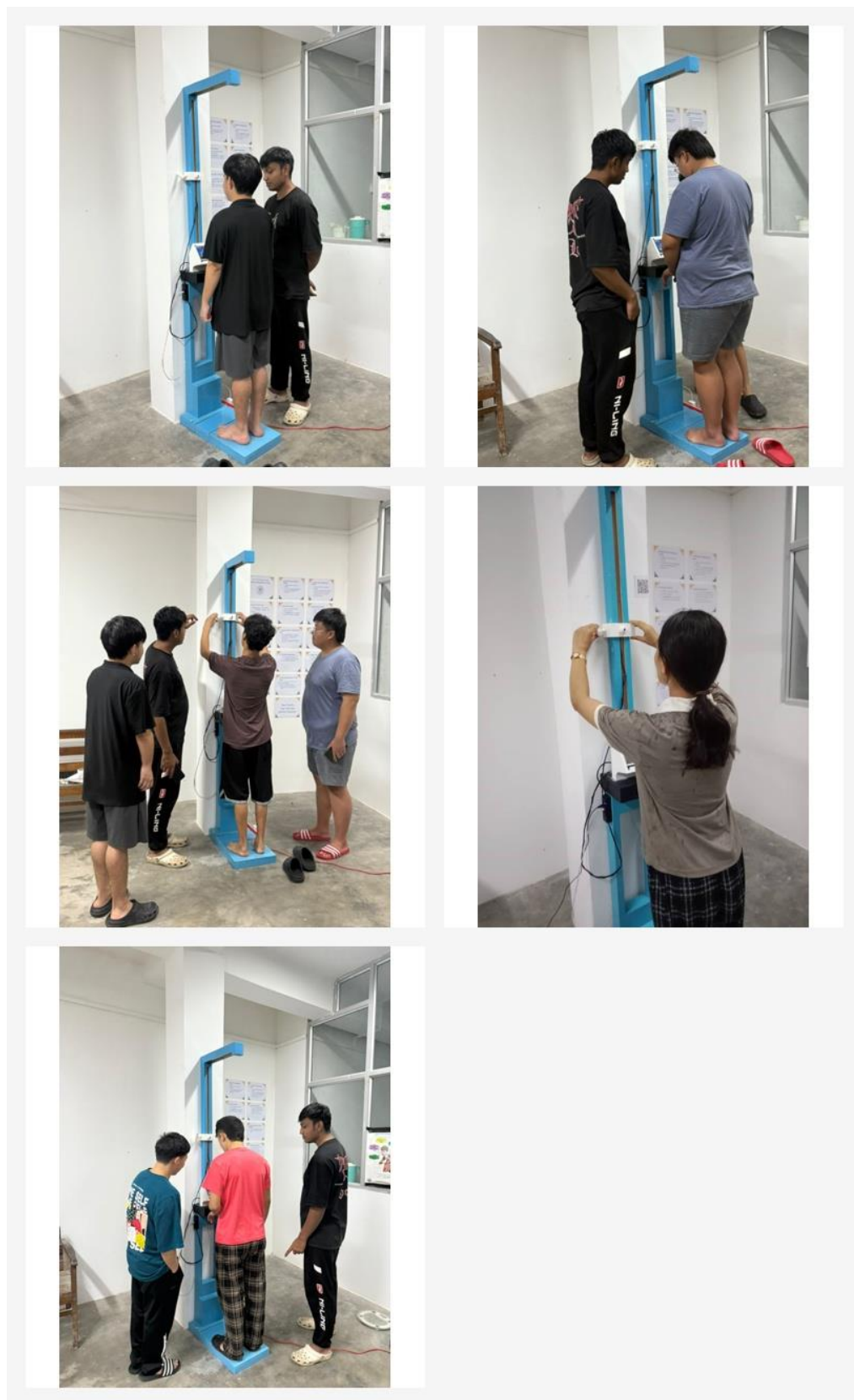


Figure 27(a): Actual user-testing photographs showing different participants using the BHMS for measurement, observation, and height-check procedures.

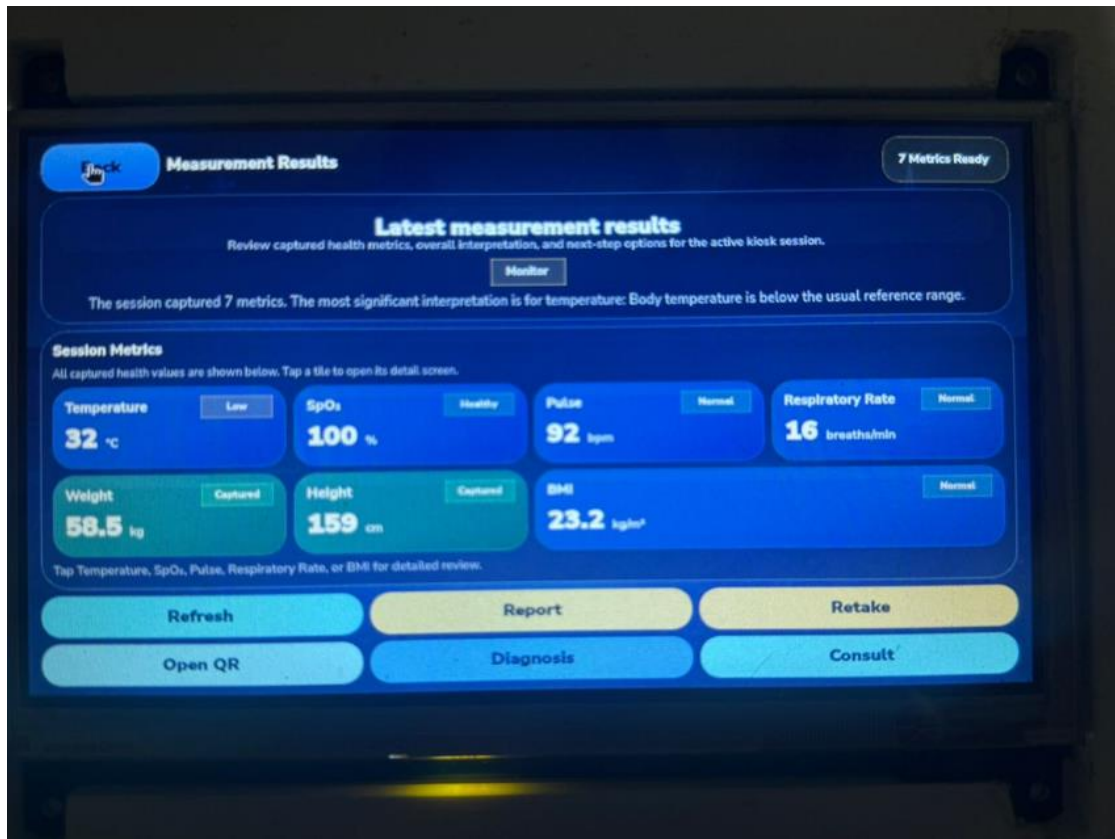


Figure 27(b): Measurement-results screen captured during integrated system testing of the complete BHMS application.

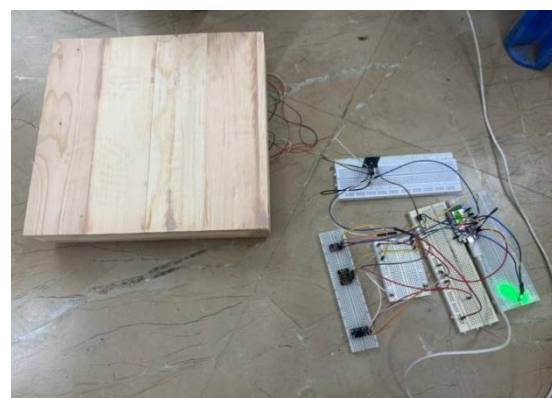
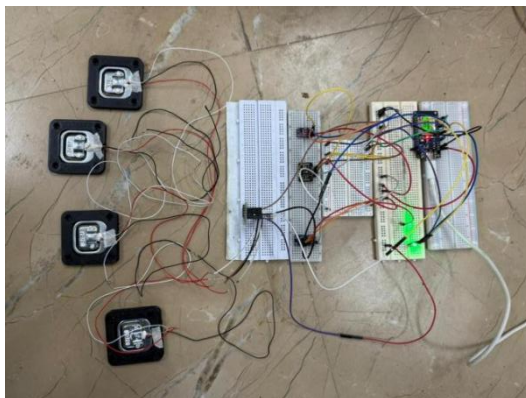


Figure 25: Integrated System Testing

5.4 Uploading of APP to the Capacitive Touch Display Screen.

After developing the health monitoring application, it was uploaded and installed into the Raspberry Pi that was connected to the 7-inch capacitive touch display screen. After installing the application into the Raspberry Pi, all necessary Python libraries and PyQt dependencies were installed to make sure that everything would work properly with regards to the graphical user interface. Once installation was completed successfully, the application was executed in the touchscreen and checked to confirm its suitability in terms of screen usage, touch capabilities, navigation through buttons, and alignment of the display screen. The successful upload of the application ensured that everything worked as expected such as the welcome screen, instructions screen, measurement screen, result screen, diagnosis screen, QR report, and administration screens.

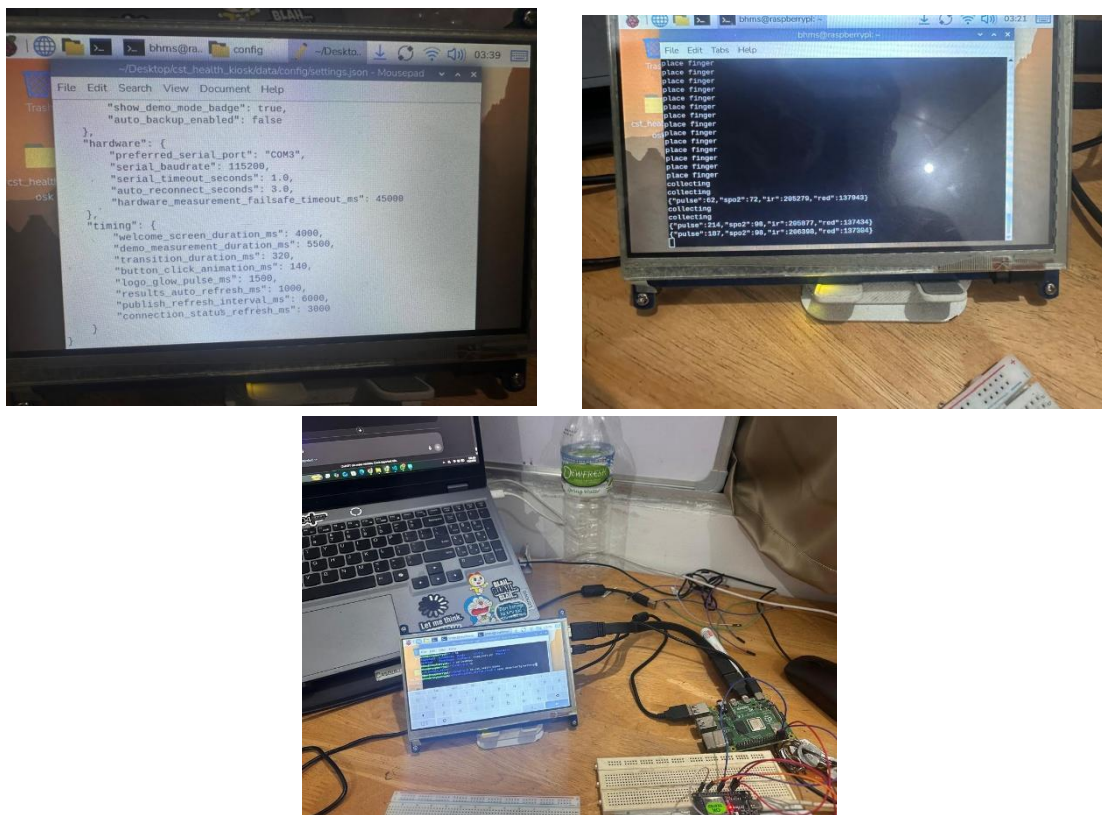


Figure 26.1: Uploading of APPLICATIONS

CHAPTER 6: SOFTWARE DESIGN AND IMPLEMENTATION

6.1 Introduction

This chapter discusses the software design and development processes of building the Basic Health Monitoring Station (BHMS). These include the ESP32 sensor data acquisition module, the Raspberry Pi processing module, touchscreen graphical user interface, QR code report generator, voice guidance, and storage system. The software is responsible for converting sensor data from the sensors into understandable health information that can be understood by non-technical college clinic patients.

ESP32 manages all the real-time sensor data gathering, and the Raspberry Pi 4B handles the graphical user interface, computation, health diagnosis, QR code report generator, and storage. The touchscreen graphical user interface will be designed for a 7-inch capacitive touch screen display, which will give detailed user guidelines using minimal supervisions.

Various software packages were utilized in designing and developing BHMS. Embedded programming, designing the applications on the Python, and graphical user interface design contributed significantly to making the BHMS into a health station.

6.2 ESP32 and Raspberry Pi Programming

The program for ESP32 is designed to receive data provided by various biomedical and sensing devices. Data acquisition is done via ESP32 microcontroller, including acquiring weight with the help of the HX711 amplifier, distance using the VL53L0X sensor, temperature using MLX90614 sensor, and also obtaining information about the pulse (number of pulses, SpO2%, respiratory frequency) using MAX30102 and pulse sensors. After preprocessing the received data, ESP32 will send packets of information with labels enabling the Raspberry Pi to correctly classify the data and further perform actions based on them. ESP32 maintains the serial communication with the Raspberry Pi and implements stepwise measurements in accordance with the application algorithm.

Raspberry Pi 4B computes BMI, classifies health status, updates GUI screens, generates reports and QR codes, and also stores all data locally. For this computer, Python with PyQt graphical library was used for implementing the software. Its functionality included transition of screens, processing of buttons, creation of the diagnosis report, generation of the report, printing QR codes, sound triggering, and administrative actions.

6.3 GUI Design

The following screen shots will be able to help in illustrating all of the functions available in the graphical user interface more effectively. Functions like navigation, report storage, configuration, administrator protection, interpretation, diagnosis, and consulting are just some of the functions that have been shown below.

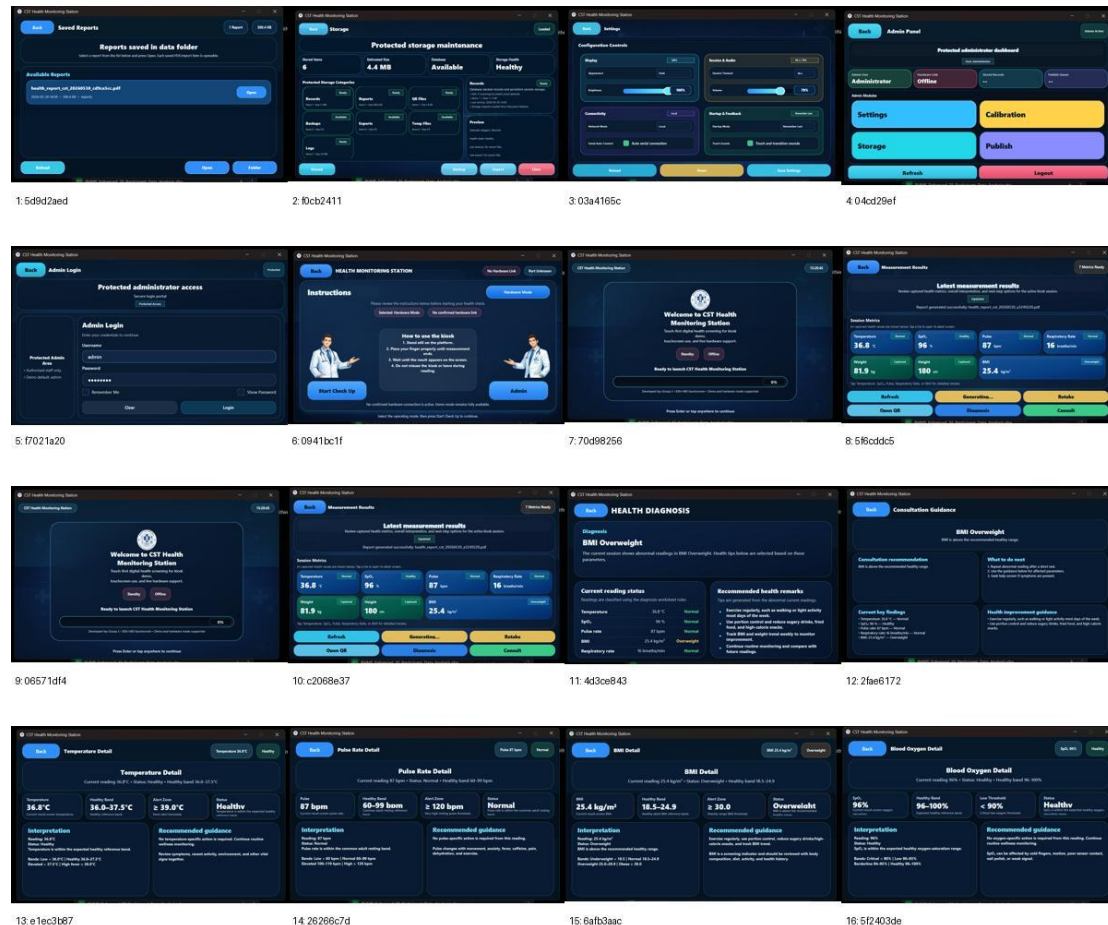


Figure 6.1: Additional GUI screen gallery showing the main operational and administrative screens of the BHMS application.

GUI had been designed in such a way that it could easily be operated by any user without the need for any technical skills. Large buttons, label names, separated metric cards, and indicator lights have been incorporated into the GUI design. Below are some screenshot examples of the GUI windows developed for BHMS application.

Sequential Screen Flow of the Final BHMS Application

To provide clarity in the final report, the screens of the application have been included below in the order that the user encounters them during the process. This ordering has been done according to the final developed GUI interface and the contexts used in the PowerPoint presentation. The sequence starts from the welcome page, proceeds to measurement guidance and results interpretation, and ends with report access and administration screens.

The aim of showcasing the sequence of the screens is to show that BHMS is not just a sensor device but also an entire system for guided digital screening. Each individual screen plays a role in the process, ranging from user onboarding to measurement management, results presentation, health interpretation, consultation management, report generation, QR code access, and system administration.

Step	Screen	Main Purpose
1	Welcome Screen	Launches the application and introduces the BHMS interface.
2	Instruction and Mode Selection Screen	Allows the user to choose Start Check Up or Admin and read operating instructions.
3	Measurement Result Screen	Displays the seven measured parameters and action buttons.
4	Health Diagnosis Screen	Shows the overall diagnosis and health remarks.
5	Consultation Guidance Screen	Provides next-step advice and consultation guidance.
6	Parameter Detail Screens	Shows detailed interpretation for temperature, SpO2, pulse and BMI.
7	PDF Report and QR Viewer	Makes the result available as a printable report and phone-based viewer.
8	Administrative Screens	Supports login, settings, storage and saved-report management.

Table 6.A: Sequential flow of the main BHMS application screens.

Step 1: Welcome Screen

The welcome screen will be displayed when the user runs the kiosk software, being the first screen presented to the user. The welcome screen presents the CST Health Monitoring Station and the loading bar along with the stand-by/offline indicators. Furthermore, this screen will ask the user to continue from the beginning.

In actual practice, the purpose of the welcome screen is also to provide confirmation that the system has run correctly. Moreover, it works in collaboration with the voice assistance module that says, “Welcome to CST Basic Health Monitoring Station.”

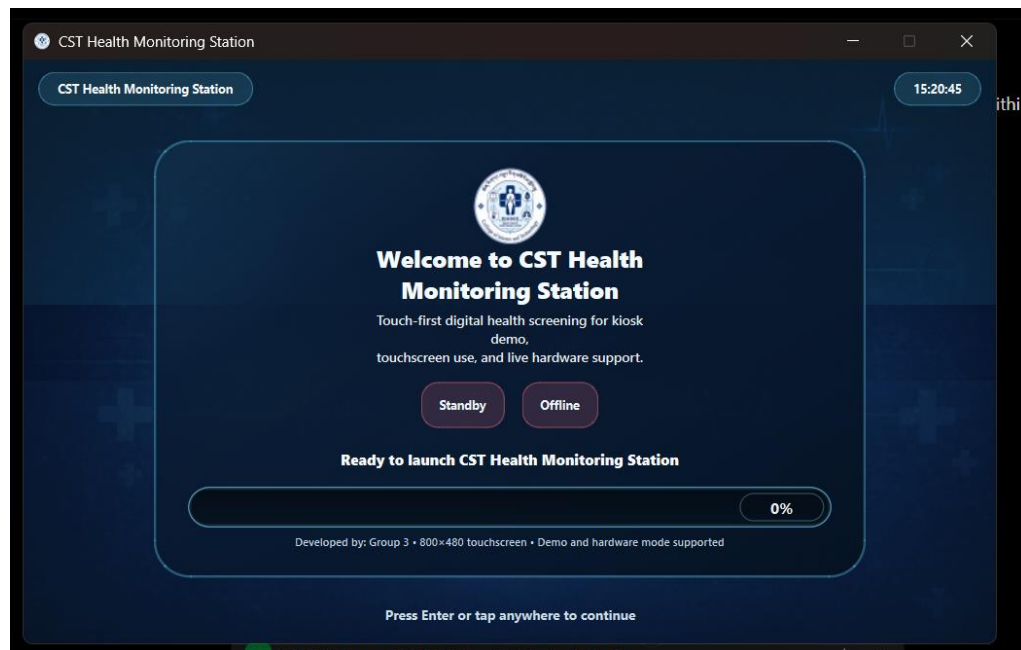


Figure Welcome Screen: Welcome Screen used in the final BHMS application.

Button / Area	Functionality
Standby / Offline	Displays the operating status of the kiosk.
Progress bar	Shows readiness of the application at startup.
Tap anywhere / Press Enter	Moves the user to the next screen.

Table Welcome Screen: Main functions of the welcome screen.

Step 2: Instruction and Mode Selection Screen

Once past the welcome page, the user arrives at the instruction and mode page. Here, the user is guided with basic instructions on operating the device and can select between the two modes available – Start Check Up and Admin modes.

The importance of this page lies in the fact that it differentiates between user mode and admin mode. Furthermore, it tells the user about the correct standing position and

placement of finger while performing the measurements. Here, the measurement process is almost about to commence and the voice guide directs the user towards the sensor action.

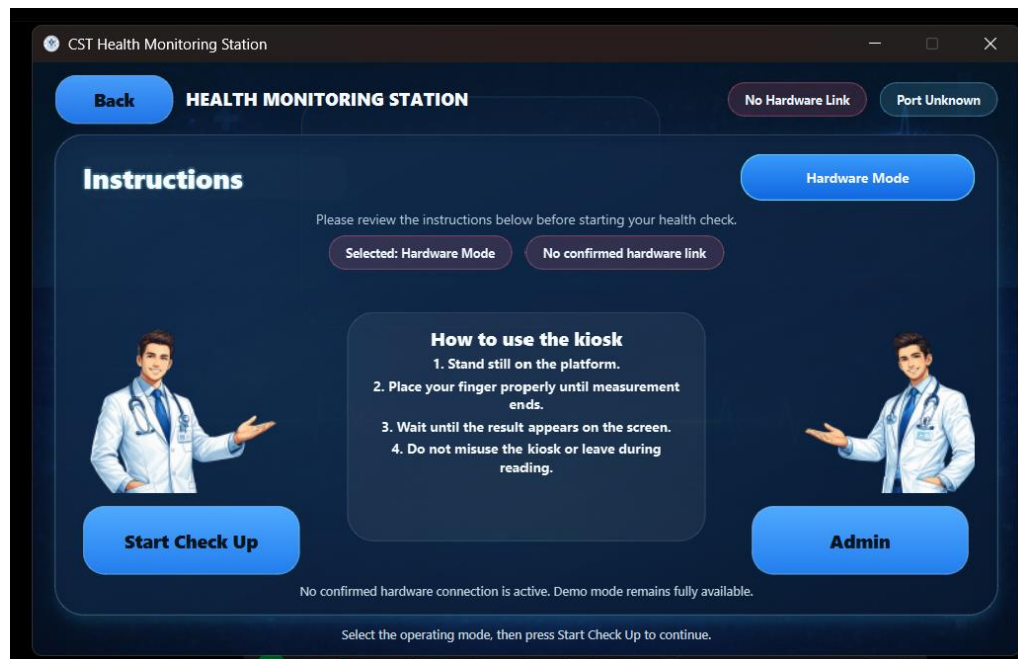


Figure Instruction and Mode Selection Screen: Instruction and Mode Selection Screen used in the final BHMS application.

Button / Area	Functionality
Start Check Up	Starts the guided health screening process.
Admin	Opens the administrator login screen.
Hardware Mode	Indicates hardware-based measurement mode.
Back	Returns to the previous screen.

Table Instruction and Mode Selection Screen: Main functions of the instruction and mode selection screen.

Step 3: Measurement Result Screen

After all measurements have been taken, the result page shows the seven parameters that are measured, which include temperature, SpO₂, pulse rate, breathing rate, weight, height, and BMI. This is the main result page of the BHMS app.

The layout of the card style helps make the result page more readable, especially when viewed on a touch screen. The user can refresh, redo the test, create a report, get the QR code, see the diagnosis, or get consultation advice. This is why the result page functions as an intermediary between raw data and meaningful health data.

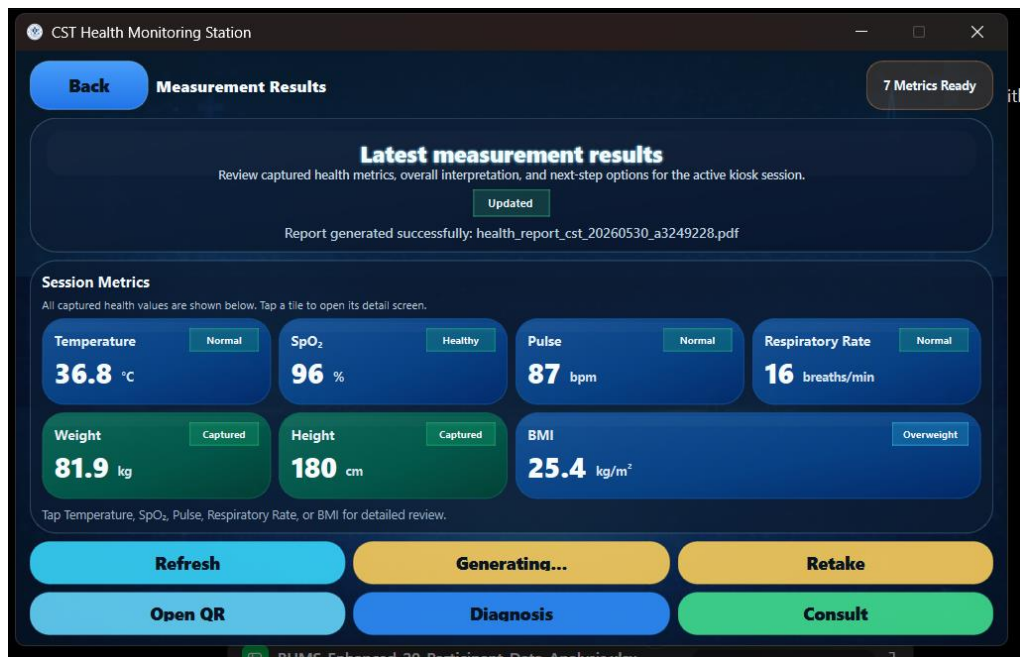


Figure Measurement Result Screen: Measurement Result Screen used in the final BHMS application.

Button / Area	Functionality
Refresh	Reloads the current result view.
Report	Generates the PDF health report.
Retake	Starts the measurement again.
Open QR	Opens the QR code for mobile report access.
Diagnosis	Opens the diagnosis screen.
Consult	Opens the consultation guidance screen.

Table Measurement Result Screen: Main functions of the measurement result screen.

Step 4: Health Diagnosis Screen

The Diagnosis screen displays a composite health assessment that is derived from the measured data. The parameter falling outside the range is determined, and the result displayed in a simple message like BMI Overweight along with a chart showing the status of the measurement readings.

This display serves well when using the kiosk in a college clinic environment since the screen provides information on what the measurements mean to the individual. Rather than displaying numeric data, it shows a key finding and gives remarks based on the parameters.

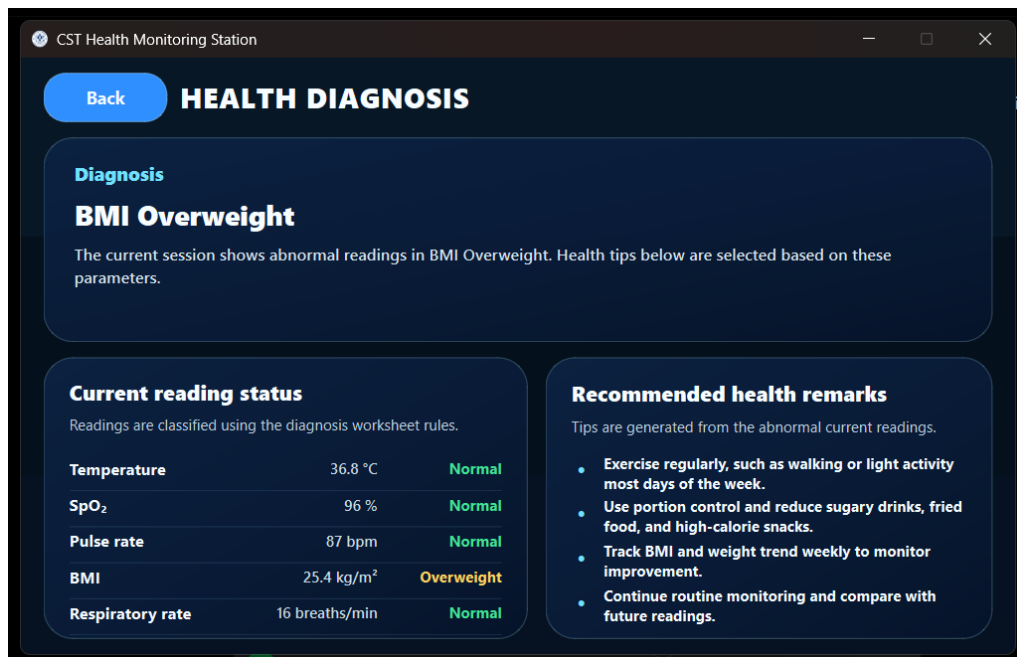


Figure Health Diagnosis Screen: Health Diagnosis Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the previous screen.
Current reading status	Shows each parameter and its classification.
Recommended health remarks	Lists tips linked to the identified condition.

Table Health Diagnosis Screen: Main functions of the health diagnosis screen.

Step 5: Consultation Guidance Screen

The Consultation Guidance screen will provide advisory information for the subsequent step. In other words, this screen will contain key information about the results that were obtained by the user as well as suggestions on what needs to be done next. It will contain health improving advice as well.

This screen will function as a supportive element for the user but not as an alternative to medical consultation. The screen will stimulate the user to make decisions related to responsible behavior.

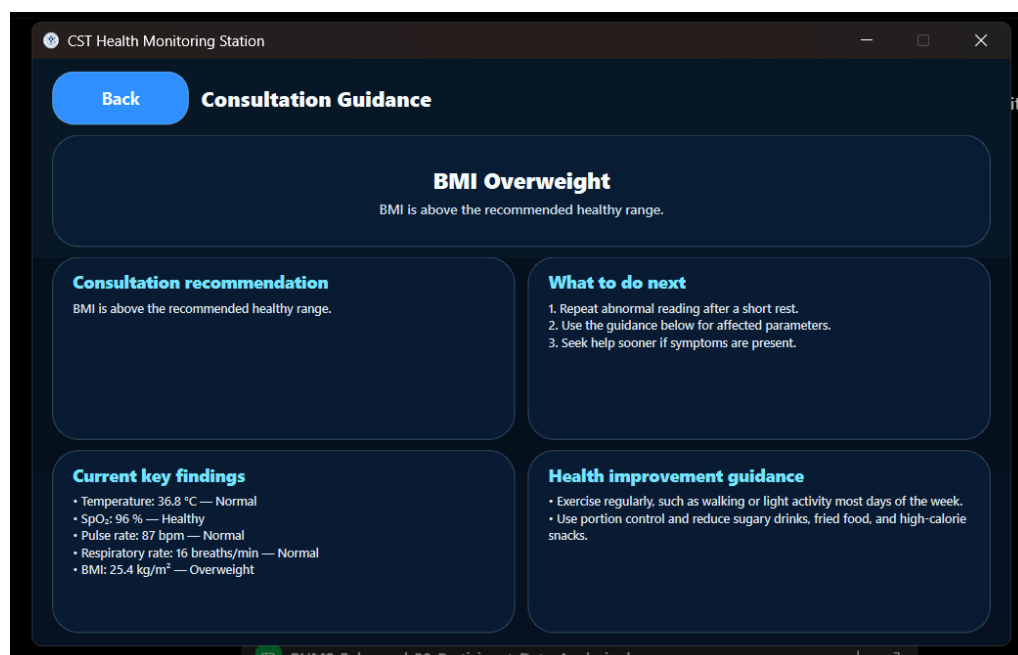


Figure Consultation Guidance Screen: Consultation Guidance Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the previous screen.
Consultation recommendation	Summarizes the consultation advice.
What to do next	Lists the suggested next actions.
Health improvement guidance	Provides simple preventive recommendations.

Table Consultation Guidance Screen: Main functions of the consultation guidance screen.

Interpretation Screens for Specific Parameters

After seeing the results display, the user may tap on any parameter card to see its more detailed interpretation screen. Such screens will provide information about threshold values, the current condition of the parameter, and advice for a specific parameter only.

Step 6: Temperature Detail Display Screen

In the temperature detail screen, information about the current body temperature, the healthy range, the alert value, and the final interpretation will be visible. This screen could prove helpful in determining the legitimacy of the contactless temperature reading.

This screen also tells the user that the temperature should be considered together with possible symptoms and the environment. The implemented version will make it clear that there is no need to do anything if everything is fine, except the alert value for a fever.

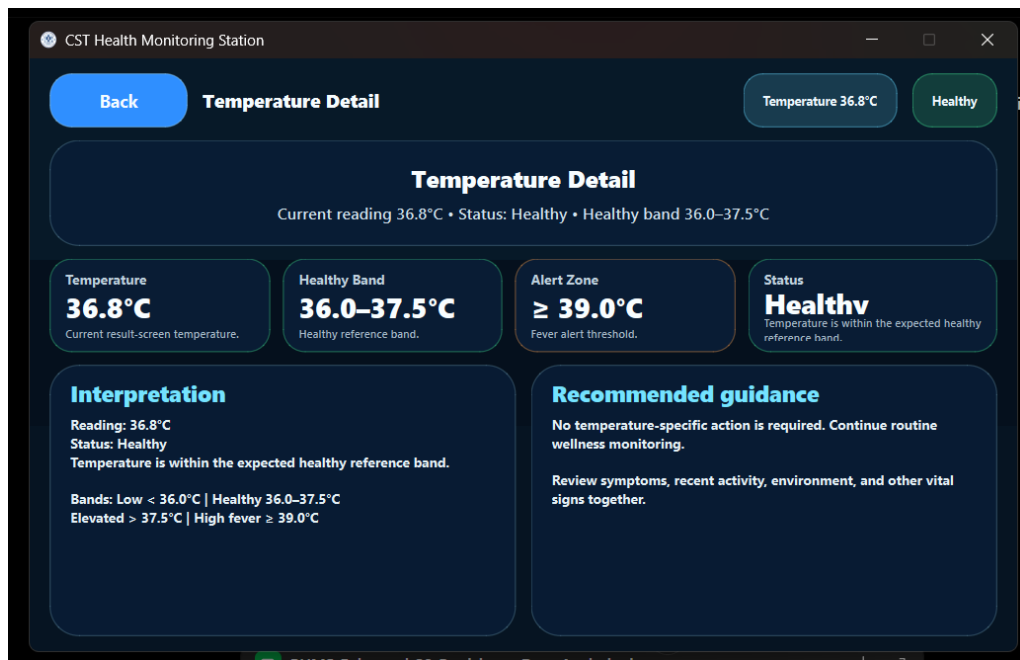


Figure Temperature Detail Screen: Temperature Detail Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the result screen.
Temperature	Shows the measured temperature value.
Healthy band	Shows the normal temperature range.
Alert zone	Shows the threshold for fever alert.
Recommended guidance	Displays temperature-related interpretation.

Table Temperature Detail Screen: Main functions of the temperature detail screen.

Step 7: Blood Oxygen Detail View

The blood oxygen detail view displays the SpO2 reading along with the normal range of readings and low threshold of SpO2 levels. This enables the user to understand the saturation value without having to read the entire report.

This display further informs the user about any possible interference factors like finger movement, poor sensor contact, nail polish, or signal quality issues.

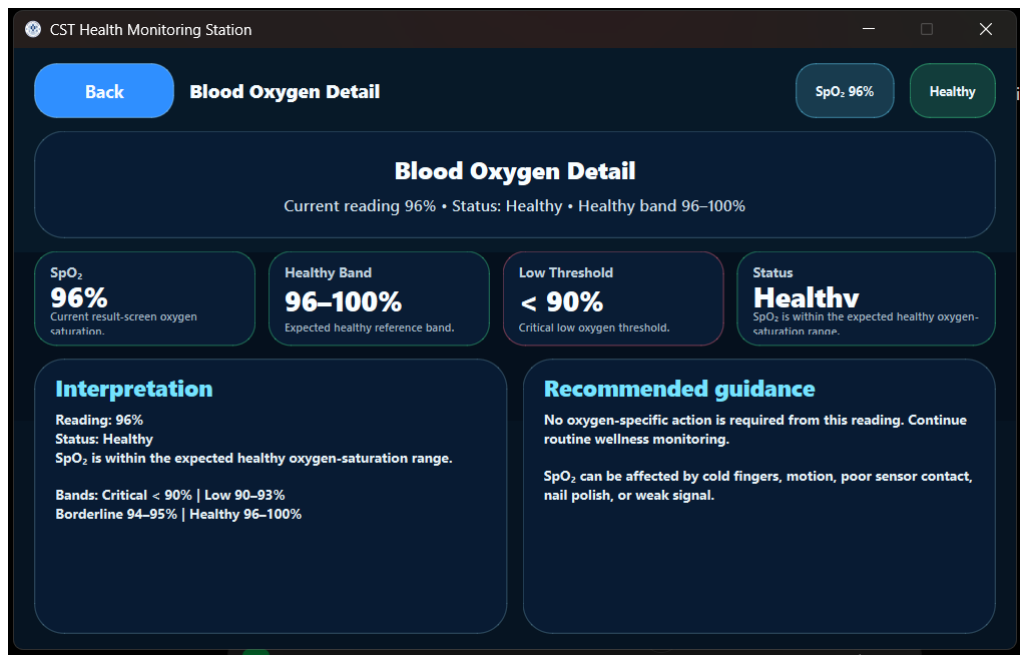


Figure Blood Oxygen Detail Screen: Blood Oxygen Detail Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the result screen.
SpO ₂ value	Shows the measured oxygen saturation.
Healthy band	Displays the expected saturation range.
Low threshold	Shows the critical lower limit.
Recommended guidance	Provides oxygen-saturation interpretation.

Table Blood Oxygen Detail Screen: Main functions of the blood oxygen detail screen.

Step 8: Pulse Rate Detail Screen

This screen highlights the pulse rate figure, comparing it against the average resting rate of adults. The other information highlighted in this screen is the highest level of alertness for abnormal pulse rates.

This screen is necessary since pulse rate may vary depending on various factors like physical exercise, caffeine intake, anxiety, fever, and pain.

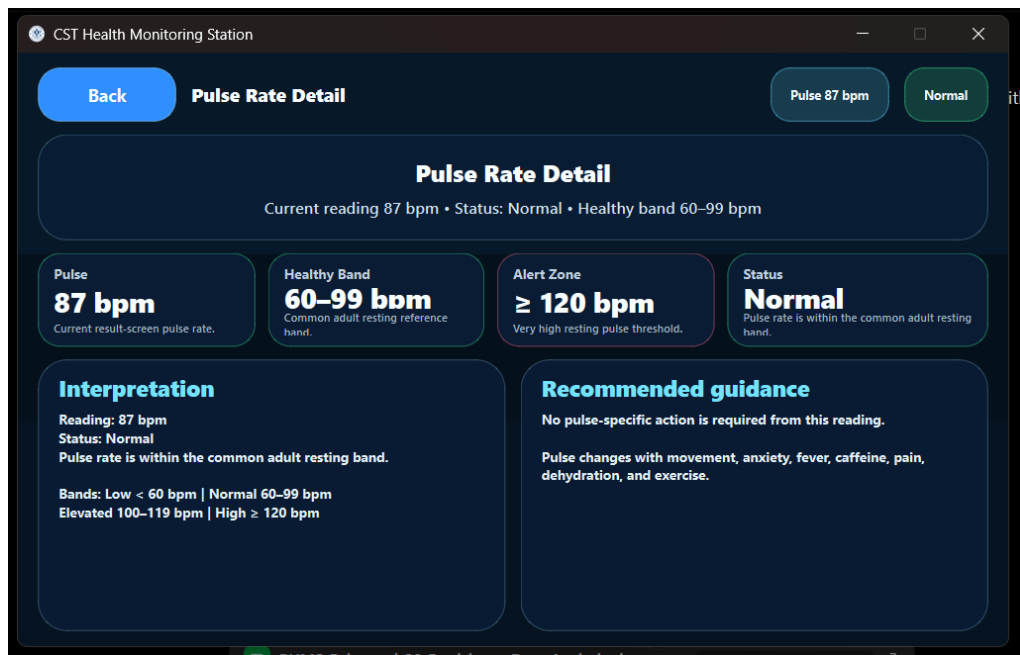


Figure Pulse Rate Detail Screen: Pulse Rate Detail Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the result screen.
Pulse	Shows the measured bpm.
Healthy band	Shows the common adult resting range.
Alert zone	Highlights the high threshold.
Recommended guidance	Provides pulse-related remarks.

Table Pulse Rate Detail Screen: Main functions of the pulse rate detail screen.

Step 9: Detail Screen for BMI

The detail screen for BMI contains both the value computed along with the range for normal, overweight, and obesity levels. The detail screen is critical since BMI is a calculated value based on physical measurements and not a directly measured parameter.

In the actual BHMS software, the detail screen includes explanation regarding the class and actions that should be taken to ensure that the BMI is within the normal range.

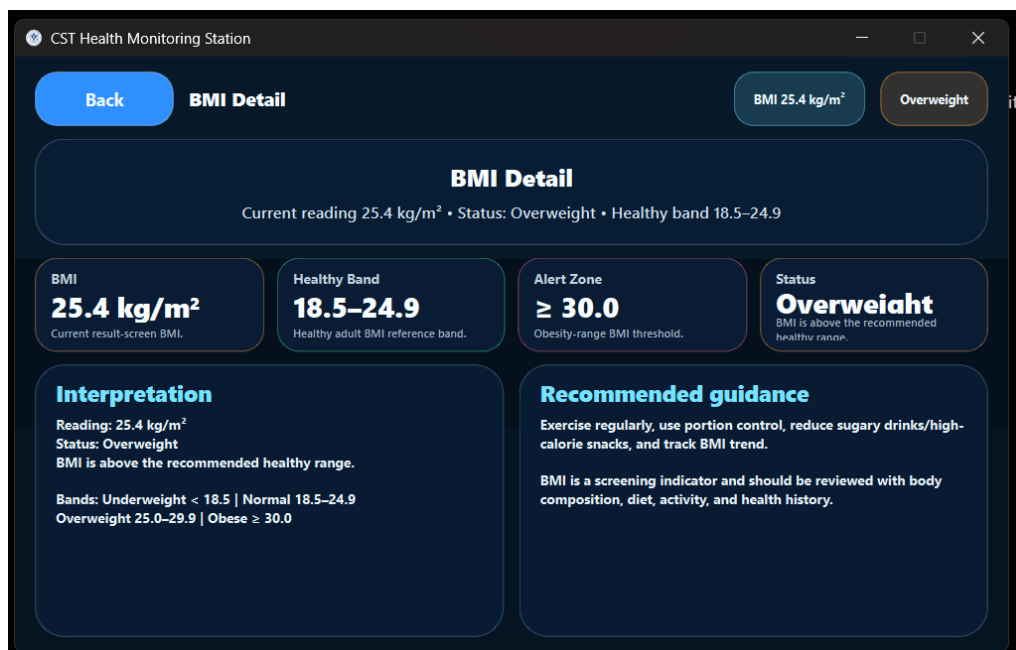


Figure BMI Detail Screen: BMI Detail Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the result screen.
BMI value	Shows the calculated BMI.
Healthy band	Shows the normal BMI interval.
Alert zone	Shows the obesity threshold.
Recommended guidance	Provides BMI-related recommendations.

Table BMI Detail Screen: Main functions of the bmi detail screen.

Step 10: PDF Report and QR-enabled Session Viewer

Following data collection and analysis, the BHMS creates a PDF report and links it with the QR code for session viewing. In this way, the BHMS expands its role from a mere local measuring instrument to a report generator that provides access to the report through either the screen itself or a scanning option using Wi-Fi connectivity to get session information on any mobile device.

The report is comprised of session information, all measurements, diagnosis summary, combination of diagnosis, recommendation, and advice in case of consultation/emergency. On the other hand, the mobile viewer offers similar information in a condensed form.



Health Measurement Report
CST Health Monitoring Station

Overall Status	All parameters are normal	Health Index: 95
Summary	All key readings are within the normal kiosk range for this session.	Excellent

Session Information

Session ID	CST-20280531-96C9FFBA
Measurement Mode	hardware
Overall Status	All parameters are normal
Date	31 May 2026
Time	01:28 AM

Measured Health Parameters

Metric	Value	Category / Interpretation
Weight	56.1 kg	-
Height	167.2 cm	-
BMI	20.1	Normal
Temperature	36.7 °C	Normal
SpO ₂	98 %	Normal
Pulse Rate	90 bpm	Normal
Respiratory Rate	15.0 breaths/min	Normal

Diagnosis Summary
All key readings are within the normal kiosk range for this session.
No major issues were flagged in the current session.

Combined Diagnosis Impression
Likely health remark: No clear issue seen
Temperature is Normal, SpO₂ is Normal, Pulse rate is Normal, BMI is Healthy, Respiration rate is Normal.
All key readings are within the normal kiosk range for this session.
• Continue your normal routine and healthy habits.

CST Health Monitoring Station • 1.0.0 Generated: 2026-05-31T01:28:59

Figure 6.B: Example of the generated BHMS PDF health measurement report.

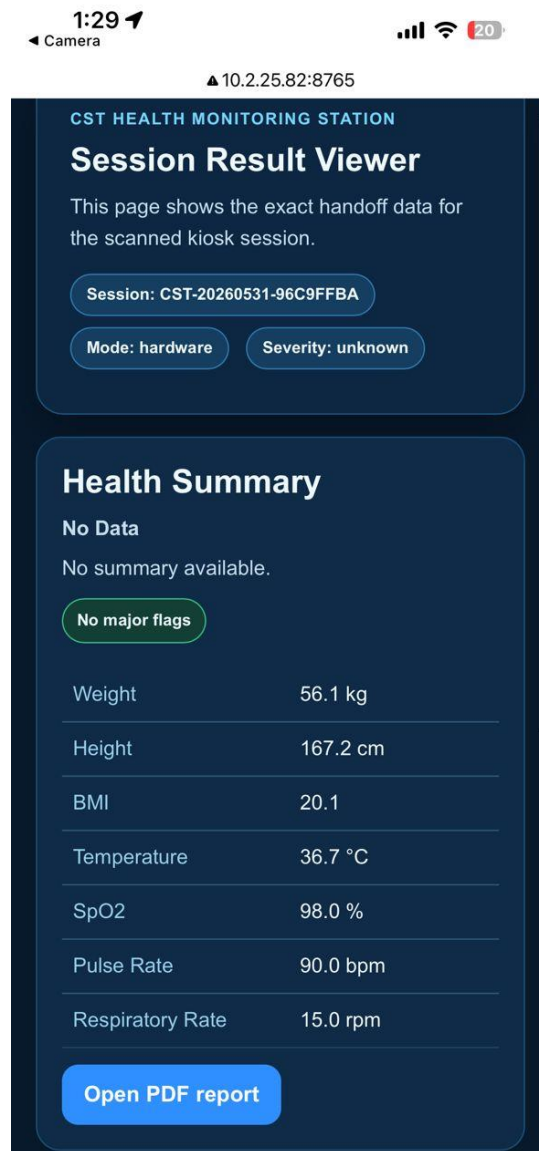


Figure 6.C: Mobile QR session result viewer showing the handoff data.

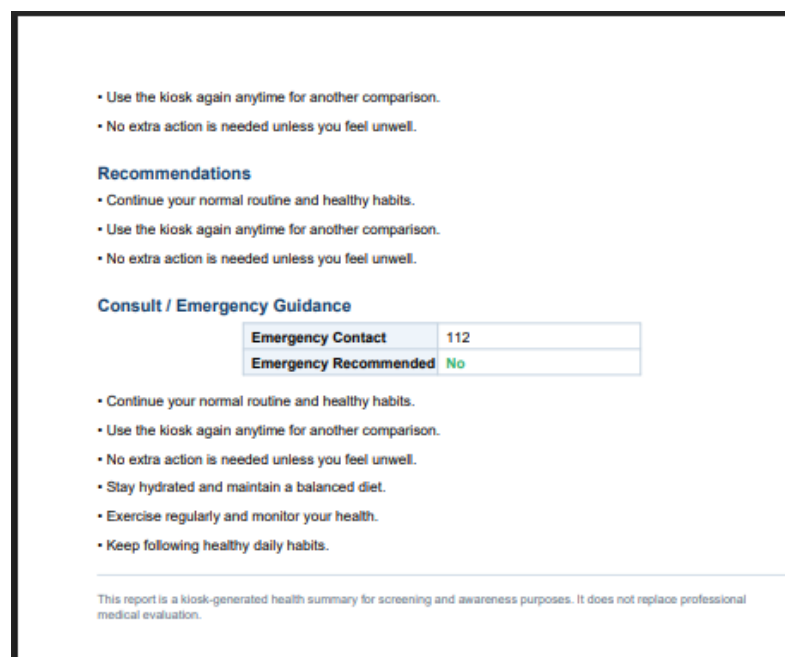


Figure 6.D: Recommendation and emergency guidance section in the generated report.

Report component	Purpose
Session Information	Identifies the measurement session, date, time and mode.
Measured Health Parameters	Lists the final values and interpretations.
Diagnosis / Recommendations	Summarizes the finding and suggests next steps.
QR Viewer	Allows the user to view the report on a phone browser.

Table 6.B: Main components of the BHMS report and QR viewer.

Administrative and Service Support Screens

These screens have not been integrated in the normal measurement process flow, but these are important screens for maintaining and managing the system. They have been shown in the final presentation as well and are added here to fulfill the requirement of complete GUI documentation.

Step 11: Admin Login Screen

The admin login screen ensures that no unauthorized person can use the configuration and maintenance options within the system. It asks the username and password along with providing an option for remembering it and toggling off and on the visibility of the password.

It is required because the configuration and calibration options need to be managed by authorized persons only, and thus the admin login process is separated from the normal student user login process.

Figure Admin Login Screen: Admin Login Screen used in the final BHMS application.

Button / Area	Functionality
Back	Returns to the previous screen.
Username / Password	Accepts administrator credentials.
Show Password	Temporarily reveals the password.
Login	Opens the administrative dashboard.
Clear	Clears the entered credentials.

Table Admin Login Screen: Main functions of the admin login screen.

Step 12: Admin Panel

Upon successful login, the admin panel becomes the central administration hub. The module buttons that will be visible include those of the Settings, Calibration, Storage and Publish modules, along with refresh and logout buttons.

The importance of this screen lies in its role as a one-stop administrative center. Rather than hiding the administration activities in different places, it brings all these functions together on one screen.

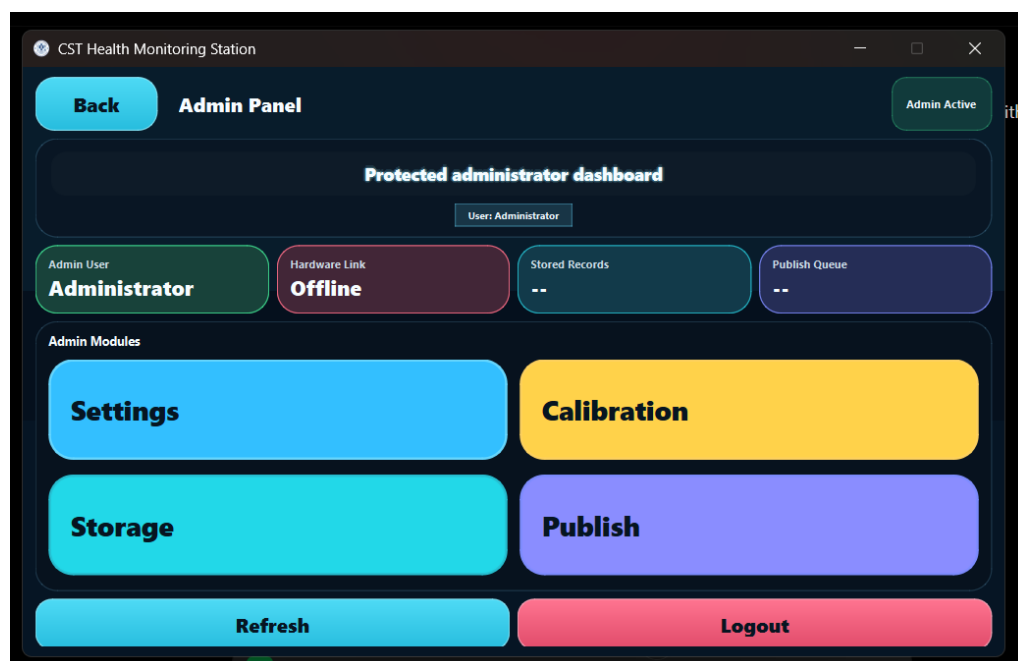


Figure Admin Panel: Admin Panel used in the final BHMS application.

Button / Area	Functionality
Settings	Opens system configuration controls.
Calibration	Opens calibration-related management functions.
Storage	Opens storage maintenance page.
Publish	Reserved for publishing or export workflow.
Refresh	Refreshes the admin dashboard.
Logout	Closes the administrative session.

Table Admin Panel: Main functions of the admin panel.

Step 13: Settings Window

Settings allow you to adjust the following parameters: screen brightness, volume levels, start-up mode, network mode, and touch-audio. In addition, there are the options "Save", "Reload", and "Reset". From the point of view of system design, this window improves usability since it enables the kiosk to adapt quickly to the new environment. For example, the screen brightness can be set according to the clinic lighting, while audio can be turned off in an infirmary environment.

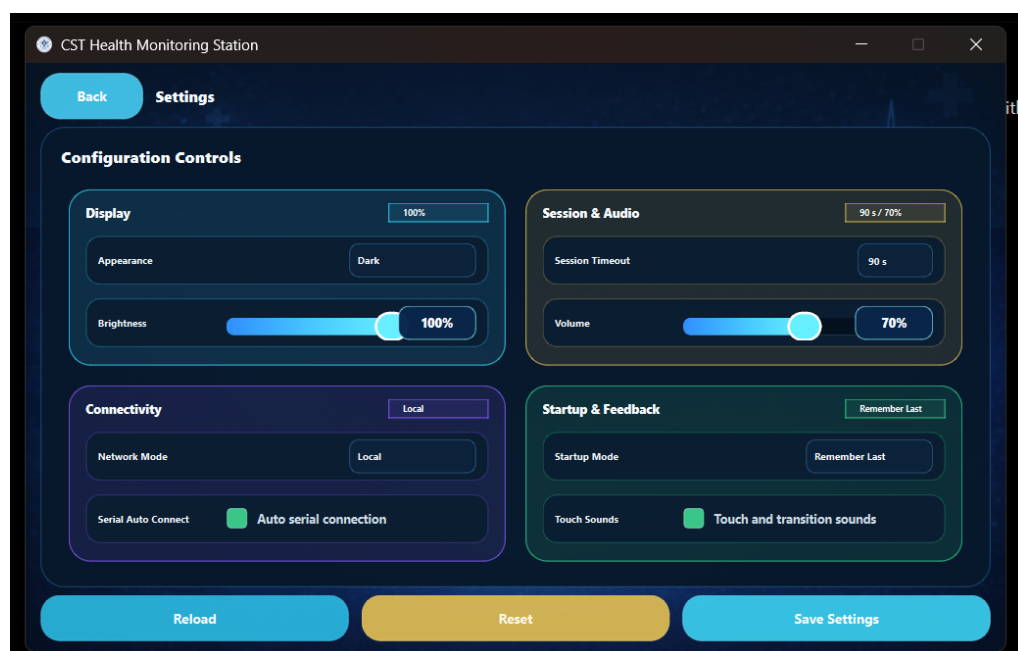


Figure Settings Screen: Settings Screen used in the final BHMS application.

Button / Area	Functionality
Brightness	Adjusts the screen brightness.
Volume	Adjusts the voice and sound level.
Network mode	Displays or sets the network setting.
Save Settings	Saves the updated settings.
Reset	Restores default values.
Reload	Reloads the current settings.

Table Settings Screen: Main functions of the settings screen.

Step 14: Storage Maintenance Screen

The storage maintenance screen will contain information about stored files, their size estimates, storage condition and state of the database. It will also categorize data by records, reports, QR codes, backup, export files, temporary files, and logs.

The screen is necessary in order to run the kiosk on a permanent basis since BHMS keeps records of session outcomes, PDF reports, and other information associated with QR codes.

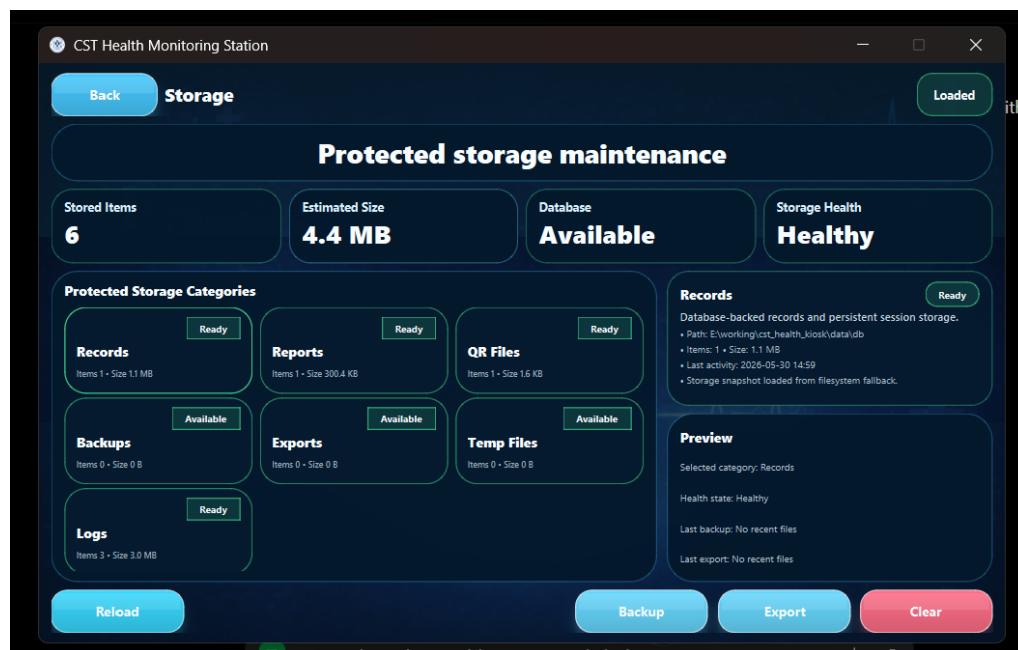


Figure Storage Maintenance Screen: Storage Maintenance Screen used in the final BHMS application.

Button / Area	Functionality
Reload	Reloads the storage summary.
Backup	Backs up stored records.

Export	Exports selected data.
Clear	Clears selected storage items.
Preview area	Shows information about selected category.

Table Storage Maintenance Screen: Main functions of the storage maintenance screen.

Screen 15: Saved Reports Screen

The screen displays the list of report file names that have been generated before, and it provides an option to open the selected report, reload the list of report files, or display the path of the folder where all the report files are stored.

In case of an infirmary, this screen becomes helpful, as the person who operates the system may need to open the generated report file again in PDF format.

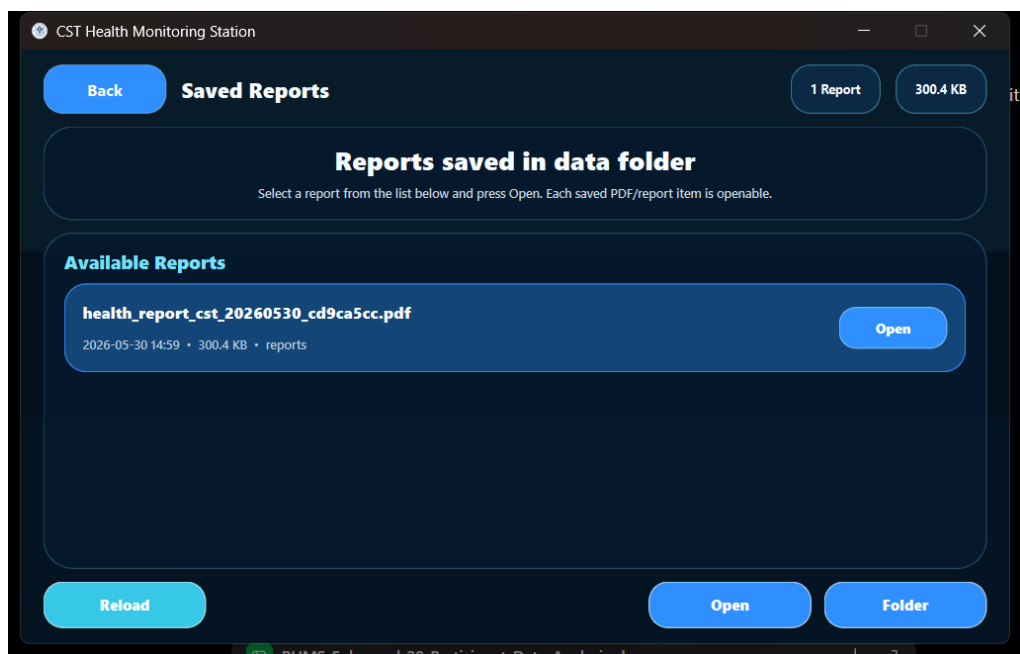


Figure Saved Reports Screen: Saved Reports Screen used in the final BHMS application.

Button / Area	Functionality
Open	Opens the selected report file.
Reload	Refreshes the report list.
Folder	Opens the reports folder location.
Available reports list	Shows all saved report files.

Table Saved Reports Screen: Main functions of the saved reports screen.

As seen in the above-mentioned sequence, the actual steps carried out by the user precede those performed for services and administration purposes. The following subsections in this chapter present the detailed screen-by-screen documentation and other software modules for BHMS.

6.5.1 Saved Reports Screen

Open, Reload, Folder buttons for report access and storage navigation.

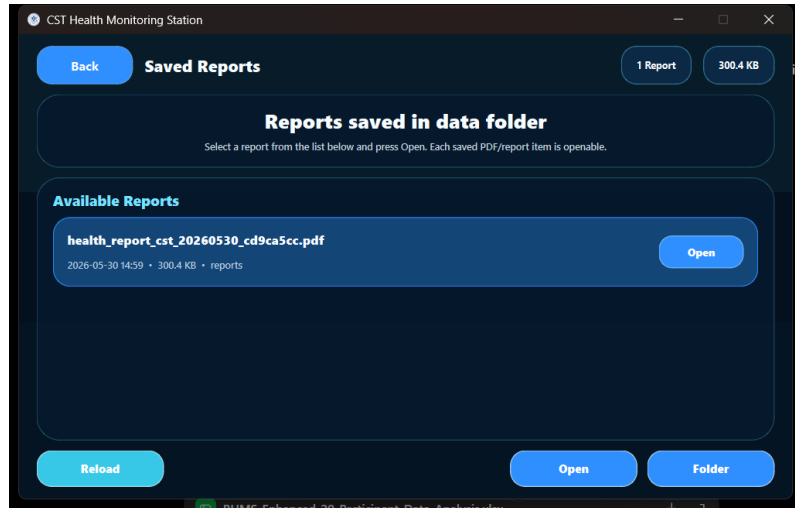


Figure 6.1: Saved Reports Screen

Table 6.1: Button and function description for Saved Reports Screen

Button / Indicator	Function
Reload	Refreshes file list or storage status.
Open	Opens selected saved report.
Folder	Opens storage directory.
Backup	Creates backup of records.
Export	Exports available data.
Clear	Deletes selected temporary or protected files after confirmation.

6.5.2 Storage Maintenance Screen

Backup, Export, Clear and Reload controls for protected storage maintenance.

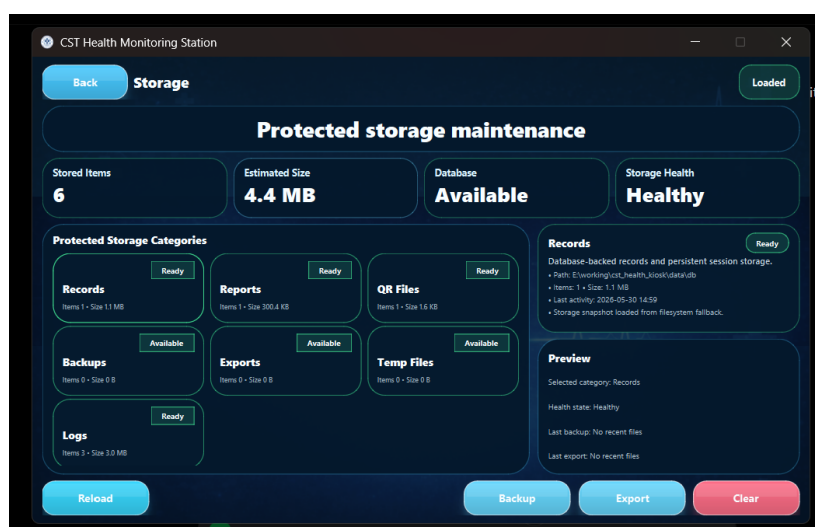


Figure 6.2: Storage Maintenance Screen

Table 6.2: Button and function description for Storage Maintenance Screen

Button / Indicator	Function
Reload	Refreshes file list or storage status.
Open	Opens selected saved report.
Folder	Opens storage directory.
Backup	Creates backup of records.
Export	Exports available data.
Clear	Deletes selected temporary or protected files after confirmation.

6.5.3 Settings Screen

Brightness, volume, appearance, network mode, touch sounds and save settings controls.

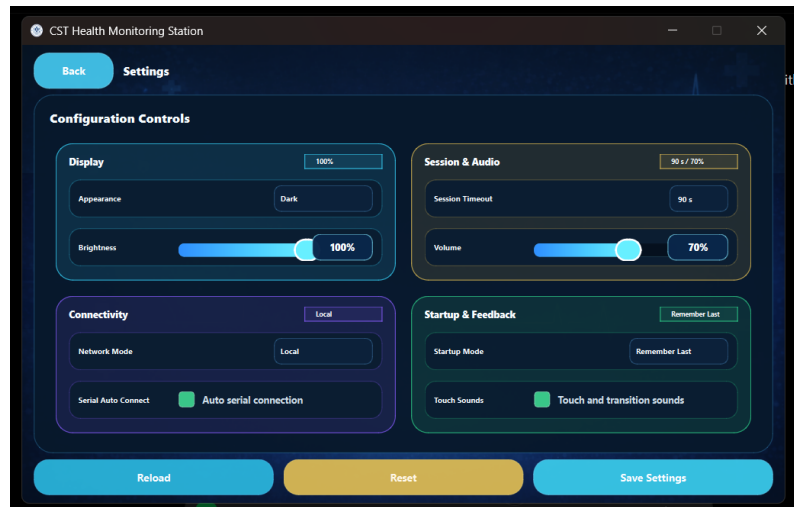


Figure 6.3: Settings Screen

Table 6.3: Button and function description for Settings Screen

Button / Indicator	Function
Appearance	Controls dark theme setting.
Brightness	Adjusts screen brightness.
Volume	Adjusts audio prompt volume.
Serial Auto Connect	Enables automatic hardware communication.
Save Settings	Stores the selected configuration.
Reset	Restores default settings.

6.5.4 Admin Panel

Admin modules for settings, calibration, storage, publish, refresh and logout.

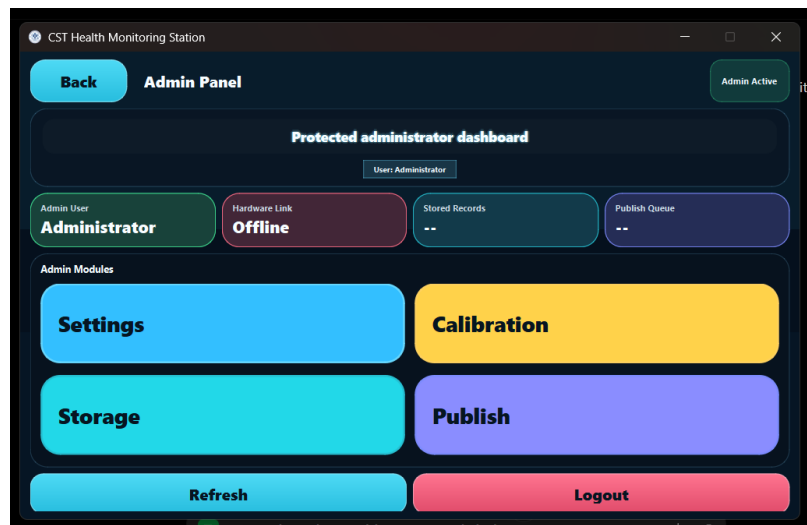


Figure 6.4: Admin Panel

Table 6.4: Button and function description for Admin Panel

Button / Indicator	Function
Settings	Opens configuration controls.
Calibration	Opens calibration-related settings.
Storage	Opens storage maintenance screen.
Publish	Manages report publishing queue.
Refresh	Reloads admin values.
Logout	Exits protected admin mode.

6.5.5 Admin Login Screen

Protected admin login with username, password, remember me, show password, clear and login controls.

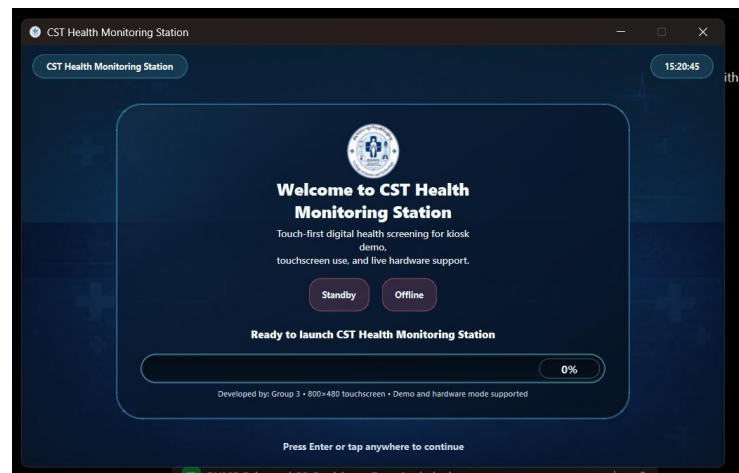


Figure 6.5: Admin Login Screen

Table 6.5: Button and function description for Admin Login Screen

Button / Indicator	Function
Back	Returns to previous screen.
Login	Authenticates admin user.
Clear	Clears entered credentials.
Remember Me	Stores user preference for login.
Show Password	Toggles password visibility.

6.5.6 Instruction and Mode Selection Screen

Start Check Up, Admin and Hardware Mode buttons with operating instructions.

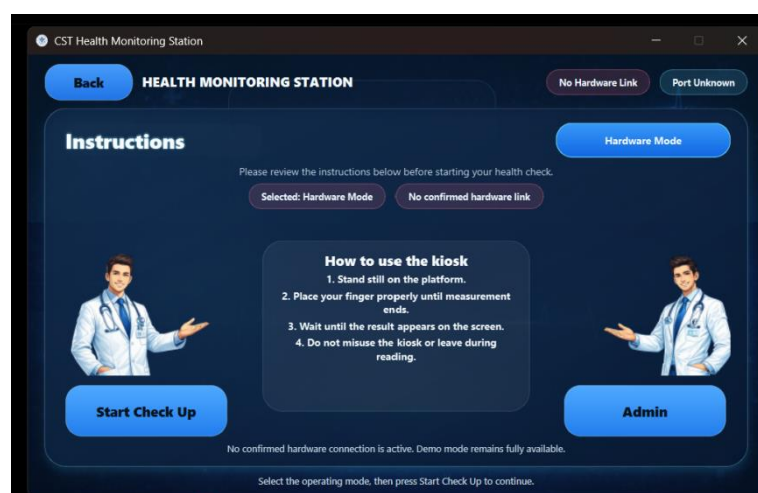


Figure 6.6: Instruction and Mode Selection Screen

Table 6.6: Button and function description for Instruction and Mode Selection Screen

Button / Indicator	Function
Start Check Up	Begins the health measurement workflow.
Admin	Opens protected admin login panel.
Hardware Mode	Selects physical sensor-based measurement mode.
Back	Returns to the previous screen.

6.5.7 Welcome Screen

Welcome interface with loading bar and standby/offline indicators.

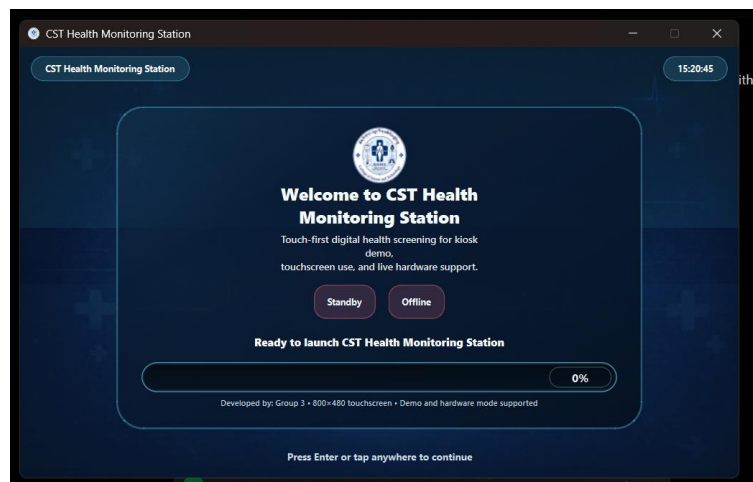


Figure 6.7: Welcome Screen

Table 6.7: Button and function description for Welcome Screen

Button / Indicator	Function
Standby / Offline	Shows current launch state before starting the kiosk application.
Loading Bar	Indicates application loading status.
Tap/Enter Instruction	Allows user to proceed to the next screen.

6.5.8 Measurement Result Screen

Health metrics include seven indicators and actions like report, QR code, diagnose, consult, refresh and retake.

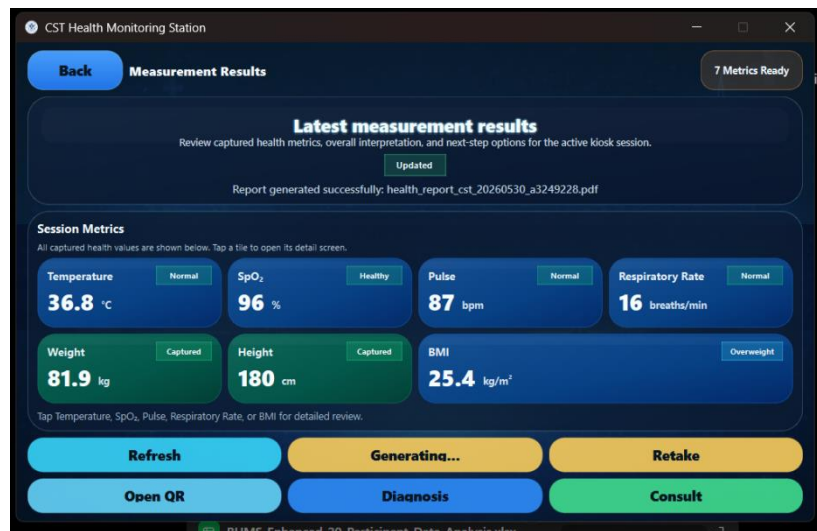


Figure 6.8: Measurement Result Screen

Table 6.8: Button and function description for Measurement Result Screen

Button / Indicator	Function
Refresh	Updates the displayed results.
Report	Generates or opens the session report.
Retake	Restarts the measurement process.
Open QR	Displays the QR page for report access.
Diagnosis	Opens health diagnosis screen.
Consult	Shows consultation guidance based on current readings.

6.5.9 Health Diagnosis Screen

Shows abnormal reading diagnosis and parameter-based health remarks.



Figure 6.9: Health Diagnosis Screen

Table 6.9: Button and function description for Health Diagnosis Screen

Button / Indicator	Function
Back	Returns to the result screen.
Current Reading Status	Lists the status of each measured parameter.
Recommended Health Remarks	Displays parameter-based tips for abnormal readings.

6.5.10 Consultation Guidance Screen

Shows current key findings, recommendation and next action instructions.

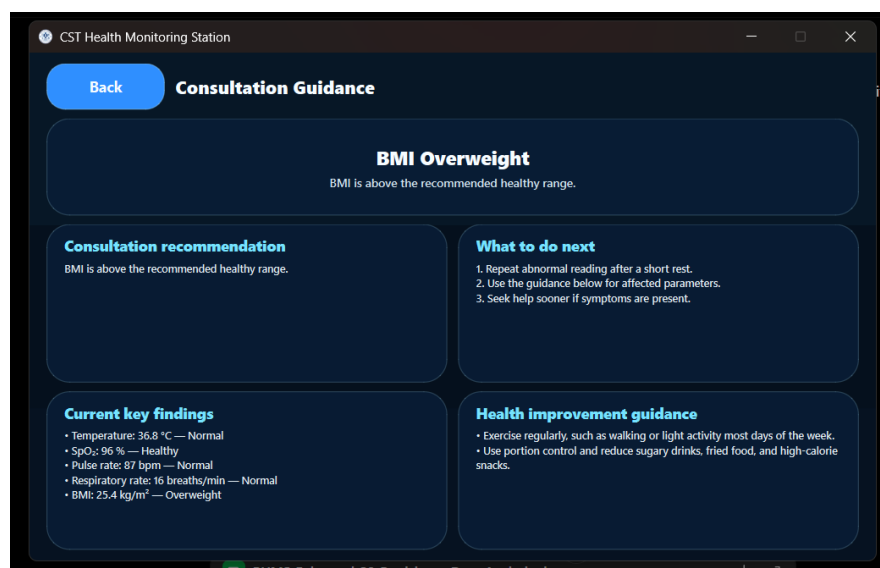


Figure 6.10: Consultation Guidance Screen

Table 6.10: Button and function description for Consultation Guidance Screen

Button / Indicator	Function
Back	Returns to diagnosis/result interface.
Consultation Recommendation	Displays main health concern.
What To Do Next	Provides next action guidance.
Health Improvement Guidance	Provides basic preventive health suggestions.

6.5.11 Temperature Detail Screen

Provides temperature value, healthy range, alert zone, status and guidance.

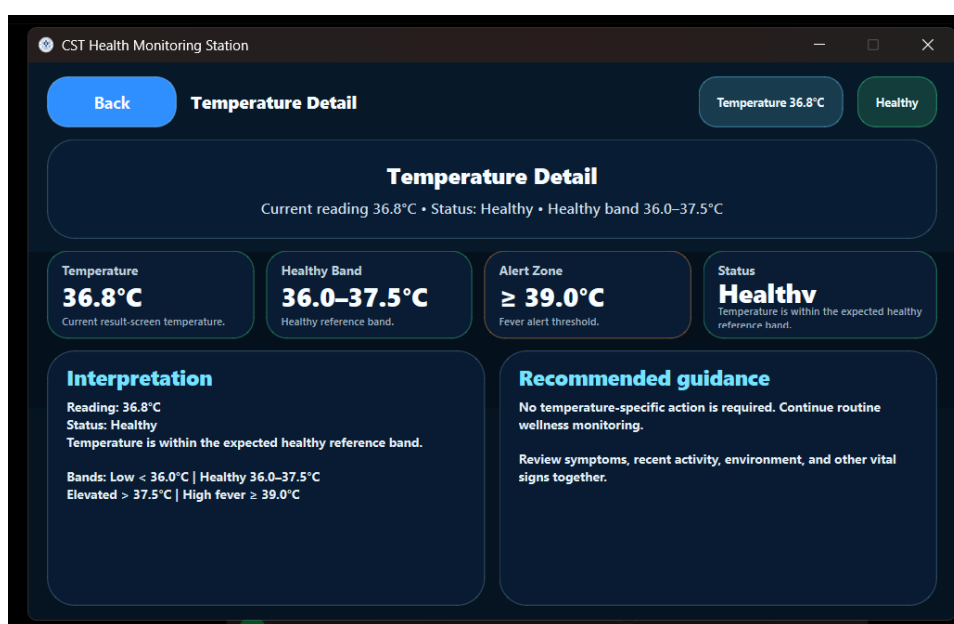


Figure 6.11: Temperature Detail Screen

Table 6.11: Button and function description for Temperature Detail Screen

Button / Indicator	Function
Back	Returns to the result screen.
Healthy Band	Shows the normal reference range.
Alert Zone	Shows threshold for abnormal classification.
Status	Shows classification result.
Recommended Guidance	Explains the meaning of the reading and user action.

6.5.12 Blood Oxygen Detail Screen

Provides SpO₂ value, healthy band, low threshold, status and guidance.

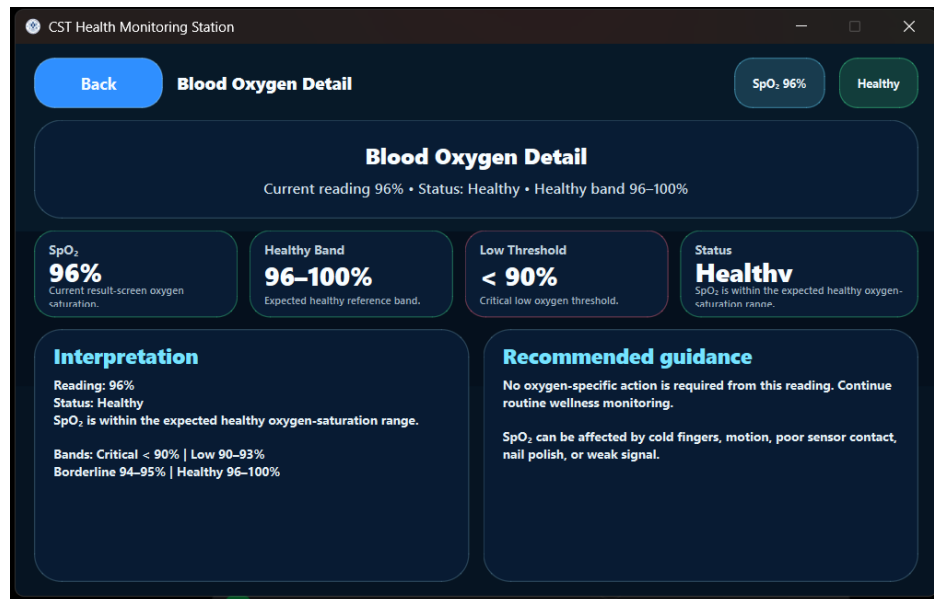


Figure 6.12: Blood Oxygen Detail Screen

Table 6.12: Button and function description for Blood Oxygen Detail Screen

Button / Indicator	Function
Back	Returns to the result screen.
Healthy Band	Shows the normal reference range.
Alert Zone	Shows threshold for abnormal classification.
Status	Shows classification result.
Recommended Guidance	Explains the meaning of the reading and user action.

6.5.13 Pulse Rate Detail Screen

Provides pulse value, normal band, alert zone, status and guidance.

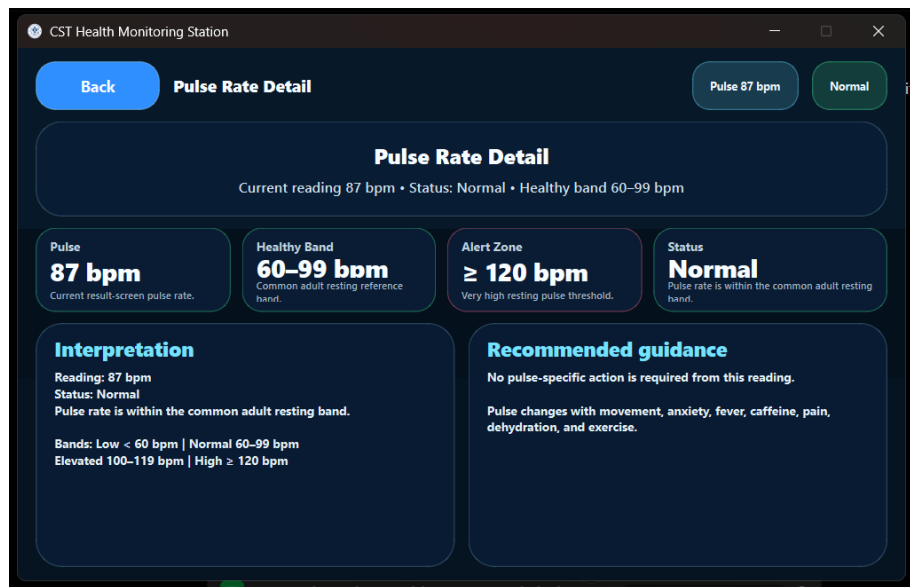


Figure 6.13: Pulse Rate Detail Screen

Table 6.13: Button and function description for Pulse Rate Detail Screen

Button / Indicator	Function
Back	Returns to the result screen.
Healthy Band	Shows the normal reference range.
Alert Zone	Shows threshold for abnormal classification.
Status	Shows classification result.
Recommended Guidance	Explains the meaning of the reading and user action.

6.5.14 BMI Detail Screen

Provides BMI value, healthy band, overweight/obese alert zones and guidance.

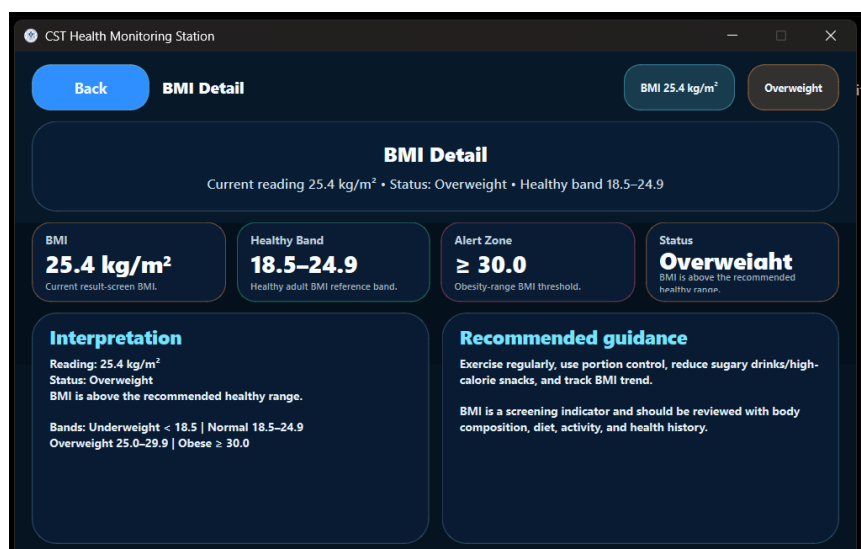


Figure 6.14: BMI Detail Screen

Table 6.14: Button and function description for BMI Detail Screen

Button / Indicator	Function
Back	Returns to the result screen.
Healthy Band	Shows the normal reference range.
Alert Zone	Shows threshold for abnormal classification.
Status	Shows classification result.
Recommended Guidance	Explains the meaning of the reading and user action.

6.6 Voice Guidance Module

The BHMS system also has an audio guidance module which makes use of DFPlayer Mini voice module along with a 3-watt speaker. Commands for triggering voice messages through the serial commands sent by ESP32 occur at every measurement step. The order of messages is as follows: (1) "welcome message"; (2) "Turn the temperature knob on your head"; (3) "temperature measurement complete, insert left index finger in left slot"; (4) "measurement complete, insert right index finger in right slot"; (5) "BPM measurement complete, please stand still for weighing"; (6) "weight measurement complete, please stand straight"; (7) "height measurement complete, check results on display and scan QR code".

Table 6.17: Voice assistance hardware used in BHMS

Component	Role	Description
Obocraze MP3-TF-16P V3.0 SD Card Module	Serial MP3 playback module	Supports FAT16/FAT32 file system, serial command control, 24-bit DAC output, dynamic range around 90 dB and SNR around 85 dB.
ND TECH 2 Inch 3W 4 Ω Full Range Speaker	Audio output device	Mini woofer loudspeaker used to play voice prompts from the MP3 module.
MicroSD Card	Audio storage	Stores recorded MP3 prompt files in the required sequence.
ESP32 / Raspberry Pi command	Trigger mechanism	Triggers the correct audio file at each

signal

measurement stage.

Table 6.18: Voice prompt sequence used during measurement workflow

Measurement Stage	Voice Prompt	Purpose
Welcome screen	“Welcome to CST Basic Health Monitoring Station.”	Played after the user enters the welcome screen.
Start measurement	“Adjust the temperature screw over your head.”	Guides user to align temperature measurement position.
Temperature completed	“Temperature measurement done. Move your left index finger and place it on the left hole in sensor box.”	Guides user to place finger for SpO ₂ and respiratory-rate measurement.
SpO ₂ /RR completed	“Reading taken. Move your right index finger and place it in right hole.”	Guides user to place finger for pulse/BPM measurement.
BPM completed	“BPM reading is done. Stand still for weight measurement.”	Prepares user for stable weight reading.
Weight completed	“Weight reading is done. Stand straight.”	Guides user to stand correctly for height measurement.
Height completed	“Height measurement is done. Now you can view the result on the screen. Use the Wi-Fi provided to scan the QR to see your report.”	Informs user that final result and QR report are available.

6.7 Serial Communication

Serial communication acts as the link between the ESP32 sensor module and the Raspberry Pi processing module. Information for each parameter is transferred in packets with appropriate labeling to ensure proper identification of height, weight, BMI, temperature, pulse, oxygen saturation levels, and respiratory rate by the receiver.

Table 6.19: Serial data labels used for ESP32 to Raspberry Pi communication

Data Label	Meaning	Example Packet
------------	---------	----------------

HEIGHT	Height value in centimetres	HEIGHT: 164.0
WEIGHT	Weight value in kilograms	WEIGHT: 57.9
BMI	Calculated BMI value	BMI: 21.5
TEMP	Temperature in degree Celsius	TEMP: 36.8
SPO2	Blood oxygen saturation	SPO2: 98
PULSE	Pulse rate in beats per minute	PULSE: 78
RR	Respiratory rate in breaths per minute	RR: 16

There were no instances of significant delay while transferring information from ESP32 to Raspberry Pi during testing. Application worked efficiently after completing measurement, showing obtained results on the screen.

6.8 Data Processing Algorithm

Data processing algorithm is used to convert information from sensors to health indicators. It detects whether information is present and filters unstable information until necessary measurements are taken.

Table 6.20: Data processing algorithm used in BHMS

Step	Processing Action
1	Initialize Raspberry Pi GUI, ESP32 and sensors.
2	Check serial connection and hardware mode status.
3	Read temperature and confirm stable value.
4	Read SpO ₂ and respiratory-rate value from optical sensing stage.
5	Read pulse value from the pulse/BPM stage.
6	Read weight from load cell after stabilization.
7	Read height from VL53L0X distance sensor.
8	Calculate BMI using height and weight.
9	Classify each parameter as low, normal, high or abnormal according to rule set.
10	Display results, generate PDF report and QR code.

6.9 Calculation of Body Mass Index

Body Mass Index can be calculated automatically using the input values for weight and height. This helps minimize errors in calculating the body mass index and also ensures that classification is done instantaneously.

$$\text{BMI} = \text{Weight (kg)} / [\text{Height(m)}]^2$$

For instance, if the weight is 58 kg while the height is 1.64 meters, then the Body Mass Index will be $58 / (1.64 \times 1.64) = 21.6 \text{ kg/m}^2$.

6.10 Health Classification Logic

Table 6.21: Health classification logic used in BHMS

Parameter	Classification Range	Use in System
Temperature	Low < 36.0 °C; Normal 36.0–37.5 °C; Elevated > 37.5 °C; High fever ≥ 39.0 °C	Used for temperature detail screen and diagnosis.
SpO ₂	Critical < 90%; Low 90–93%; Borderline 94–95%; Healthy 96–100%	Used for blood oxygen detail and report.
Pulse Rate	Low < 60 bpm; Normal 60–99 bpm; Elevated 100–119 bpm; High ≥ 120 bpm	Used for pulse detail and guidance.
BMI	Underweight < 18.5; Normal 18.5–24.9; Overweight 25.0–29.9; Obese ≥ 30.0	Used for BMI detail and diagnosis.
Respiratory Rate	Low, normal and high ranges based on adult resting breathing classification	Used for respiratory-rate status and diagnosis.

This classification logic is meant for preliminary testing purposes only. It is not an alternative to a clinical diagnosis. Abnormal test results are meant to suggest repeat testing or consulting a medical expert.

6.11 QR Report Creation

Once the testing process has been done by the user, a PDF report of all the test parameters, including classification and health information, along with date and time, will be created by the Raspberry Pi. A QR code will then be created to allow users to access the report on their mobile devices after connecting to the Wi-Fi network.

The reporting feature will facilitate documentation, as results will not have to be manually documented. This feature also enables the user to keep an electronic copy of the results for future comparison purposes.

The functionality of the report and QR code viewer was tested. The report is composed of session ID, values of parameters, statuses, diagnosis, recommendation, and advice on

emergencies. On the other hand, the mobile QR code viewer provides the same information as the session handover data but with a more suitable layout for phones.

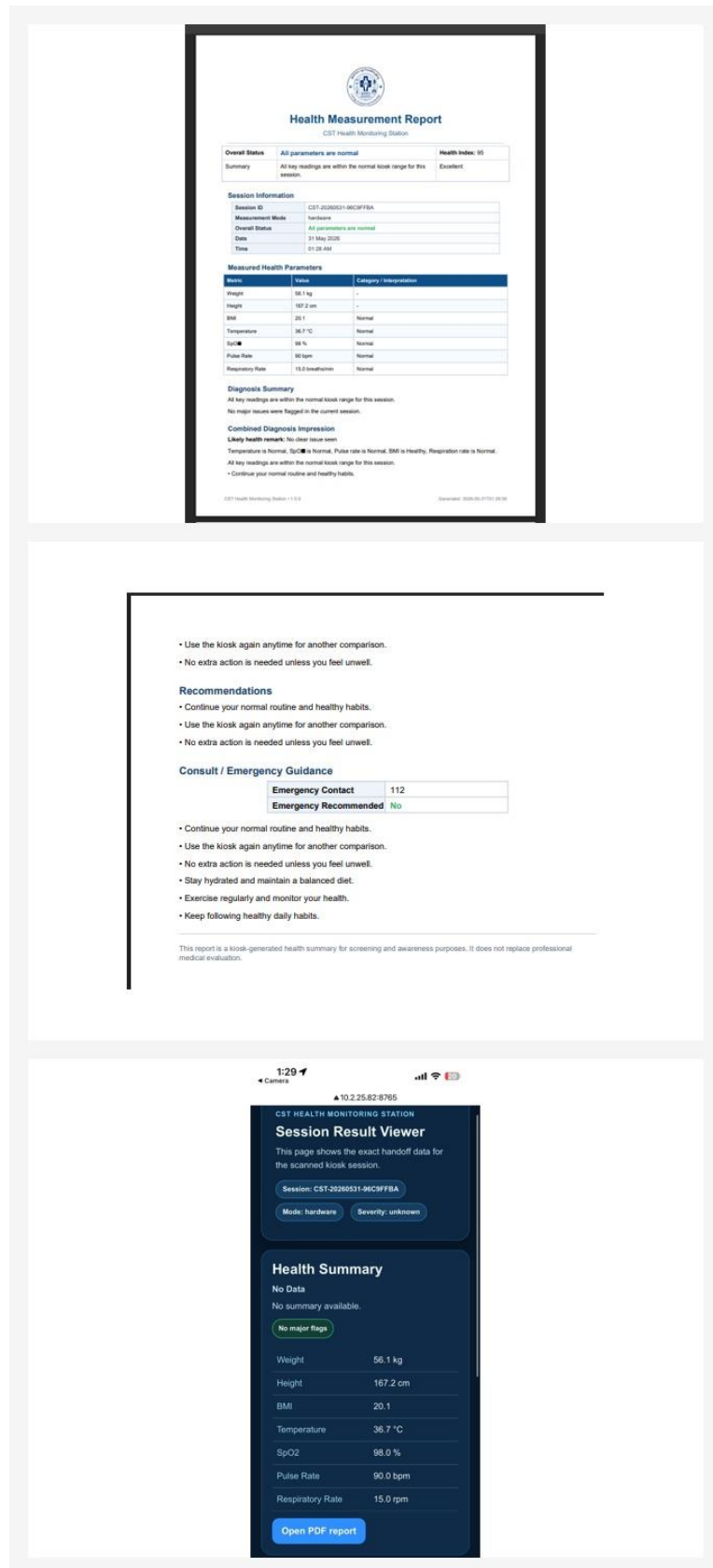


Figure 6.15(a): Generated PDF report pages and the mobile QR-linked session-result viewer produced by the BHMS after measurement completion.

These additional figures have been included from the slides and testing documentation for the project to ensure that the report reflects the actual process followed during the development and testing of the Basic Health Monitoring Station.



Figure 6.15: Generated PDF health measurement report from the BHMS.

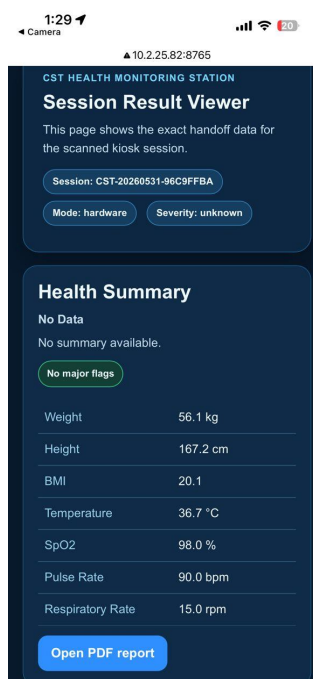


Figure 6.16: Mobile QR session result viewer showing the scanned report data.

• Use the kiosk again anytime for another comparison.
• No extra action is needed unless you feel unwell.

Recommendations

- Continue your normal routine and healthy habits.
- Use the kiosk again anytime for another comparison.
- No extra action is needed unless you feel unwell.

Consult / Emergency Guidance

Emergency Contact	112
Emergency Recommended	No

- Continue your normal routine and healthy habits.
- Use the kiosk again anytime for another comparison.
- No extra action is needed unless you feel unwell.
- Stay hydrated and maintain a balanced diet.
- Exercise regularly and monitor your health.
- Keep following healthy daily habits.

This report is a kiosk-generated health summary for screening and awareness purposes. It does not replace professional medical evaluation.

Figure 6.17: Emergency guidance and recommendation section included in the report.

CHAPTER 7: TESTING WITH NUMBER OF SAMPLES

7.1 Introduction

Tests were performed to confirm individual parts, the functionality of the system as well as the measuring consistency of the Basic Health Monitoring Station. These tests consisted of an individual component test, calibration, integration, and data acquisition based on samples. Finally, the collected data set was increased to 20 individuals to ensure better basis for analysis.

The goal of this chapter is to introduce tests that have been implemented to validate ESP32, Raspberry Pi, sensors, display, voice module, and measurement process. Calibration tables as well as offsets are introduced here to increase measuring accuracy.

7.2 Individual Component Test

Before connecting all parts together, the individual components had to be tested first. It ensured correct connections as well as eliminated any errors such as unstable sensor value readings or communication problems.

Table 7.1: Individual component testing summary

Component	Testing Method	Observation
ESP32 WROOM-32	Serial monitor and sensor reading test	Successfully collected labelled readings from sensors.
Raspberry Pi 4B	Operating system and GUI launch test	Successfully ran the PyQt touchscreen application.
7-inch Touchscreen Display	Touch response and screen layout test	Buttons and screen transitions responded properly.
HX711 + Load Cells	Known weight testing	Weight values were corrected through offset calibration.
MAX30102	Finger placement test	SpO ₂ and pulse-related signals were obtained after calibration.
MLX90614	Reference thermometer comparison	Non-contact temperature readings required offset correction.
VL53L0X	Known distance testing	Distance offset was applied for height calculation.
Voice Module and Speaker	Audio file playback test	Voice prompts played according to the measurement workflow.

7.3 Calibration Process

Calibration was done by comparing sensor readings with reference or real values. Afterward, offsets were determined and applied in order to minimize the gap between sensor and reference values. The following are the processes involved in the calibration of each sensor.

7.3.1 MAX30102 Calibration

The MAX30102 sensor was calibrated by comparing the raw SpO₂ values with reference values from an oximeter. The raw sensor values were below the reference values. As such, a positive offset of +1.67% was used.

Table 7.2: MAX30102 SpO₂ calibration before correction

Reading	Medical Device SpO ₂ (%)	MAX30102 Raw SpO ₂ (%)	Offset (%)
1	98.00	96.31	-1.69

2	97.00	95.42	-1.58
3	99.00	97.26	-1.74
4	96.00	94.28	-1.72
5	98.00	96.37	-1.63
6	97.00	95.29	-1.71
7	99.00	97.41	-1.59
8	98.00	96.25	-1.75
9	97.00	95.36	-1.64
10	96.00	94.35	-1.65

Table 7.3: MAX30102 SpO₂ calibration after correction

Reading	Medical Device SpO ₂ (%)	MAX30102 Raw (%)	Correction Applied (%)	Corrected MAX30102 (%)
1	98	96.46	+1.67	98.13
2	97	95.21	+1.67	96.88
3	99	97.34	+1.67	99.01
4	96	94.41	+1.67	96.08
5	98	96.18	+1.67	97.85
6	97	95.47	+1.67	97.14
7	99	97.12	+1.67	98.79
8	98	96.52	+1.67	98.19
9	97	95.33	+1.67	97.00
10	96	94.26	+1.67	95.93

MAX30102 calibration

The application of the correction factor made the MAX30102 sensor readings close to the reference values and acceptable for screening at a prototype level.

7.3.2 MLX90614 Calibration

The MLX90614 infrared sensor was calibrated using reference temperature values. Raw data from the MLX90614 infrared sensor was slightly lower than the reference value. Correction factor of +2.30 °C was added.

Table 7.4: MLX90614 temperature calibration

Reading	Actual Temperature (°C)	Sensor Raw Temperature (°C)	Offset Applied (°C)	Corrected Temperature (°C)
1	36.30	34.02	+2.30	36.32
2	36.10	33.79	+2.30	36.09

3	35.90	33.61	+2.30	35.91
4	36.40	34.08	+2.30	36.38
5	36.20	33.95	+2.30	36.25

7.3.3 Calibration of Pulse Sensor

The pulse sensors were compared with a BPM measuring tool. There is an adjustment factor of -0.30 BPM which can be applied to improve the performance of the pulse sensors. The accuracy of pulse detection depends on finger movement and pressure.

Table 7.5: Pulse rate calibration

Reading	Actual Pulse Rate (BPM)	Sensor Raw Pulse Rate (BPM)	Offset Applied (BPM)	Corrected Pulse Rate (BPM)
1	77.00	78.40	-0.30	78.10
2	78.00	76.90	-0.30	76.60
3	76.00	76.80	-0.30	76.50
4	79.00	79.60	-0.30	79.30
5	75.00	75.20	-0.30	74.90

7.3.4 VL53L0X Sensor Calibration

Calibration of the VL53L0X sensor was done through comparison between the measured distance and the actual distances. The offset used was around -40 mm. Height is measured using the fixed position of the sensor and then subtracted from the measured distance, so positioning matters.

Table 7.6: VL53L0X height sensor calibration

Reading	Sensor Raw Distance (mm)	Calibrated Distance (mm)	Offset Applied (mm)	Calculated Person Height (cm)
1	404.00	363.55	-40.45	163.60
2	405.00	364.45	-40.55	163.60
3	405.60	364.99	-40.61	163.50
4	392.60	353.30	-39.30	164.70
5	395.80	356.18	-39.62	164.40
6	400.40	360.31	-40.09	164.00

7.3.5 Calibration of Load Cell

Load cell calibration involved adding different standard loads ranging from 0 kg to 200 kg in steps of 20 kg. The reading was taken, and then any offsets that might exist were accounted for. For example, at 200 kg the raw reading was 196.30 kg, and an offset of +3.70 kg was needed.

Table 7.7: Load cell calibration using standard weights

Reading No.	Standard Weight Applied (kg)	Raw Load Cell Reading (kg)	Offset / Correction (kg)	Corrected Reading (kg)
0	0.00	0.00	0.00	0.00
1	20.00	19.62	0.38	20.00
2	40.00	39.28	0.72	40.00
3	60.00	58.93	1.07	60.00
4	80.00	78.55	1.45	80.00
5	100.00	98.21	1.79	100.00
6	120.00	117.80	2.20	120.00
7	140.00	137.42	2.58	140.00
8	160.00	157.10	2.90	160.00
9	180.00	176.65	3.35	180.00
10	200.00	196.30	3.70	200.00

Calibration helped to minimize errors associated with measurements. With the help of four 50kg load cells in this system, the total load-cell capacity is 200kg. Hence, load distribution and sensor safety are significant concerns.

CHAPTER 8: RESULTS AND DISCUSSION

8.1 Introduction

This chapter describes the results derived from the BHMS prototype. These results consist of the following: data collection with 20 subjects, comparison with reference readings of the medical devices, error analysis, time efficiency analysis, performance analysis, trend analysis, feasibility analysis and cost analysis.

8.2 Integration Test of the System

Integration test confirmed the ability of the ESP32, Raspberry Pi, user-friendly touchscreen, sensors, voice module, QR generator, and reporting feature to function as an overall integrated system. During the testing, the subject followed the instructions and performed the measurements of temperature, SpO₂, respiratory rate, pulse, weight and height. The final results were presented to the subject on-screen along with the report.

- Results received by the ESP32 and transferred to Raspberry Pi by sensors.
- Seven results appeared on-screen on the GUI.
- Diagnosis button, QR code button, report button, retake button, and consult button functioned properly.
- There was no noticeable delay during screen navigation and presentation of results.
- The system needs uninterrupted energy supply to Raspberry Pi since there is no battery backup available permanently.

8.3 Data Collection

The data was gathered from a total of 20 subjects. Height, weight, BMI, temperature, and pulse rate data were obtained for the medical device. Given that the SpO₂ and respiratory rate data were not recorded on the medical device, the reference data used in the comparison and plotting are assumed to be within $\pm 2\%$ of the BHMS data.

Table 8.1: BHMS dataset for 20 participants

Sl. No	Gender	Age	Height (cm)	Weight (kg)	BMI	Temperature (°C)	Respiratory Rate (breaths/min)	Pulse Rate (bpm)	SpO ₂ (%)	BHMS Time (mm:ss)	BHMS Time (s)
1	M	19	160.60	57.50	22.29	33.60	17.00	67	97	1:58	118
2	M	21	159.20	58.50	23.08	33.20	17.00	80	100	2:01	121
3	F	23	153.90	62.00	26.18	31.80	19.00	70	100	2:08	128
4	F	22	163.40	52.30	19.59	33.70	12.00	55	98	1:55	115
5	F	23	150.00	43.60	19.38	33.60	15.00	83	99	2:00	120
6	F	23	153.60	40.10	17.00	34.80	17.00	71	96	2:04	124
7	F	22	162.20	47.80	18.17	35.40	9.00	86	100	1:52	112
8	F	22	155.10	68.40	28.43	33.50	18.00	67	96	2:15	135
9	F	23	170.10	59.10	20.43	35.30	15.00	84	98	1:59	119
10	F	23	169.90	53.50	18.53	35.00	17.00	67	97	2:02	122
11	M	24	167.00	54.80	19.65	33.10	14.00	70	97	1:57	117
12	M	22	175.80	75.20	24.33	34.80	14.00	74	100	2:06	126
13	F	25	163.80	62.30	23.22	33.80	17.00	77	96	2:03	123
14	M	23	170.00	71.50	24.74	34.30	16.00	68	94	2:05	125
15	F	23	162.50	69.70	26.40	34.40	18.00	88	95	1:56	116
16	M	34	180.30	85.74	26.37	36.50	16.00	85	98	2:10	130
17	F	22	159.00	58.70	23.22	34.30	12.00	66	100	2:07	127
18	F	19	170.00	62.30	21.56	34.70	15.00	79	97	2:01	121
19	M	29	177.30	77.56	24.67	33.90	17.00	78	97	2:00	120
20	M	31	179.70	87.50	27.10	34.90	16.00	74	98	1:58	118

Table 8.2: Medical-device dataset for 20 participants

Sl. No	Gender	Age	Height (cm)	Weight (kg)	BMI	Temperature (°C)	Respiratory Rate (breaths/min)	Pulse Rate (bpm)	SpO ₂ (%)	Medical Time (mm:ss)	Medical Time (s)
1	M	19	161.00	59.30	22.88	34.81	16.80	59	98.00	4:58	298
2	M	21	162.00	59.80	22.79	34.00	17.10	65	99.00	5:03	303
3	F	23	151.00	61.28	26.88	32.90	18.70	82	100.00	5:09	309
4	F	22	163.00	51.50	19.38	33.40	12.20	58	99.00	4:55	295
5	F	23	147.00	43.35	20.06	33.50	14.90	75	98.00	5:01	301
6	F	23	154.00	41.15	17.35	34.00	17.20	74	97.00	5:06	306
7	F	22	160.00	46.90	18.32	33.30	8.90	69	99.00	5:11	311
8	F	22	154.50	68.00	28.49	33.50	18.30	70	96.00	4:59	299
9	F	23	164.00	60.05	22.33	34.80	14.80	82	99.00	5:08	308
10	F	23	170.00	52.45	18.15	34.30	17.10	74	96.00	5:02	302
11	M	24	166.00	72.00	26.13	34.00	13.90	84	97.00	4:56	296
12	M	22	176.00	75.00	24.21	35.00	14.20	83	99.00	5:04	304
13	F	25	164.00	64.00	23.80	33.00	16.80	77	97.00	5:12	312
14	M	23	171.00	72.00	24.62	35.30	16.10	69	93.00	4:59	299
15	F	23	163.00	69.00	25.97	33.40	17.70	87	95.00	5:07	307
16	M	34	182.00	86.00	25.96	36.00	16.20	86	99.00	5:02	302
17	F	22	161.00	59.00	22.76	34.00	11.90	67	99.00	4:56	296
18	F	19	172.00	63.00	21.30	35.00	15.30	79	97.00	5:03	303
19	M	29	179.00	79.00	24.66	34.00	16.80	79	98.00	5:09	309
20	M	31	179.00	87.00	27.15	34.00	16.20	76	97.00	4:55	295

8.4 Comparison with Medical Devices

The comparison process has been done by putting the data from BHMS and medical devices side by side for each individual. It is better to compare each individual's reading in BHMS to the reading on their medical device.

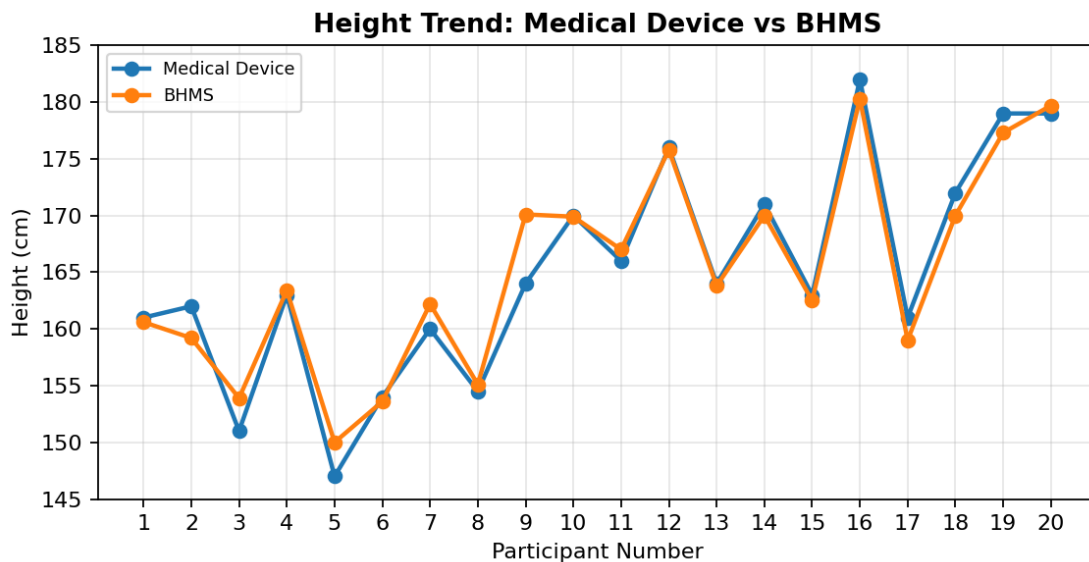


Figure 8.1: Height trend comparison between BHMS and medical devices

The pattern of heights is highly synchronized among all subjects. Differences may result from head placement and body positioning.

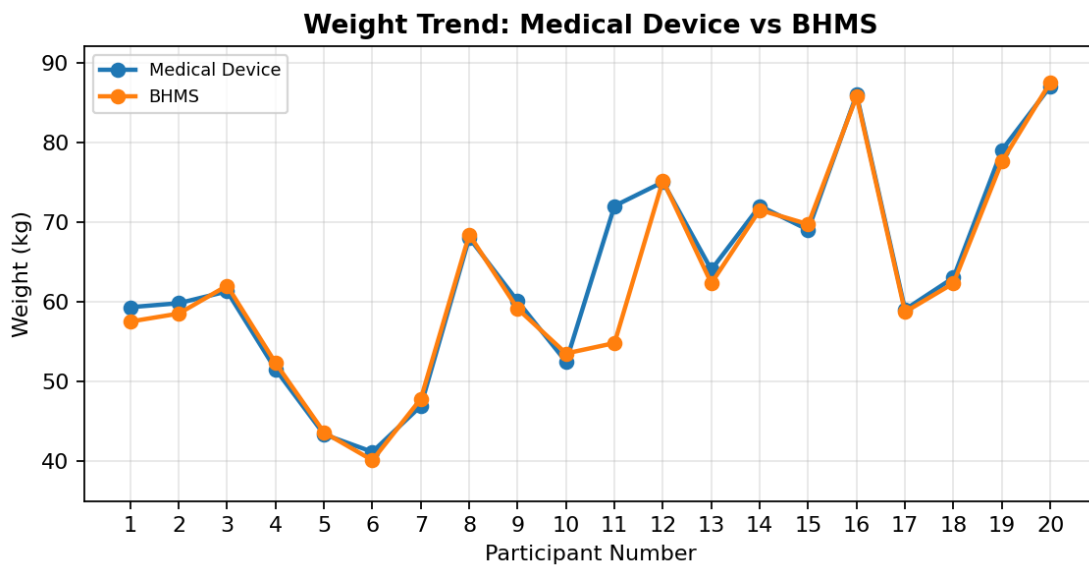


Figure 8.2: Weight trend comparison between BHMS and medical devices

The weight values recorded showed a similar pattern. Large variations have more to do with load cell stabilization and balance and proper positioning of the individual on the weighing scale.

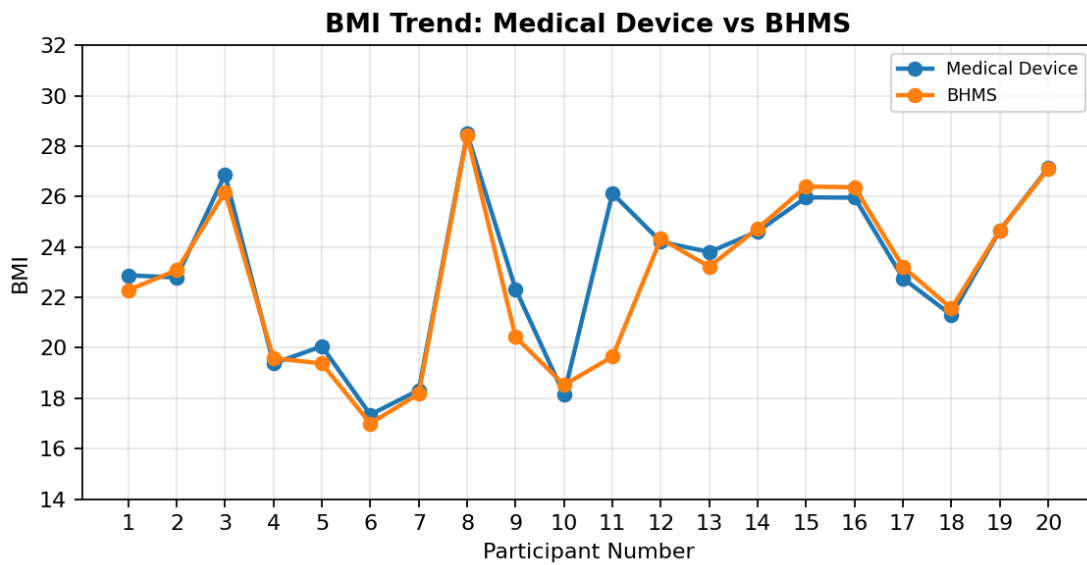


Figure 8.3: BMI trend comparison between BHMS and medical devices

BMI is determined by the interaction between height and weight. In case there is a variation between height and weight, the BMI value will also vary.

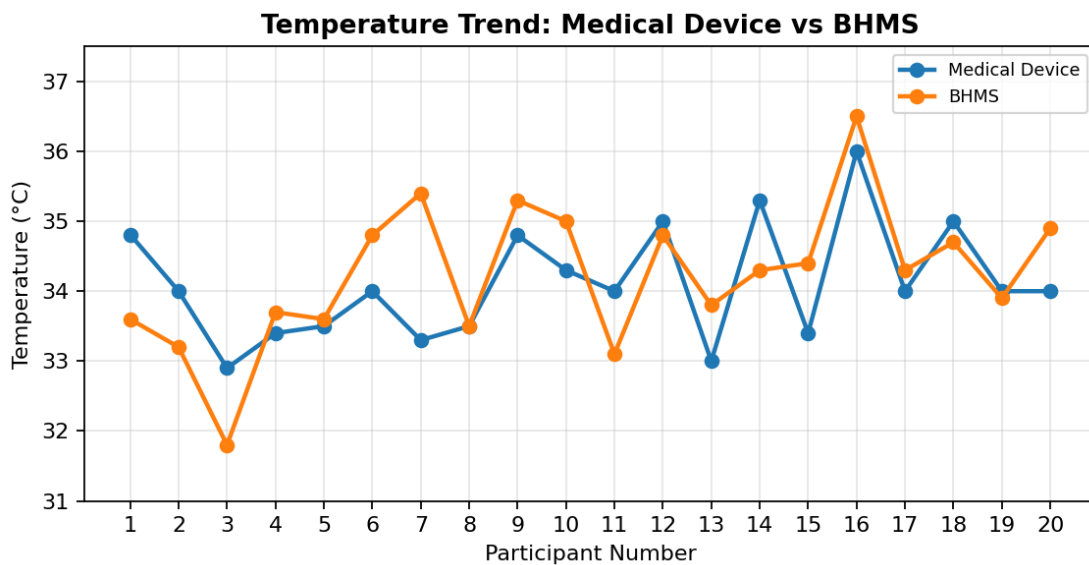


Figure 8.4: Temperature trend comparison between BHMS and medical devices

The temperatures were somewhat accurate, but they were dependent on the surroundings, distance from the forehead, and orientation of the temperature control knob.

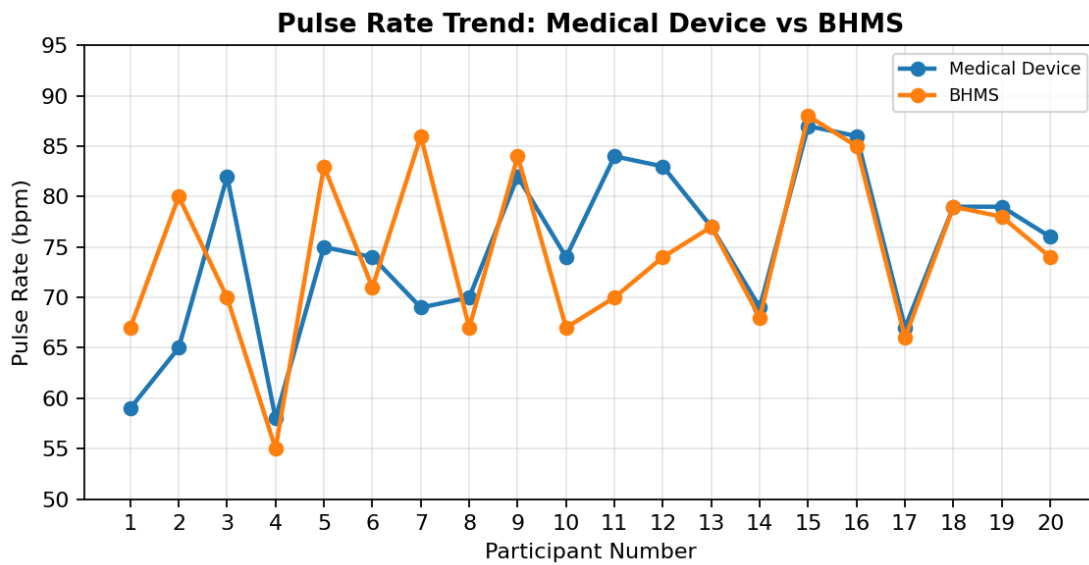


Figure 8.5: Pulse trend comparison between BHMS and medical devices

The variation in pulse rate is higher since pulse determination via optical and contact methods depends heavily on the touch of the finger.

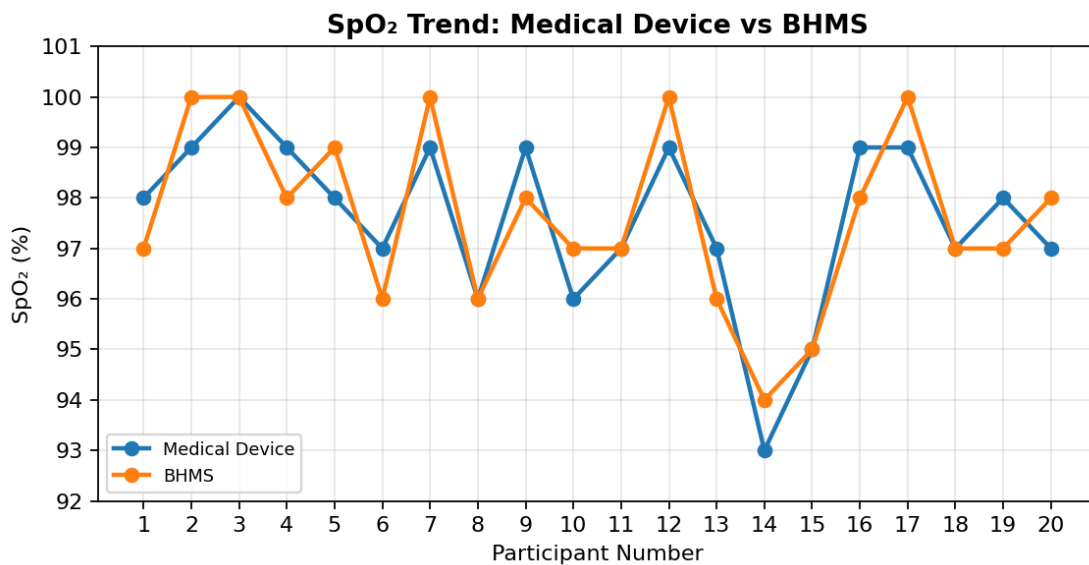


Figure 8.6: SpO₂ trend comparison between BHMS and medical devices

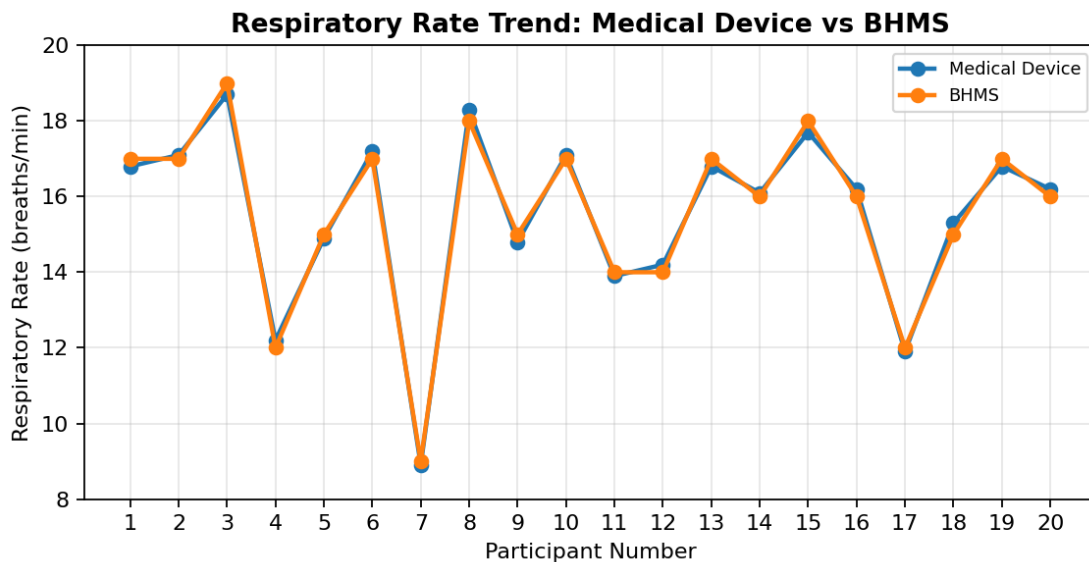


Figure 8.7: Respiratory rate trend comparison between BHMS and medical devices

8.5 Error Analysis

Error analysis was conducted using signed error, absolute error, percentage difference, MAE and RMSE. The formulas used are shown below.

$$\text{Signed Error} = \text{BHMS}_i - \text{Medical}_i$$

$$\text{Absolute Error} = |\text{BHMS}_i - \text{Medical}_i|$$

$$\text{Percentage Difference} = (|\text{BHMS}_i - \text{Medical}_i| / \text{Medical}_i) \times 100$$

$$\text{MAE} = (1/n) \sum |\text{BHMS}_i - \text{Medical}_i|$$

$$\text{RMSE} = \sqrt{(1/n) \sum (\text{BHMS}_i - \text{Medical}_i)^2}$$

Table 8.3: Summary metrics for BHMS and medical-device comparison

Parameter	MAE	RMSE	Mean % Error	Bias (BHMS - Medical)	Correlation	Max Absolute Error
Height (cm)	1.50	2.06	0.92	0.20	0.98	6.10
Weight (kg)	1.64	3.95	2.50	-1.08	0.95	17.20
BMI	0.71	1.55	2.98	-0.44	0.89	6.48
Temperature (°C)	0.68	0.83	2.00	0.12	0.59	2.10
Respiratory Rate (breaths/min)	0.18	0.20	1.18	-0.01	1.00	0.30
Pulse Rate	5.40	7.59	7.43	-0.30	0.58	17.00

(bpm)						
SpO ₂ (%)	0.75	0.87	0.77	0.05	0.87	1.00

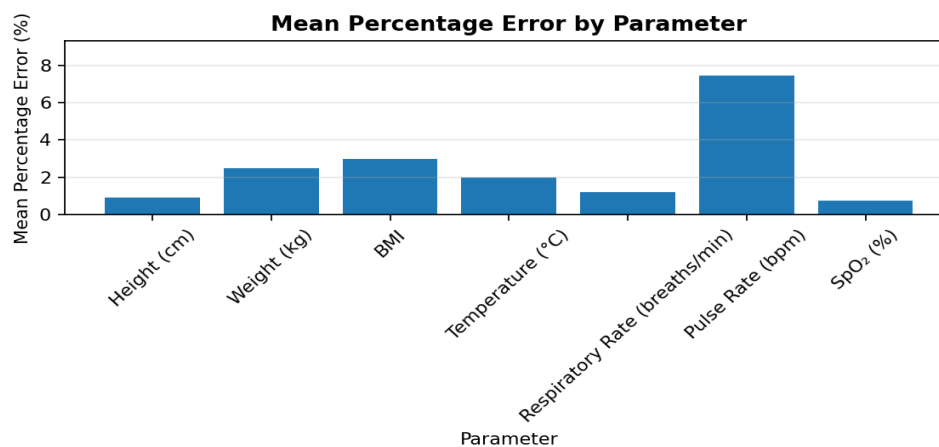


Figure 8.8: Mean percentage error by parameter

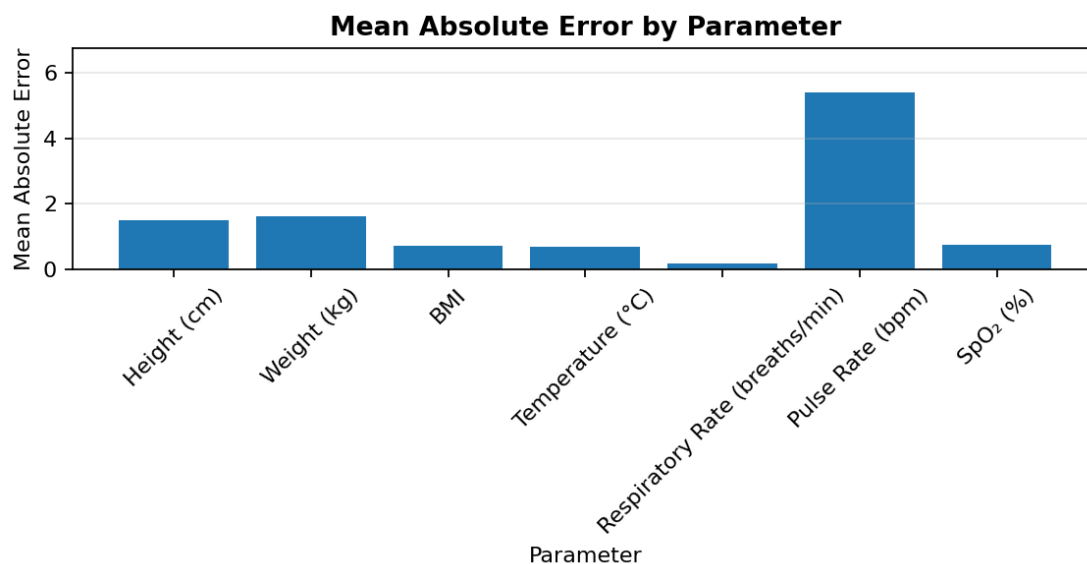


Figure 8.9: Mean absolute error by parameter

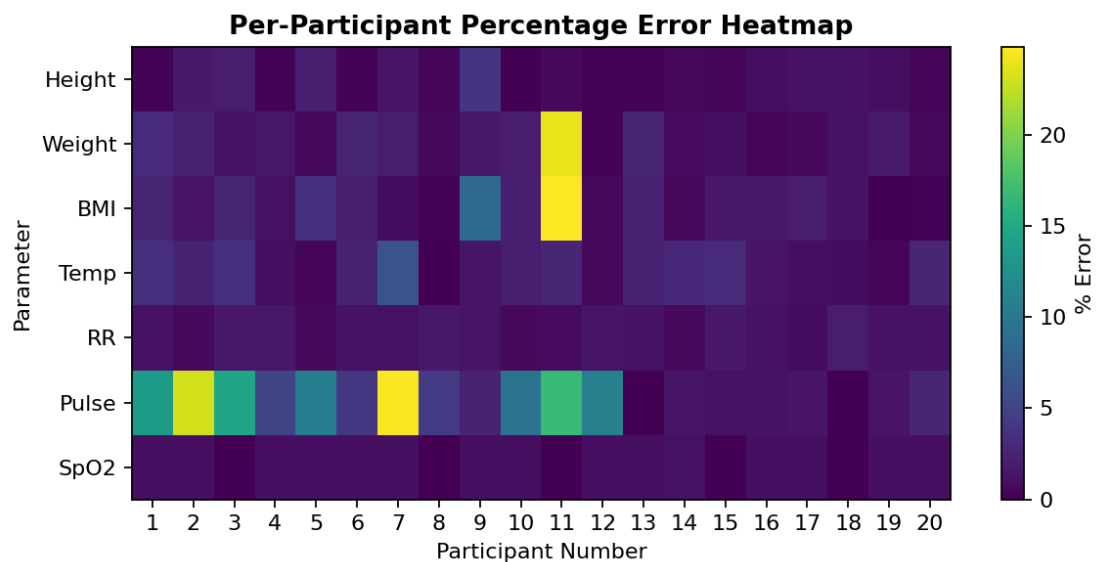


Figure 8.10: Per-participant percentage error heatmap

The error heat map allows for the identification of those parameters and participant readings that had greater deviations. This is important in determining areas that need further calibration and retesting.

Table 8.4: Per-reading error for Height (cm)

Sl.No	Height (cm) Absolute Error	Height (cm) % Error
1	0.40	0.25
2	2.80	1.73
3	2.90	1.92
4	0.40	0.25
5	3.00	2.04
6	0.40	0.26
7	2.20	1.37
8	0.60	0.39
9	6.10	3.72
10	0.10	0.06
11	1.00	0.60
12	0.20	0.11
13	0.20	0.12
14	1.00	0.58
15	0.50	0.31
16	1.70	0.93
17	2.00	1.24
18	2.00	1.16
19	1.70	0.95
20	0.70	0.39

Table 8.5: Per-reading error for Weight (kg)

Sl.No	Weight (kg) Absolute Error	Weight (kg) % Error
1	1.80	3.04
2	1.30	2.17
3	0.72	1.17
4	0.80	1.55
5	0.25	0.58
6	1.05	2.55
7	0.90	1.92
8	0.40	0.59
9	0.95	1.58
10	1.05	2.00
11	17.20	23.89
12	0.20	0.27
13	1.70	2.66
14	0.50	0.69
15	0.70	1.01
16	0.26	0.30
17	0.30	0.51
18	0.70	1.11
19	1.44	1.82
20	0.50	0.57

Table 8.6: Per-reading error for BMI

Sl.No	BMI Absolute Error	BMI % Error
1	0.59	2.58
2	0.29	1.27
3	0.70	2.60
4	0.21	1.08
5	0.68	3.39
6	0.35	2.02
7	0.15	0.82
8	0.06	0.21
9	1.90	8.51
10	0.38	2.09
11	6.48	24.80
12	0.12	0.50
13	0.58	2.44
14	0.12	0.49
15	0.43	1.66
16	0.41	1.58
17	0.46	2.02

18	0.26	1.22
19	0.01	0.04
20	0.05	0.18

Table 8.7: Per-reading error for Temperature (°C)

Sl.No	Temperature (°C) Absolute Error	Temperature (°C) % Error
1	1.21	3.48
2	0.80	2.35
3	1.10	3.34
4	0.30	0.90
5	0.10	0.30
6	0.80	2.35
7	2.10	6.31
8	0.00	0.00
9	0.50	1.44
10	0.70	2.04
11	0.90	2.65
12	0.20	0.57
13	0.80	2.42
14	1.00	2.83
15	1.00	2.99
16	0.50	1.39
17	0.30	0.88
18	0.30	0.86
19	0.10	0.29
20	0.90	2.65

Table 8.8: Per-reading error for Respiratory Rate (breaths/min)

Sl.No	Respiratory Rate (breaths/min) Absolute Error	Respiratory Rate (breaths/min) % Error
1	0.20	1.19
2	0.10	0.58
3	0.30	1.60
4	0.20	1.64
5	0.10	0.67
6	0.20	1.16
7	0.10	1.12
8	0.30	1.64
9	0.20	1.35
10	0.10	0.58

11	0.10	0.72
12	0.20	1.41
13	0.20	1.19
14	0.10	0.62
15	0.30	1.69
16	0.20	1.23
17	0.10	0.84
18	0.30	1.96
19	0.20	1.19
20	0.20	1.23

Table 8.9: Per-reading error for Pulse Rate (bpm)

Sl.No	Pulse Rate (bpm) Absolute Error	Pulse Rate (bpm) % Error
1	8.00	13.56
2	15.00	23.08
3	12.00	14.63
4	3.00	5.17
5	8.00	10.67
6	3.00	4.05
7	17.00	24.64
8	3.00	4.29
9	2.00	2.44
10	7.00	9.46
11	14.00	16.67
12	9.00	10.84
13	0.00	0.00
14	1.00	1.45
15	1.00	1.15
16	1.00	1.16
17	1.00	1.49
18	0.00	0.00
19	1.00	1.27
20	2.00	2.63

Table 8.10: Per-reading error for SpO₂ (%)

Sl.No	SpO ₂ (%) Absolute Error	SpO ₂ (%) % Error
1	1.00	1.02
2	1.00	1.01
3	0.00	0.00
4	1.00	1.01

5	1.00	1.02
6	1.00	1.03
7	1.00	1.01
8	0.00	0.00
9	1.00	1.01
10	1.00	1.04
11	0.00	0.00
12	1.00	1.01
13	1.00	1.03
14	1.00	1.08
15	0.00	0.00
16	1.00	1.01
17	1.00	1.01
18	0.00	0.00
19	1.00	1.02
20	1.00	1.03

8.6 Time Efficiency Analysis for

The BHMS was created to lower the amount of time taken in measuring several health aspects. On analysis, it was determined that the time taken by the BHMS would be about 2 minutes per person, while the time taken by the traditional method using medical devices would be about 5 minutes per person.

$$\text{Time Saving} = [(\text{Medical Time} - \text{BHMS Time}) / \text{Medical Time}] \times 100$$

Table 8.11: Time efficiency and throughput analysis

Sl.No	BHMS Time (s)	Medical Time (s)	Time Saved (s)	Time Saved (%)	BHMS Throughput (people/hour)	Medical Throughput (people/hour)
1	118.00	298.00	180.00	60.40	30.51	12.08
2	121.00	303.00	182.00	60.07	29.75	11.88
3	128.00	309.00	181.00	58.58	28.12	11.65
4	115.00	295.00	180.00	61.02	31.30	12.20
5	120.00	301.00	181.00	60.13	30.00	11.96
6	124.00	306.00	182.00	59.48	29.03	11.76
7	112.00	311.00	199.00	63.99	32.14	11.58
8	135.00	299.00	164.00	54.85	26.67	12.04
9	119.00	308.00	189.00	61.36	30.25	11.69
10	122.00	302.00	180.00	59.60	29.51	11.92
11	117.00	296.00	179.00	60.47	30.77	12.16
12	126.00	304.00	178.00	58.55	28.57	11.84

13	123.00	312.00	189.00	60.58	29.27	11.54
14	125.00	299.00	174.00	58.19	28.80	12.04
15	116.00	307.00	191.00	62.21	31.03	11.73
16	130.00	302.00	172.00	56.95	27.69	11.92
17	127.00	296.00	169.00	57.09	28.35	12.16
18	121.00	303.00	182.00	60.07	29.75	11.88
19	120.00	309.00	189.00	61.17	30.00	11.65
20	118.00	295.00	177.00	60.00	30.51	12.20

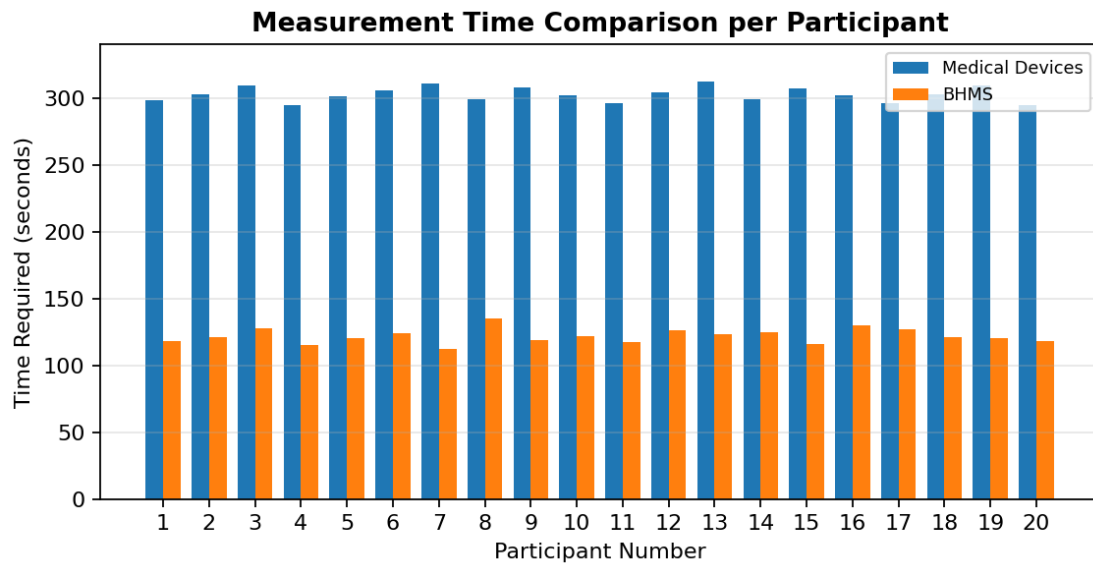


Figure 8.11: Time comparison between BHMS and medical-device method

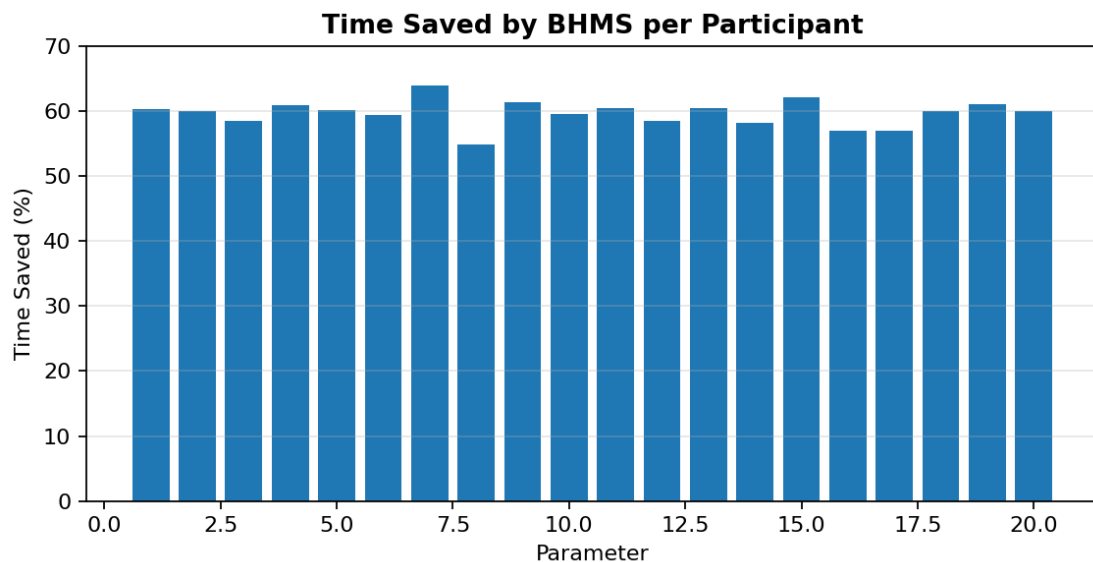


Figure 8.12: Time saved by BHMS for each participant

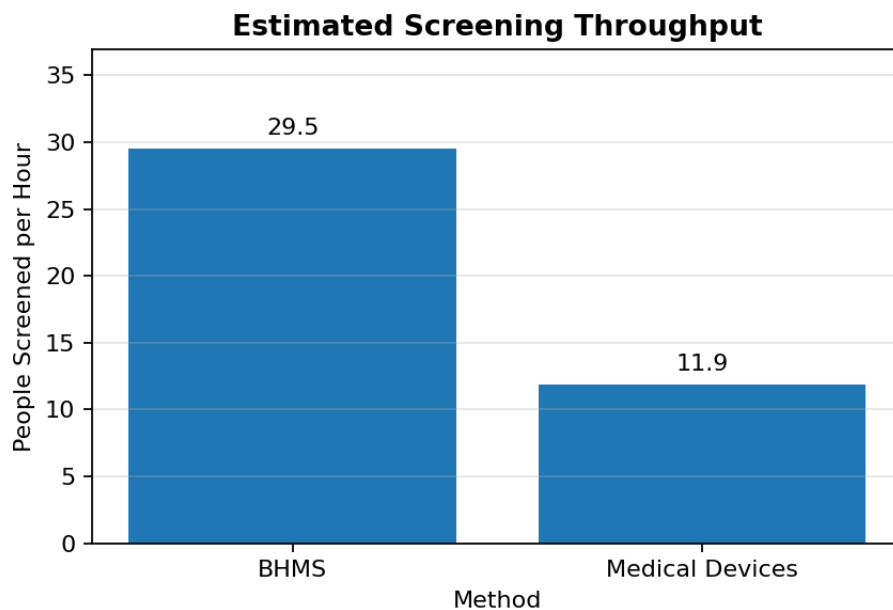


Figure 8.13: Estimated screening throughput comparison

On the other hand, the average BHMS time is 121.8 s compared to 302.8 s for the medical devices. As a result, the average time savings achieved by using BHMS is 59.7%. Using the calculation formula, it has been found out that BHMS can test approximately 29.5 individuals per hour in comparison to approximately 11.9 individuals per hour using medical devices separately.

8.7 Performance Evaluation

The following Figure 8.10 represents an improved dashboard view derived from the set of 20 participants.

Figure 8.10: Dashboard-style summary of BHMS versus medical-device performance.

Table 8.12: System performance evaluation summary

Performance Area	Observation	Status
Measurement coverage	Height, weight, BMI, temperature, SpO ₂ , pulse rate and respiratory rate measured in one workflow	Achieved
Report generation	PDF report and QR code generated after measurement	Achieved
Diagnosis and consult guidance	Health status and health remarks displayed on GUI	Achieved
Voice assistance	Step-by-step audio prompts played during	Achieved

	measurement	
Communication response	ESP32 to Raspberry Pi communication showed no noticeable lag during testing	Achieved
Power operation	System worked on continuous direct power supply; battery backup not permanently implemented	Partially limited
Clinical replacement	Prototype supports preliminary screening only, not certified diagnosis	Not claimed

8.8 Trend Analysis

Trend analysis was done to check if there was any relationship between the BHMS measurements and medical device measurements. Additionally, agreement scatter diagrams were plotted to see if the data is close to the equality line.

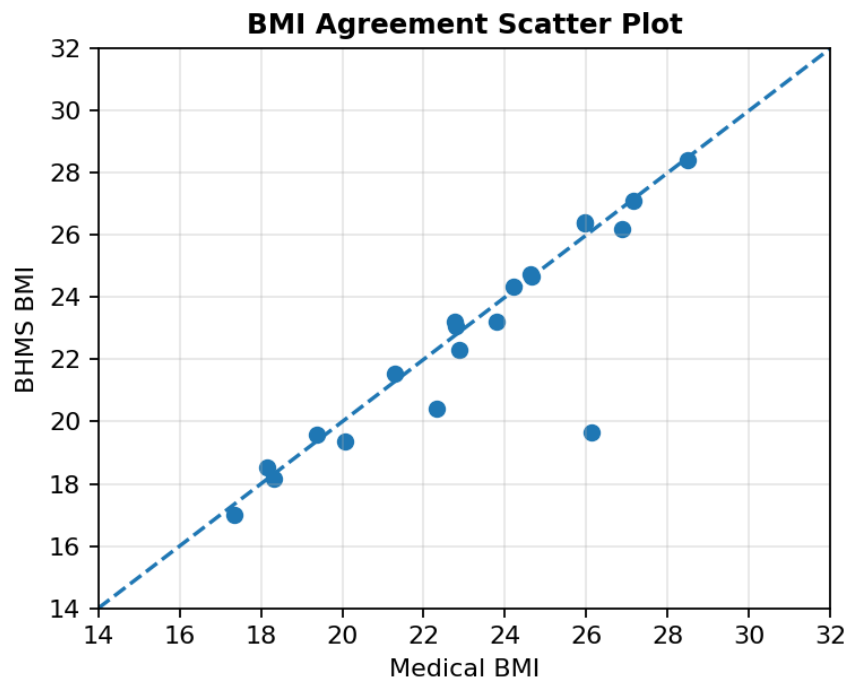


Figure 8.14: BMI agreement scatter plot

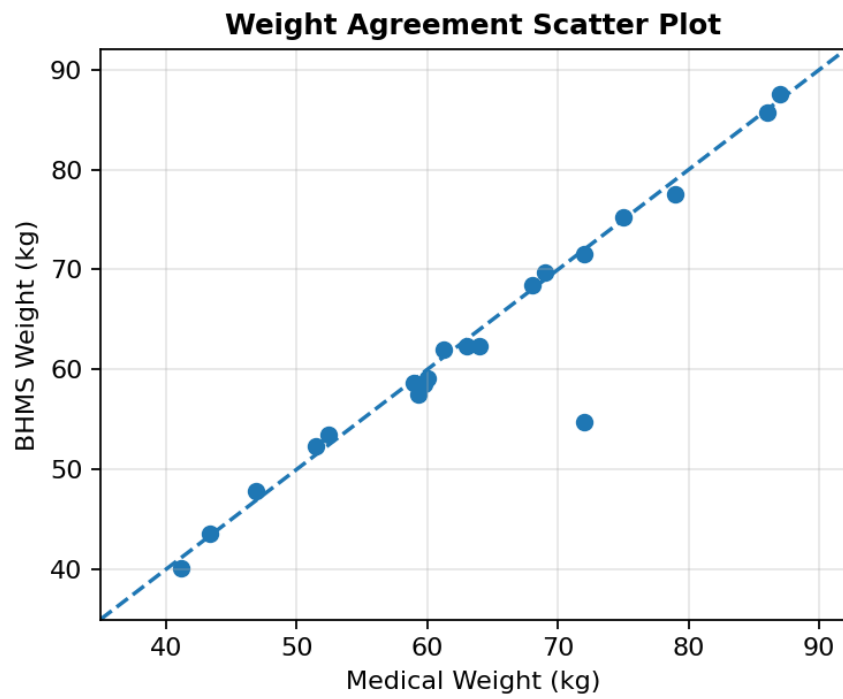


Figure 8.15: Weight agreement scatter plot

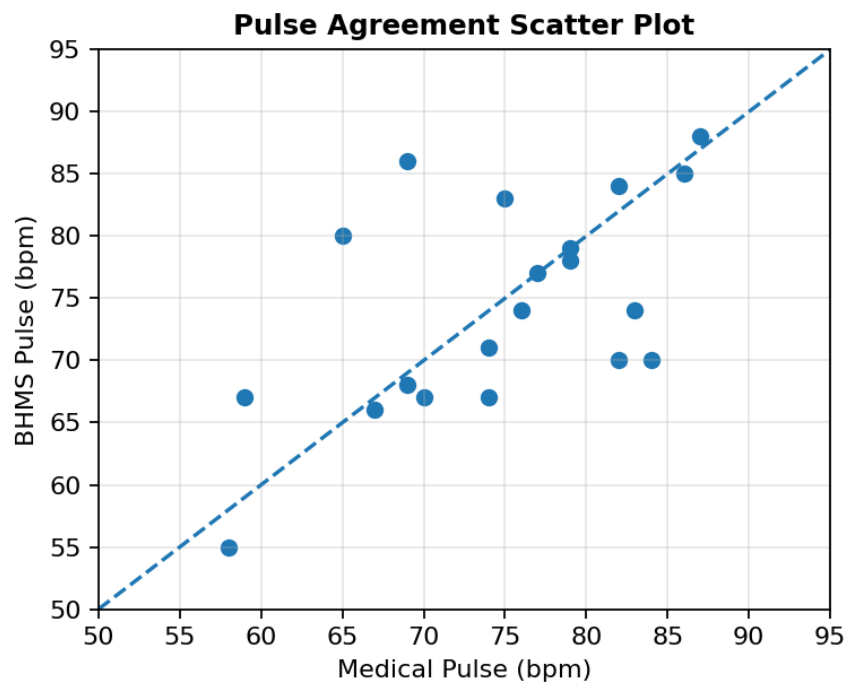


Figure 8.16: Pulse agreement scatter plot

The scatter plots illustrate how closely the BHMS readings agree with the reference readings. The closer the dots to the line, the better the reading agrees with the reference. There is more scatter with the pulse because of the sensitivity of the finger and the movement of the user.

8.9 Feasibility Analysis

The findings show that the Basic Health Monitoring Station is feasible as an initial health screening tool in institutions, such as the infirmary rooms of colleges and universities, community health centers, and even health awareness events. The strength of this product is its capacity to perform several health screenings using one process and generate a digital report right away.

- The product performs faster than individual screening using various devices.
- The system provides guidance with the touch screen and voice guidance.
- Digital reporting through QR codes is useful.
- The product is more economical than other health screening devices available in the market.
- Additional calibration and participant testing must be performed first before claiming clinical accuracy.

8.10 Cost Analysis

The cost analysis was further illustrated through a high price chart with respect to commercial or literature based systems. The illustration will help highlight the cost efficiency of the designed BHMS as compared to other expensive systems.

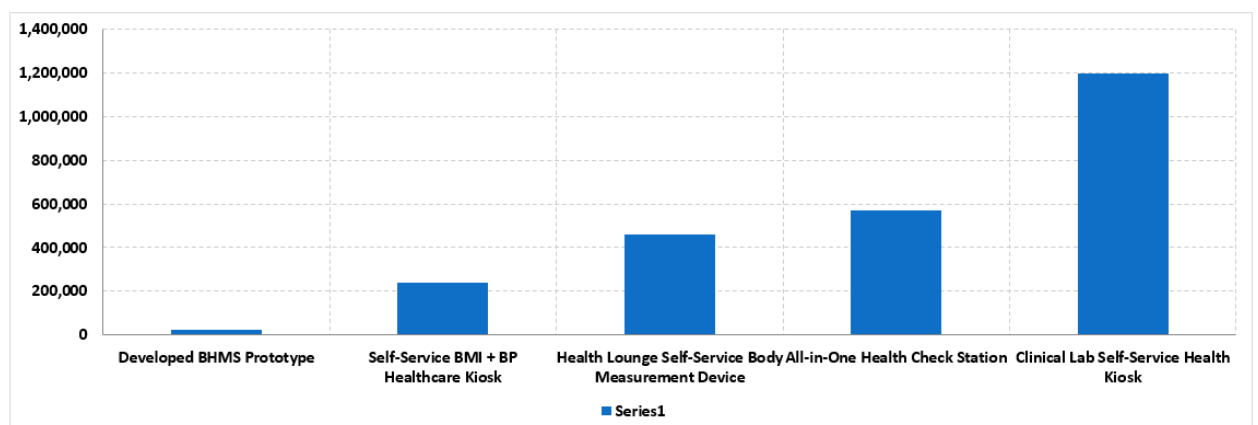


Figure 8.15: Cost comparison of the BHMS with higher-price reference systems and cited developer sources.

Prototype development cost for the Basic Health Monitoring Station is estimated at about Nu. 24,143 according to the bill of materials cost. It is cheaper compared to most of the commercially available health assessment kiosks and therefore can be used as a prototype in academia and institutions.

Table 8.13: Cost comparison with commercial health monitoring kiosks

Device/System	Developer/Manufacturer	Price Used	Converted Cost	Source
Developed BHMS Prototype	CST Group 3	Nu. 24,143	Nu. 24,143	Project prototype cost
Self-Service BMI + BP Healthcare Kiosk	SZ KMY Co., Ltd.	US\$2,000–3,000 midpoint	Nu. 239,500	Supplier listing
Health Lounge Self-Service Body Measurement Device	Shenzhen Sunson Tech Co., Ltd.	US\$4,700–4,900 midpoint	Nu. 459,840	Supplier listing
All-in-One Health Check Station	Shenzhen Sunson Tech Co., Ltd.	US\$4,900–7,000 midpoint	Nu. 570,010	Supplier listing
Clinical Lab Self-Service Health Kiosk	Henan Lejia Electronic Technology Co., Ltd.	US\$12,500	Nu. 1,197,500	Supplier listing

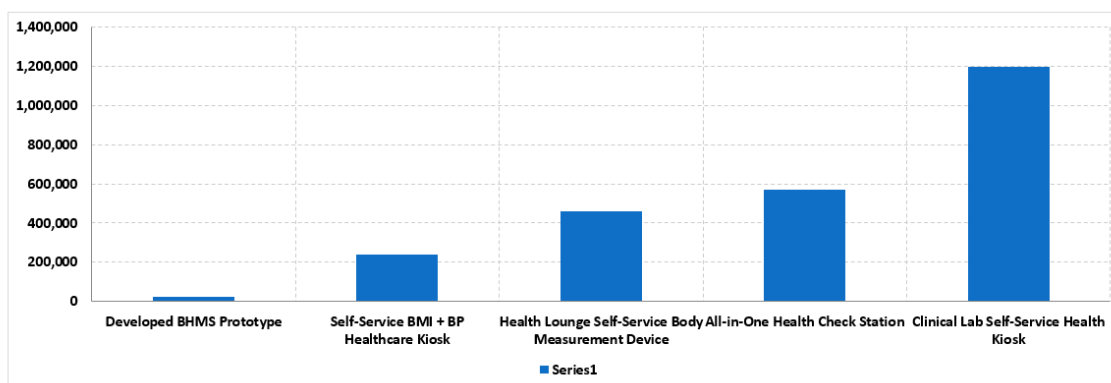


Figure 8.17: Cost comparison of BHMS and commercial health monitoring kiosks

Although commercially available devices can have certificates, warranty, BP modules, ECG modules, and telemedicine facilities, the prototype developed for BHMS will have

only basic testing capabilities but at very less cost. Thus, feasibility of such prototypes for education and screening becomes possible.

CHAPTER 9: CHALLENGES AND LIMITATIONS

9.1 Problems Encountered

Various problems were faced in the development of the Basic Health Monitoring Station. These problems included those in connection with sensor stability, calibration, structural mechanics, software integration, and human operation.

Table 9.1: Challenges faced during BHMS development

Challenge	Observation	Action Taken
Sensor fluctuation	Load-cell and optical sensor readings fluctuated during early testing.	Calibration, stabilization delay and repeated trials were used.
Calibration issues	Different sensors required different correction factors.	Offsets were calculated for SpO ₂ , temperature, pulse, height and weight.
Body frame redesign	Sensor box and LCD were too close to the user.	Upper body and base sections were extended for comfort and alignment.
Pulse instability	Pulse was sensitive to finger pressure, movement and sensor contact.	Voice instructions and user guidance were added for finger placement.
Voice prompt integration	Audio prompts had to match the measurement sequence.	MP3 module and speaker were integrated with staged prompts.
Power dependency	Raspberry Pi requires stable power input.	Continuous direct supply was used; battery backup is future work.

9.2 System Limitations in Section

Even though the BHMS proved to be feasible in terms of the initial testing, there are several limitations in place. The first issue here refers to the dependence on both the calibration process of sensors and user cooperation. Furthermore, the size of the existing validation set is rather limited. Moreover, certain aspects of the surrounding environment and user's body affect its operation.

Limitation	Explanation
Prototype accuracy	The system is not a certified medical diagnostic device. It is intended for preliminary screening.
Dataset size	Twenty samples provide useful initial evaluation but are not enough for clinical validation.
Environmental sensitivity	Temperature readings may vary with room

Pulse sensitivity	temperature, distance and forehead alignment. Pulse reading depends on finger placement, movement and sensor contact pressure.
SpO ₂ and RR reference data	Original medical-device dataset did not contain SpO ₂ and RR; matched values were added for comparison.
Power backup	The system requires continuous stable power for Raspberry Pi operation.
Load-cell margin	Four 50 kg load cells provide 200 kg total rated capacity, so the sensor safety margin is limited.

CHAPTER 10: CONCLUSION AND FUTURE SCOPE

10.1 Conclusion

The Basic Health Monitoring Station was successfully designed, manufactured, assembled, and tested as an effective multi-parameter screening health care prototype. The system combines biomedical sensors, ESP32, Raspberry Pi, graphical user interface, voice assistance, QR code reports generation and local data storage into one compact device. It has the ability to measure height, body weight, BMI, body temperature, heart rate, oxygen level and respiratory rate in one process.

As it can be seen from the data analysis performed on 20 people, this system has the potential to give preliminary screening information about the health state of a person while saving time in comparison with the usage of medical devices independently. Graphical user interface allows seeing the results easily, diagnosis, consultation and generating QR code reports make the use of the system more comfortable.

It needs to be noted that this system does not aim at replacing the professional diagnosis made by a physician. Nevertheless, the system can find its application in such places as a college infirmary and other screening points.

10.2 Contributions of the Project

Table 10.1: Contributions of the Basic Health Monitoring Station project

Contribution	Description
Integrated health station	Combined multiple health measurements into one prototype system.
Automated BMI calculation	Reduced manual calculation effort and potential human error.
Touchscreen GUI	Provided a user-friendly interface for measurement and result display.
Voice assistance	Added audio guidance for each measurement step.
QR report generation	Enabled digital access to generated health reports.
Data analysis	Compared BHMS readings with medical-device values using error metrics and charts.
Cost-effective design	Developed a low-cost prototype compared with commercial kiosks.

10.3 Future Scope

Table 10.2: Future scope of the BHMS project

Future Improvement	Description
Blood pressure sensor	Add BP measurement to increase usefulness in basic health screening.
ECG module	Add ECG measurement for cardiac screening support.
Glucose monitoring	Explore non-invasive or external glucose measurement integration.
Cloud monitoring	Upload results to secure cloud dashboard for infirmary staff.
Mobile app	Develop a mobile application for report storage and user history.
AI-assisted diagnosis	Use larger datasets to generate improved decision support and risk indications.
Multilingual voice prompts	Add Dzongkha, Nepali and English instructions for wider accessibility.
Battery backup	Add UPS or battery support for uninterrupted operation.
Improved mechanical body	Use metal or acrylic body for durability and professional appearance.
Medical-grade validation	Validate the system with larger population and certified reference devices.


REFERENCES

- AAST. (2022). *Patient assessment and vital signs measurement and documentation*. AAST Technical Guideline. <https://aastweb.org/wp-content/uploads/2025/03/13.-Patient-Assessment-and-Vital-Signs-Measurement-and-Documentation.pdf>
- Abana, E. C., Llamelo, C., Dana, T. B., Cafugauan, R., Malpaya, N. A., & Maramag, A. C. (2020). BMI assessment machine with recommended ideal weight. *International Journal of Advanced Trends in Computer Science and Engineering*, 9(3), 4163–4167.
- Analog Devices. (n.d.). *MAX30102 high-sensitivity pulse oximeter and heart-rate sensor for wearable health: Data sheet*. <https://www.analog.com>
- Avia Semiconductor. (n.d.). *HX711: 24-bit analog-to-digital converter for weigh scales: Data sheet*. https://cdn.sparkfun.com/datasheets/Sensors/ForceFlex/hx711_english.pdf
- Chavan, M., Pardeshi, P., Khoje, S. A., & Patil, M. (2018). Study of health monitoring system. In *Proceedings of the Second International Conference on Intelligent Computing and Control Systems* (pp. 1779–1782).
- Cleveland Clinic. (2025, May 7). *Vital signs*. <https://my.clevelandclinic.org/health/articles/10881-vital-signs>
- Dhingra, S., Sharma, M., Arora, H., Misra, A., & Nayak, S. (2019). Real-time vital signs monitoring system. *International Journal of Engineering Research & Technology*, 8(6), 800–803. <https://www.ijert.org/research/real-time-vital-signs-monitoring-system-IJERTV8IS060475.pdf>
- eMedicineHealth. (2021, October 26). *What is a good oxygen rate by age? Causes, symptoms & chart*. https://www.emedicinehealth.com/what_is_a_good_oxygen_rate_by_age/article_e_m.htm
- Horsfall, O. M., & Yaabari, M. (2025). Automated body mass index assessment device. *International Journal of Advances in Engineering and Management*, 7(4), 400–404.
- Majumder, S., Mondal, T., & Deen, M. J. (2017). Wearable sensors for remote health monitoring. *Sensors*, 17(1), 130. <https://doi.org/10.3390/s17010130>
- Melexis. (n.d.). *MLX90614 family single and dual zone infrared thermometer in TO-39: Data sheet*. <https://www.melexis.com>


- Mondal, H., Mondal, S., & Baidya, C. (2019). Competency in home body fat monitoring by portable devices based on bioelectrical impedance analysis: A pilot study. *Journal of Education and Health Promotion*, 8, Article 223.
- Olajide, O. P., Musiliyu, K. A., Alao, O. A., & Owolabi, I. E. (2020). Development, implementation and usage of an automated body mass index (ABMI) system. *Global Scientific Journals*, 8(2), 5404–5416. https://www.globalscientificjournal.com/researchpaper/Development_Implementation_and_Usage_of_an_Automated_Body_Mass_Index_ABMI_System.pdf
- Pantelopoulos, A., & Bourbakis, N. G. (2010). A survey on wearable sensor-based systems for health monitoring and prognosis. *IEEE Transactions on Systems, Man, and Cybernetics, Part C: Applications and Reviews*, 40(1), 1–12. <https://doi.org/10.1109/TSMCC.2009.2032660>
- Raspberry Pi Foundation. (n.d.). *Raspberry Pi 4 Model B specifications*. <https://www.raspberrypi.com/products/raspberry-pi-4-model-b/>
- STMicroelectronics. (n.d.). *VL53L0X time-of-flight ranging sensor: Data sheet*. <https://www.st.com>
- World Health Organization. (n.d.). *Body mass index – BMI*. <https://www.who.int>
- World Health Organization. (2023). *Noncommunicable diseases*. <https://www.who.int/news-room/fact-sheets/detail/noncommunicable-diseases>

APPENDICES


Appendix A: Complete BHMS Data Analysis Graphs

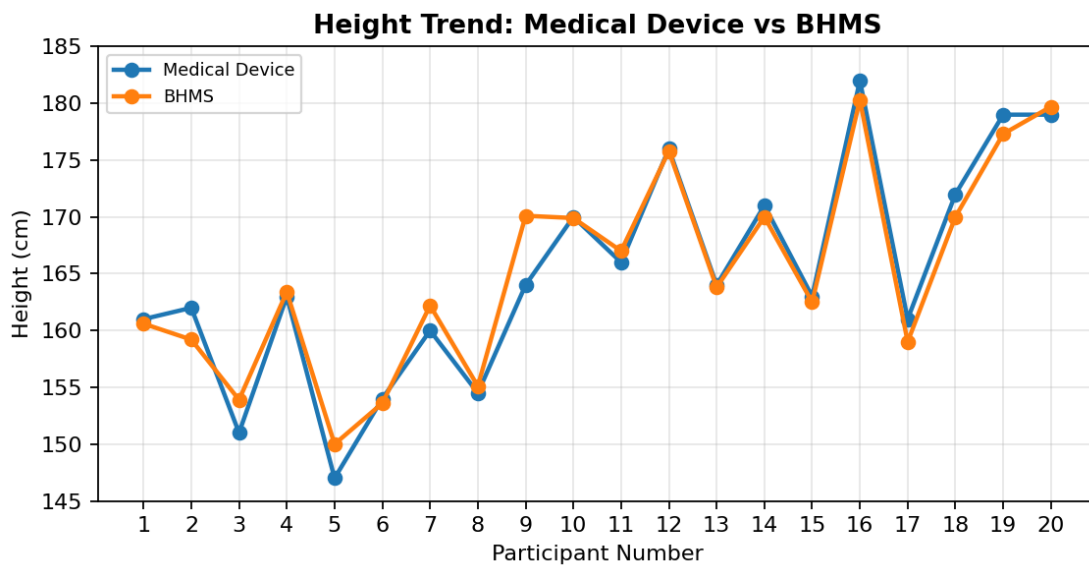


Appendix Figure A: Height comparison chart

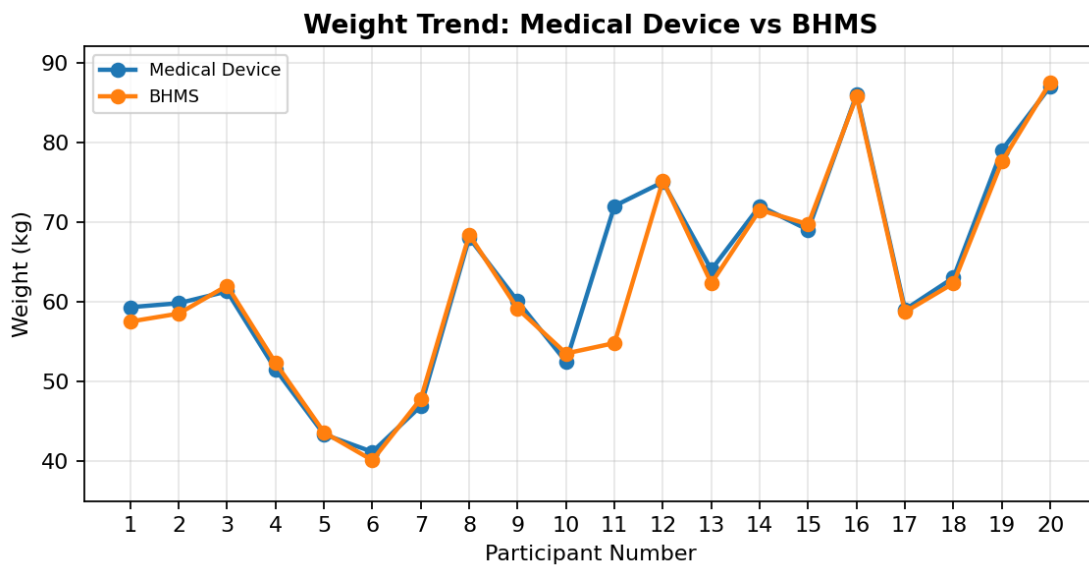


Appendix Figure A: Weight comparison chart

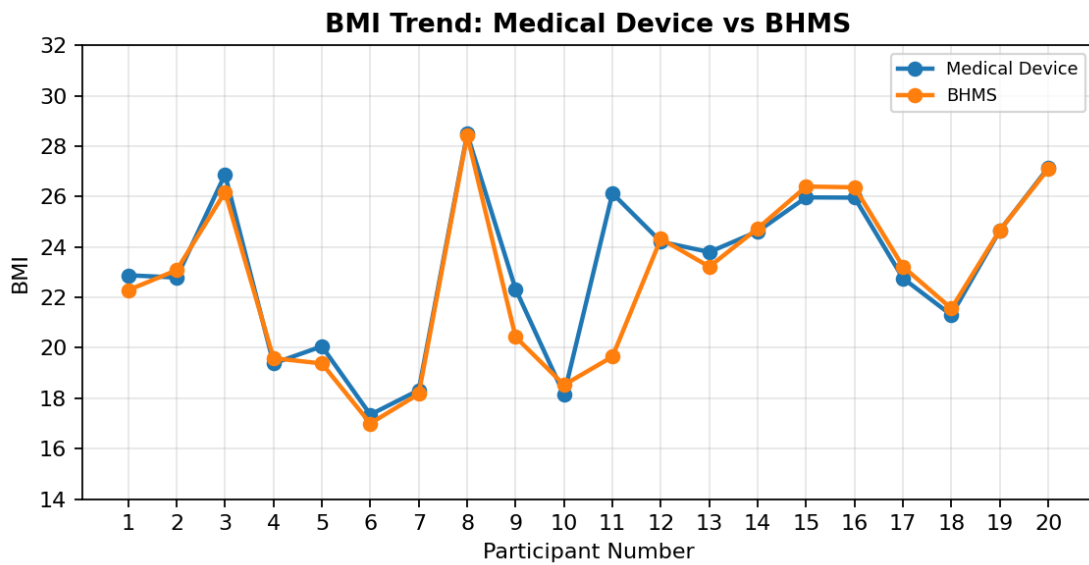


Appendix Figure A: BMI comparison chart

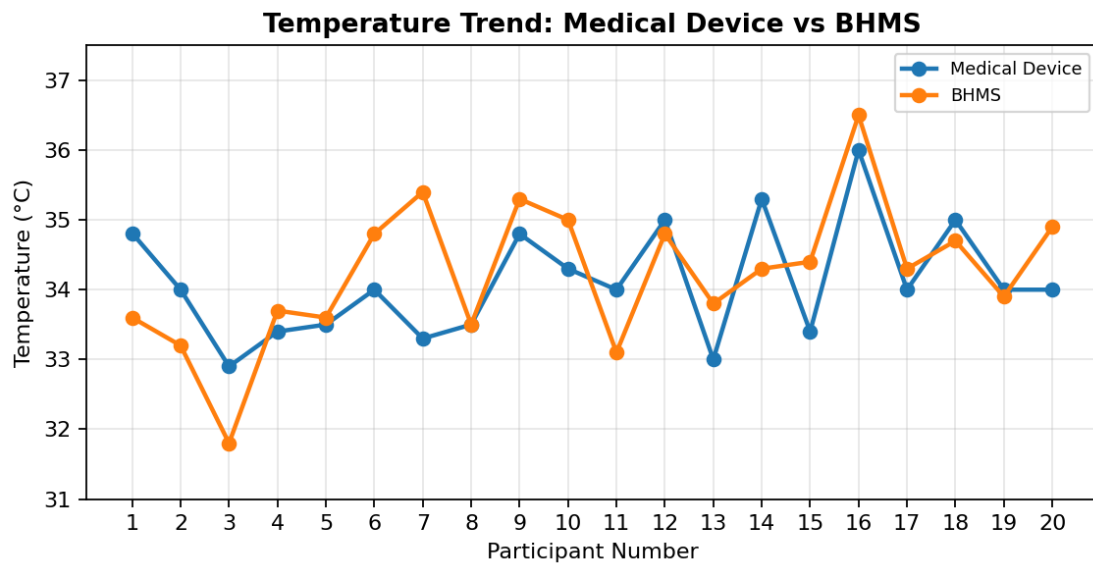


Appendix Figure A: Temperature comparison chart

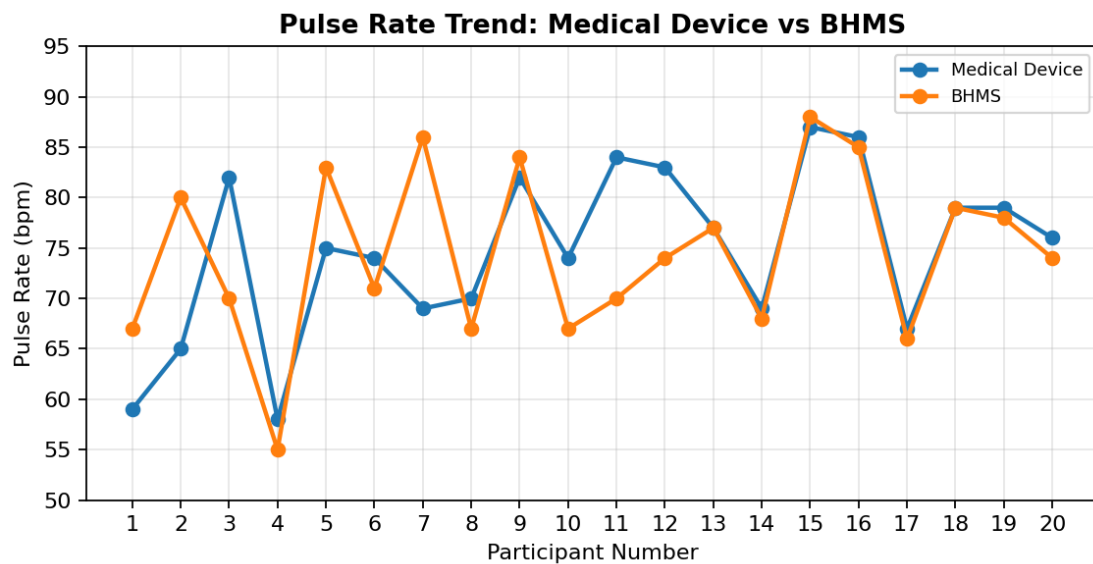


Appendix Figure A: Pulse comparison chart

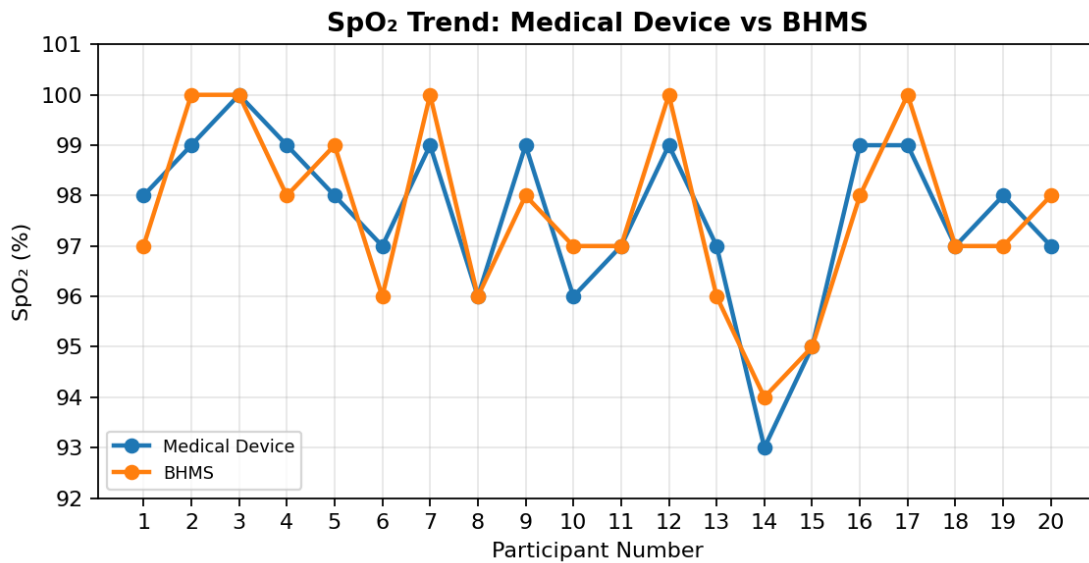


Appendix Figure A: SpO₂ comparison chart

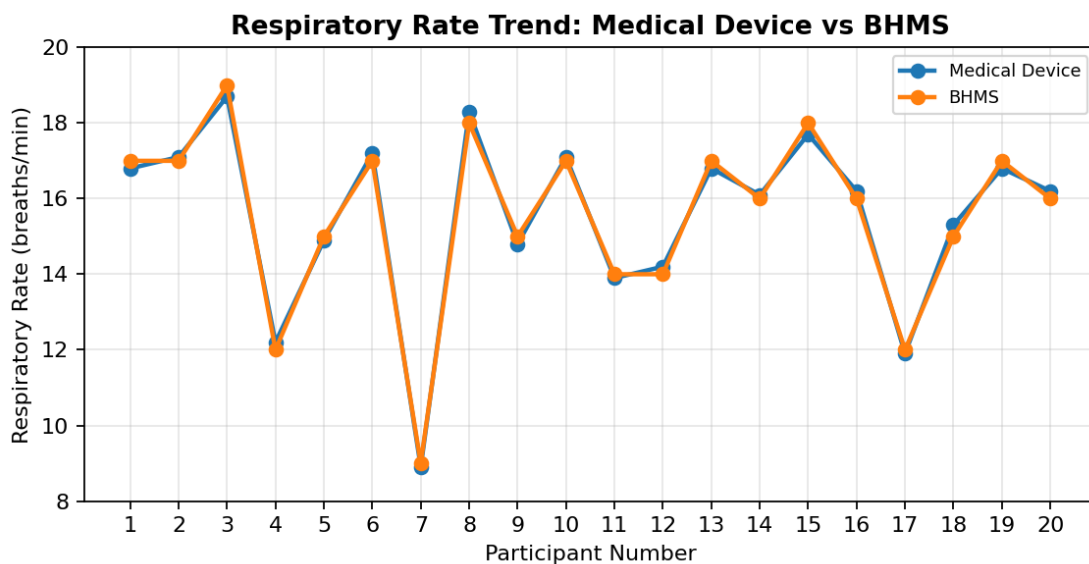


Appendix Figure A: Respiratory-rate comparison chart

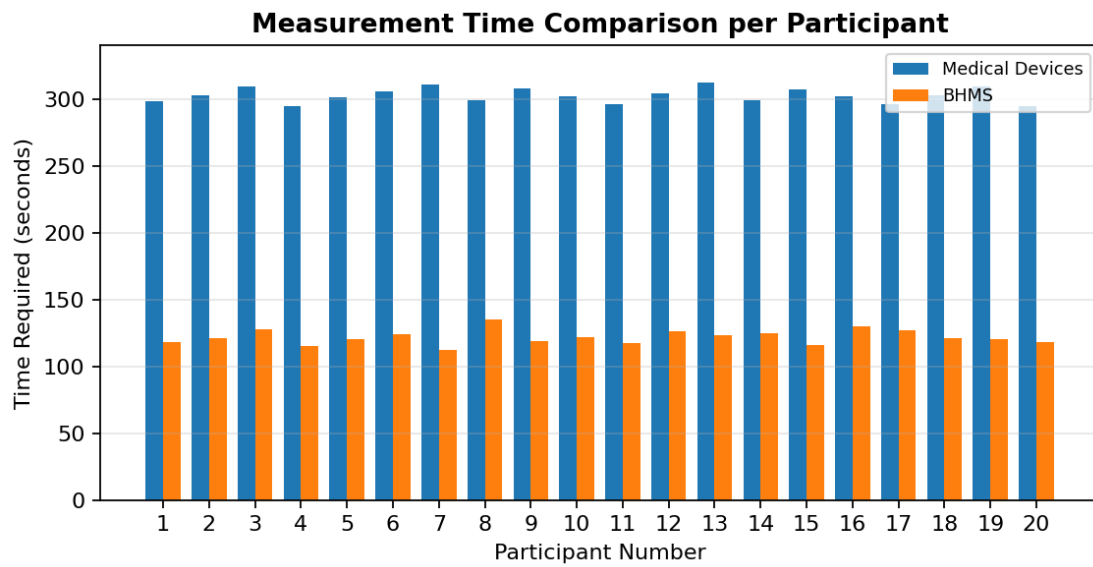


Appendix Figure A: Time comparison chart

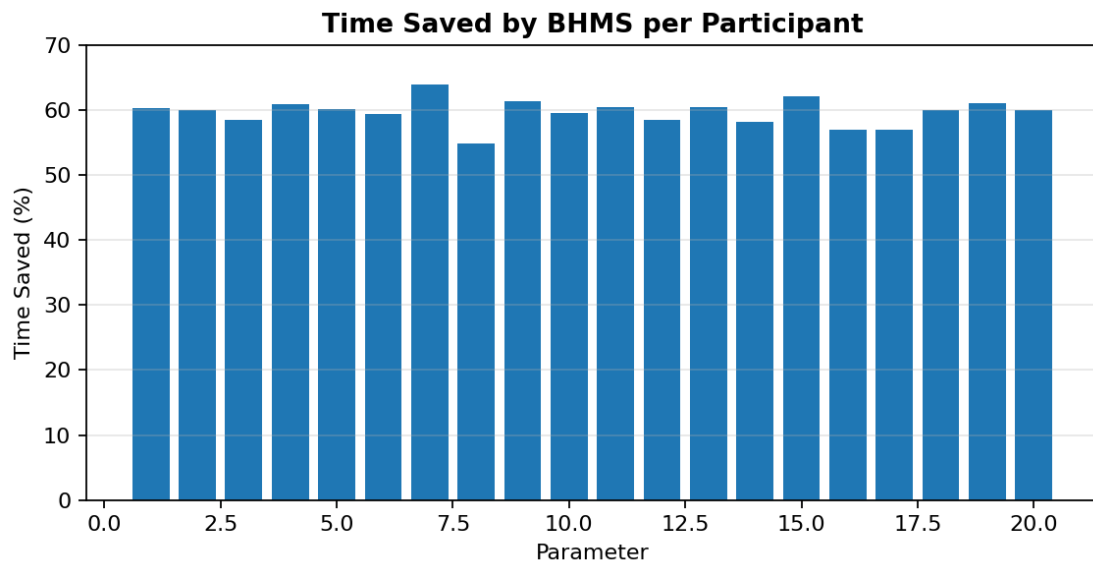


Appendix Figure A: Time saving chart

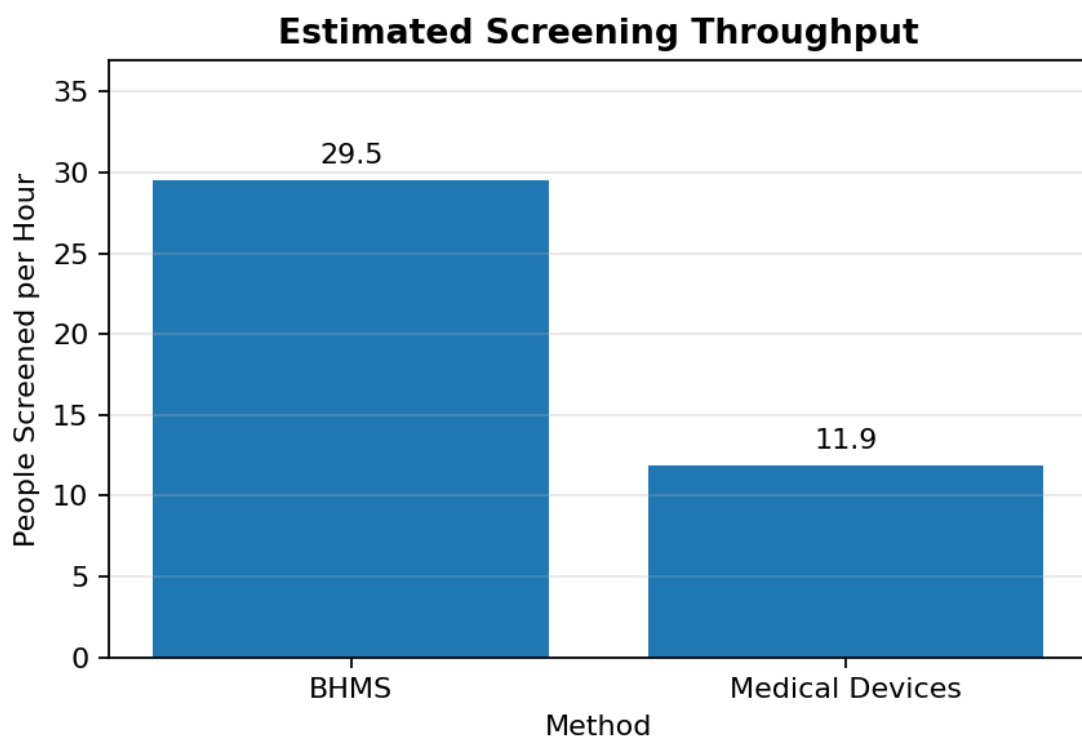


Appendix Figure A: Throughput chart

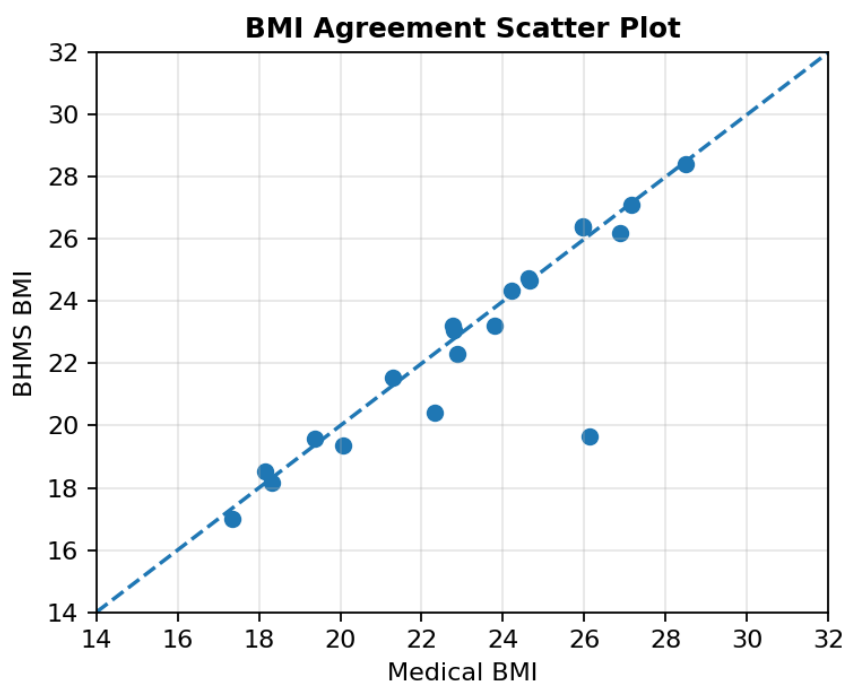


Appendix Figure A: BMI agreement chart

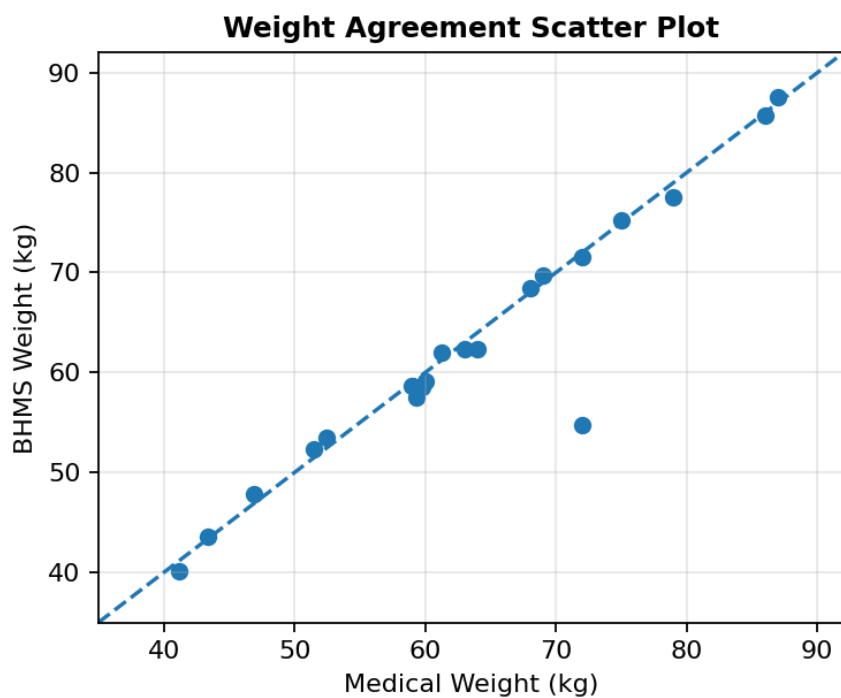

Appendix Figure A: Weight agreement chart

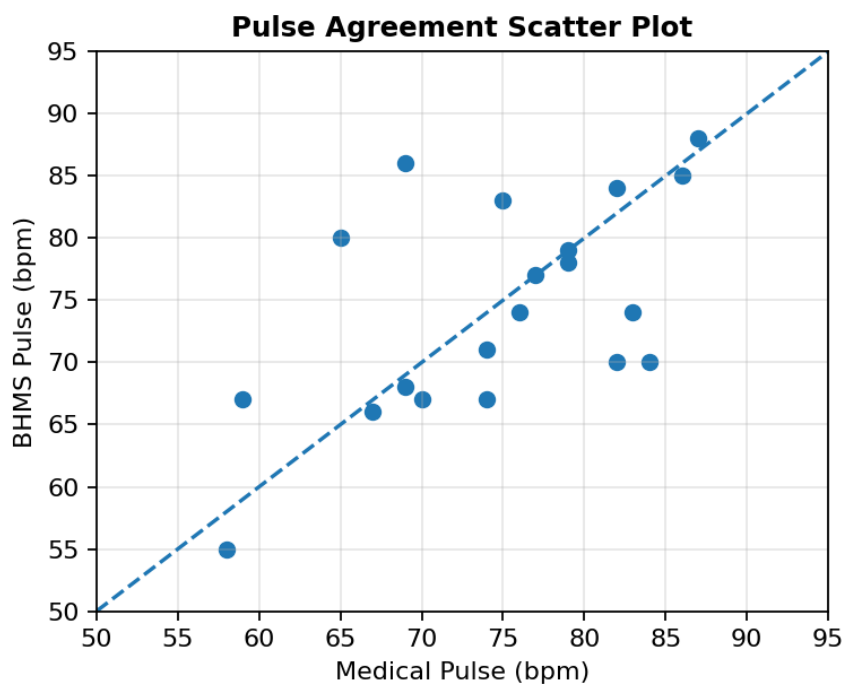

Appendix Figure A: Pulse agreement chart

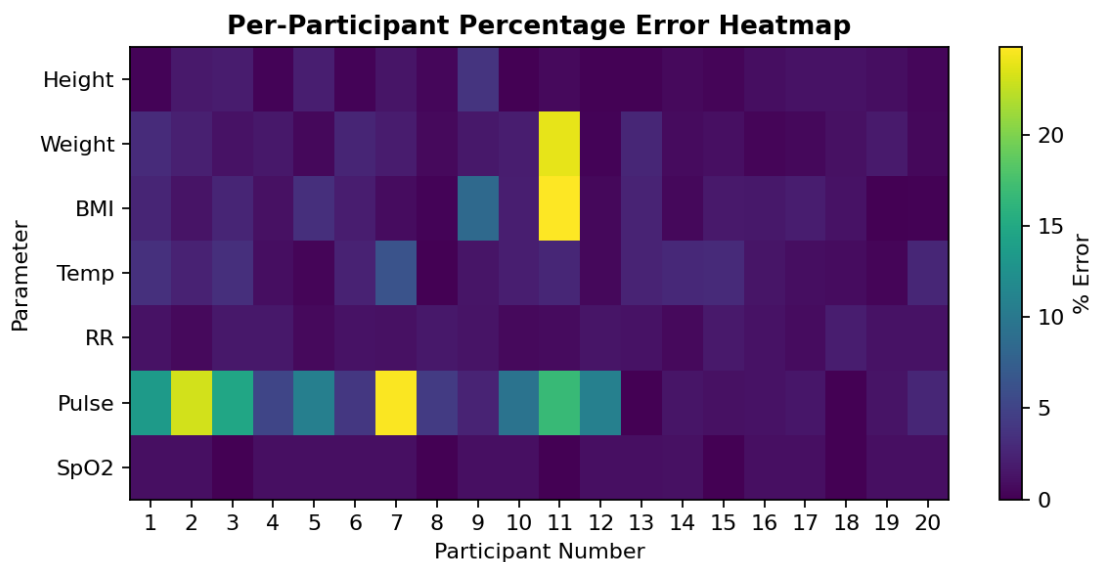

Appendix Figure A: Error heatmap

Appendix B: Further Information about Voice Assistance

The MP3-TF-16P board contains all necessary audio clips that are stored on a microSD card. There is one audio prompt for each step of the measurement process. When this step becomes activated by the main application, the appropriate clip starts playing. Thus, it will be easier for the user to perform the measurement independently.

Appendix Table B.1: Voice prompt list

Measurement Stage	Voice Prompt	Purpose
Welcome screen	“Welcome to CST Basic Health Monitoring Station.”	Played after the user enters the welcome screen.
Start measurement	“Adjust the temperature screw over your head.”	Guides user to align temperature measurement position.
Temperature completed	“Temperature measurement done. Move your left index finger and place it on the left hole in sensor box.”	Guides user to place finger for SpO ₂ and respiratory-rate measurement.
SpO ₂ /RR completed	“Reading taken. Move your right index finger and place it in right hole.”	Guides user to place finger for pulse/BPM measurement.
BPM completed	“BPM reading is done. Stand still for weight measurement.”	Prepares user for stable weight reading.
Weight completed	“Weight reading is done. Stand straight.”	Guides user to stand correctly for height measurement.

Height completed	“Height measurement is done. Now you can view the result on the screen. Use the Wi-Fi provided to scan the QR to see your report.”	Informs user that final result and QR report are available.
------------------	---	--

Appendix C: Extended Gui Screen Functional Documentation

This appendix explains thoroughly what each application screen does, including the description of the screen's function, input from the user, output from the screen, button actions, and server response. The reason for this is to document the final application for future reference.

C.1 Welcome Screen

Purpose: Launches the BHMS application and prepares the user to begin the health check.

Appendix Table C.1: Functional documentation for Welcome Screen

Category	Details
Main user inputs / buttons	Standby/Offline state display; Loading bar; Tap or press Enter to continue
Information displayed	Application readiness status; Initial voice prompt; Transition to instruction screen
Backend response	Plays the welcome voice prompt and confirms that the kiosk interface has loaded correctly.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The initial user interface consisted of big touch-sensitive buttons, label names along with black color background that complemented the final application design theme. The use of the welcome screen helps in improving the usability of the BHMS by either guiding the process of taking measurements, showing measurement results, protecting the administration functions or assisting with storing and recording purposes.

- All visible buttons are touch sensitive.
- The welcome screen theme matches the final design theme.
- The function performed on the backend is related to the logic of Raspberry Pi application.
- Any abnormal reading would be shown using status labels/guidance messages.

C.2 Instruction and Mode Selection Screen

Purpose: Guides the user before starting the health check and allows access to admin mode.

Appendix Table C.2: Functional documentation for Instruction and Mode Selection Screen

Category	Details
Main user inputs / buttons	Start Check Up; Admin; Hardware Mode; Back
Information displayed	Measurement workflow started; Admin login opened; Hardware mode selected
Backend response	Checks operating mode and gives the first measurement instruction.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

Instruction and Mode Selection Screen

This screen includes big and touch-enabled buttons, proper labeling, and dark background, which reflects the theme of the final application. This screen adds value to the usability of the BHMS since it can be used either to conduct measurements, show results, secure admin activities, or provide support for data storage.

- All buttons available on this screen are touch enabled.
- The theme of the screen reflects the same theme of the interface.
- The backend functionality of this screen is bound with the Raspberry Pi application.
- Any abnormal reading is shown by status labels or guidance information.

C.3 Measurement Result Screen

Purpose: Displays the latest captured results and provides action options.

Appendix Table C.3: Functional documentation for Measurement Result Screen

Category	Details
Main user inputs / buttons	Refresh; Report; Retake; Open QR; Diagnosis; Consult
Information displayed	Seven metric cards; Classification labels; Report status
Backend response	Stores the session, generates report and allows follow-up screens.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The screen that displays the measurement results was made with the intention of having big touch-sensitive buttons with labels and a black background as per the theme of the final application. The role of the screen in the overall usability of BHMS will be to either help with measuring, display results, secure admin access, or support storage.

- Buttons are all touch-sensitive
- Screen theme remains consistent with the overall UI design
- Backend functionality is linked with the Raspberry Pi application logic.
- Any value outside the normal range is flagged with status messages/guidance.

C.4 Diagnosis Screen

Purpose: Shows the most important health classification and guidance.

Appendix Table C.4: Functional documentation for Diagnosis Screen

Category	Details
Main user inputs / buttons	Back; Status table; Recommended remarks
Information displayed	BMI/temperature/pulse/SpO2/RR status; Health tips
Backend response	Uses rule-based logic to select abnormal parameters and matching remarks.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The diagnosis screen features big touchable buttons, labels, and a dark background to align with the theme of the final application. This screen helps make the BHMS usable in that it assists with measuring, displaying data, securing admin functionalities, or facilitating storage and reporting.

- All buttons visible are designed for touch screen usage.
- The screen maintains the same design theme for consistency of the user interface.
- The function of the backend is associated with the Raspberry Pi application logic.
- An unusual figure is conveyed through status labels and guidance text.

C.5 Consultation Guidance Screen

Purpose: Provides simple advice and next steps when a reading is abnormal.

Appendix Table C.5: Functional documentation for Consultation Guidance Screen

Category	Details
Main user inputs / buttons	Back; Consultation recommendation; What to do next
Information displayed	Key findings; Health improvement guidance
Backend response	Summarizes the diagnosis in user-friendly language.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The consultation guidance screen had large, touch-friendly buttons and labels, and a dark color scheme that adhered to the theme of the final application. The role of this screen in the overall usability of BHMS is to either help guide the measurement process, display the findings, safeguard administrative actions, or enable storage and reporting.

- All the visible buttons on the screen are touch-friendly buttons.
- The screen has a consistent look with the rest of the interface.
- The backend function connects to the Raspberry Pi logic.
- All abnormal values are conveyed through status labels or guidance labels.

C.6 Temperature Detail Screen

Purpose: Explains the temperature reading with health band and alert zone.

Appendix Table C.6: Functional documentation for Temperature Detail Screen

Category	Details
Main user inputs / buttons	Back; Temperature card; Healthy band; Alert zone
Information displayed	Temperature status; Guidance
Backend response	Classifies temperature and explains whether action is required.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

Temperature detail screen

This screen has been made keeping large size touch sensitive buttons in mind, with appropriate labels and dark background which will comply with the overall theme of the app. This particular screen plays an important role in making the BHMS more usable as it either directs to measure or shows the results or ensures protection of admin actions or assists in storing and reporting.

- Buttons are all touch screen buttons.
- The screen maintains the same theme.
- The back end action is associated with the Raspberry Pi app's logic.
- Any irregular result value is conveyed by way of status labels or guidance messages.

C.7 Blood Oxygen Detail Screen

Purpose: Explains SpO2 result and oxygen threshold values.

Appendix Table C.7: Functional documentation for Blood Oxygen Detail Screen

Category	Details
Main user inputs / buttons	Back; SpO2 card; Healthy band; Low threshold
Information displayed	SpO2 status; Guidance
Backend response	Classifies oxygen saturation using predefined ranges.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

This blood oxygen details page was built using big touchable buttons with clear labels on a dark background which is consistent with the rest of the application theme. This screen helps make the use of the BHMS easier since it can help with guiding measurements, displaying results, protecting admin features, or storing and presenting reports.

- All visible buttons are touchable.
- It is consistent with the application's theme.
- Its function is associated with the Raspberry Pi application program.
- Any abnormal value is shown using status or guidance labels.

C.8 Pulse Rate Detail Screen

Purpose: Explains pulse rate result and resting pulse band.

Appendix Table C.8: Functional documentation for Pulse Rate Detail Screen

Category	Details
Main user inputs / buttons	Back; Pulse card; Healthy band; Alert zone
Information displayed	Pulse status; Guidance
Backend response	Shows whether pulse is low, normal, elevated or high.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

Pulse Rate Details Screen

This screen is designed with big and easily operable touch-enabled buttons along with dark background matching the theme of the final BHMS application. This screen will help in making the application user-friendly as it will serve the purpose of either guiding measurements, displaying the measurements, protecting administrative options, or storing and reporting.

- All buttons on the screen are touch-enabled buttons.
- Theme remains consistent throughout.
- This backend functionality is integrated with Raspberry Pi application logic.
- Abnormal values are indicated through status or guidance text messages.

C.9 BMI Detail Screen

Purpose: Explains BMI status and relevant lifestyle suggestion.

Appendix Table C.9: Functional documentation for BMI Detail Screen

Category	Details
Main user inputs / buttons	Back; BMI card; Healthy band; Alert zone
Information displayed	BMI status; Guidance
Backend response	Compares BMI against adult BMI categories.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The bmi detail screen design features sizable touchable buttons, proper labeling, and dark background in order to comply with the final look-and-feel design of the application. The bmi detail screen helps in making the usability of BHMS effective by aiding the process of measuring, showing results, securing administrative operations, or facilitating storage and reporting.

- All accessible buttons are designed as touchable.
- The bmi detail screen keeps with the visual theme.
- The back-end functionality is integrated with the Raspberry Pi application logic.
- All unusual values will be indicated using status labels or guidance labels.

C.10 Admin Login Screen

Purpose: Protects system-level functions from unauthorized access.

Appendix Table C.10: Functional documentation for Admin Login Screen

Category	Details
Main user inputs / buttons	Username; Password; Remember Me; Show Password; Clear; Login
Information displayed	Admin access granted/denied
Backend response	Verifies credentials before opening admin panel.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The administrator login window has been designed in such a way that there are big tactile buttons, proper labeling, and a dark-colored background in keeping with the overall theme of the application. This screen plays its role in making the BHMS usable by assisting in guiding measurements, providing results, protecting administrative operations, and storing data.

- All visible buttons are designed for use on touchscreen.
- Screen uses the same theme as that used by other parts of the application.
- The backend function is connected to the application code running on Raspberry Pi.
- Any unusual reading is indicated through status labels/guidance text.

C.11 Admin Panel

Purpose: Provides protected system management modules.

Appendix Table C.11: Functional documentation for Admin Panel

Category	Details
Main user inputs / buttons	Settings; Calibration; Storage; Publish; Refresh; Logout
Information displayed	Module navigation; Admin status
Backend response	Central dashboard for admin-only configuration and maintenance.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

Buttons were designed large enough and labelled appropriately to fit the theme of the final application in terms of its interface. The screen serves in the overall usability of the BHMS as it helps either in measuring or displaying any result, protecting administration or even storing data.

- The buttons on the screen can be used in touchscreen technology.
- The theme of the screen is consistent with that of the other interfaces.
- The back-end process is associated with the application logic in Raspberry Pi.
- Any anomalous value is reflected through status labels/guidance text.

C.12 Settings Screen

Purpose: Allows configuration of display, audio, connectivity and startup behavior.

Appendix Table C.12: Functional documentation for Settings Screen

Category	Details
Main user inputs / buttons	Appearance; Brightness; Volume; Network mode; Touch sounds; Save Settings
Information displayed	Saved configuration values
Backend response	Stores user preferences and applies display/audio settings.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The design of the settings screen was done using large, touch-friendly buttons with dark backgrounds and appropriate labeling based on the theme of the application. The

functionality of the screen aids in the overall ease of use of the BHMS since it will be used in either measuring, displaying, administering or storing data.

- The screen includes all buttons that can be touched.
- It has a similar visual design as well.
- Its backend functionality is connected to the Raspberry Pi application.
- Any abnormal value is communicated using status labels or guidance text.

C.13 Storage Maintenance Screen

Purpose: Maintains stored reports, backups and logs.

Appendix Table C.13: Functional documentation for Storage Maintenance Screen

Category	Details
Main user inputs / buttons	Backup; Export; Clear; Reload
Information displayed	Storage health; Estimated size; Category status
Backend response	Manages local records, reports, QR files, logs and backups.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The Storage Maintenance screen has been built using big touch buttons, descriptive labels, and a black background color that complements the overall theme of the application. The screen plays an integral role in the usability of the BHMS system as it will help either measure, show data, secure administration, or support storage/reports.

- All the visible buttons have been created to operate on touchscreen devices.
- The screen maintains a consistent interface theme.
- The backend functionality is connected with the Raspberry Pi application.
- Any unusual value is communicated through status labels or guiding messages.

C.14 Saved Reports Screen

Purpose: Displays generated report files and allows them to be opened.

Appendix Table C.14: Functional documentation for Saved Reports Screen

Category	Details
Main user inputs / buttons	Open; Reload; Folder
Information displayed	Saved PDF list; File size; Time stamp
Backend response	Lists files saved in the report folder and opens selected reports.
Importance in BHMS workflow	This screen supports self-service operation, reduces user confusion and improves the completeness of the measurement workflow.

The Saved Reports screen features large buttons that are touch-sensitive, clearly labeled and have a black background which reflects the overall style of the app. The Saved Reports screen helps in ensuring that the BHMS is usable, in that it measures, reports, protects administrative functions or facilitates data storage.

- All the visible buttons on the screen are touch-enabled buttons.
- The screen has the same visual theme.
- The backend functionality is tied into the application code of the Raspberry Pi.
- Any unusual number is represented using a status label or guidance text.

Appendix D: Consolidated Data Interpretation And Workflow Notes

The notes above replace the lengthy repeated explanation notes. Detailed statistical analysis, tables, graphs, etc. are provided in Chapter 8 and Appendix A. The objective of this appendix is merely to furnish a quick reference guide to be used during viva, presentation, and future improvements to the Basic Health Monitoring Station.

Area	Important Interpretation for Report and Presentation
Height	Small differences can occur due to posture, head position and sensor alignment. The user should stand straight under the VL53L0X sensor.
Weight	Load-cell readings depend on platform stability and load distribution. Tare and calibration are required before use.
BMI	BMI is calculated from height and weight, therefore any height or weight deviation directly affects the BMI result.
Temperature	Non-contact temperature is affected by forehead distance, room condition and sensor alignment. The screw adjustment improves positioning.
Pulse and SpO2	Optical readings are sensitive to finger movement, pressure and contact. Voice prompts help guide correct finger placement.
Respiratory Rate	Respiratory rate is included for preliminary screening and should be validated with matched reference instruments in future work.
Time Efficiency	The BHMS reduces the need for several separate devices and supports faster screening with digital result storage.
Power Operation	The Raspberry Pi requires continuous power because a permanent battery backup was not implemented.
Feasibility	The prototype is suitable for preliminary screening, academic demonstration and institutional health awareness. It is not a certified diagnostic device.

Table D.1: Consolidated interpretation notes for BHMS results and workflow.

D.1 Report Generation and Voice-Assisted Workflow

After capturing all the readings in total, the Raspberry Pi will save the session and then classify the results and create a PDF report with the QR code generated. For mobile applications, the session values will be the same, but the PDF report gives a summary of the readings together with the diagnosis comments as well as what actions should be taken in case of an emergency. With voice assistance, the entire process is guided, including welcoming the participant, temperature readings, left finger for SpO2 and RR, right finger for the pulse, weighing and height.

Appendix F: Source Code Summary And Key Program Excerpts

This section will give a concise version of the source code for Basic Health Monitoring Station Application. The full source code of the application can be found separately from the project source code package/repo. Meanwhile, the source code excerpts included in this report are those considered crucial in understanding sensor capturing, weight calibration, communication between sensors and computer through serial interface, BMI calculations, QR report creation, and GUI processes.

F.1 Source Code File Organization

Module / File	Main Purpose
ESP32 integrated firmware	Reads sensors, controls measurement sequence and sends labelled JSON values to Raspberry Pi.
ESP32 HX711 ESP-NOW weight sender	Receives load-cell data through HX711, applies 10-point calibration and sends weight to the receiver ESP32.
main.py	Launches the Raspberry Pi PyQt application and initializes navigation, services and screen flow.
serial_service.py	Receives and parses ESP32 serial packets containing height, weight, BMI, temperature, SpO2, pulse and respiratory-rate data.
report_service.py	Creates the PDF health report using measured values, classifications, diagnosis summary and recommendations.
qr_service.py	Generates the QR code and session viewer link for mobile access.
diagnosis_service.py	Applies rule-based health classification and creates health remarks.

screens/*.py	Implements welcome, instruction, measurement, result, QR, diagnosis, consult, settings, storage and admin screens.
data/config/*.json	Stores threshold, calibration and application configuration values.

Table F.1: Compact source-code organization for the BHMS project.

F.2 ESP32 Load-Cell Calibration and ESP-NOW Sender Excerpt

```
#include <WiFi.h>
#include <esp_now.h>
#include "HX711.h"

HX711 myScale;
uint8_t dataPin = 32;
uint8_t clockPin = 33;
uint8_t receiverAddress[] = {0xF4, 0x2D, 0xC9, 0x72, 0x84, 0x60};

const float MAX_WEIGHT_KG = 200.0;
const int NUM_POINTS = 10;
const float calWeightKg[NUM_POINTS] = {20,40,60,80,100,120,140,160,180,200};
const float calRawValue[NUM_POINTS] = {
    442526.28, 838240.94, 1244028.00, 1652952.37, 2053779.00,
    2470557.25, 2908480.25, 3327427.50, 3718852.50, 4045510.25
};

typedef struct struct_message {
    float weightKg;
    float rawValue;
    bool overload;
} struct_message;

float convertRawToKg(float raw) {
    float previousRaw = 0.0;
    float previousWeight = 0.0;
    for (int i = 0; i < NUM_POINTS; i++) {
        float nextRaw = calRawValue[i];
        float nextWeight = calWeightKg[i];
        if (raw >= previousRaw && raw <= nextRaw) {
            return previousWeight + ((raw - previousRaw) *
                (nextWeight - previousWeight) / (nextRaw - previousRaw));
        }
    }
}
```

```

    previousRaw = nextRaw;
    previousWeight = nextWeight;
}
return calWeightKg[NUM_POINTS - 1];
}

void loop() {
    float rawValue = myScale.get_value(10);
    float weightKg = convertRawToKg(rawValue);
    bool overload = weightKg > MAX_WEIGHT_KG;
    struct_message outgoingData = {weightKg, rawValue, overload};
    esp_now_send(receiverAddress, (uint8_t *)&outgoingData, sizeof(outgoingData));
    delay(500);
}

```

F.3 ESP32 BMI and JSON Data Transfer Excerpt

```

float calculateBMI(float weightKg, float heightCm) {
    float h = heightCm / 100.0;
    if (weightKg <= 0 || heightCm <= 0 || h <= 0) return 0.0;
    return weightKg / (h * h);
}

void sendFullResultPacket() {
    String json = "{";
    json += "\"type\": \"measurement\",";
    json += "\"weight\": " + String(finalWeight, 1) + ",";
    json += "\"height\": " + String(finalHeight, 1) + ",";
    json += "\"bmi\": " + String(finalBMI, 1) + ",";
    json += "\"temperature\": " + String(finalTemperature, 1) + ",";
    json += "\"spo2\": " + String(finalSpO2) + ",";
    json += "\"pulse_rate\": " + String(finalPulse) + ",";
    json += "\"respiratory_rate\": " + String(finalRespiration);
    json += "}";
    PiSerial.println(json);
    PiSerial.flush();
}

```

F.4 Raspberry Pi Serial Parsing Excerpt

```

def process_serial_packet(packet: str) -> dict:
    """Convert ESP32 serial packet into a measurement dictionary."""
    packet = packet.strip()
    if not packet:
        return {}

```

```
data = json.loads(packet) if packet.startswith("{}") else parse_key_value_packet(packet)
return {
    "height": data.get("height"),
    "weight": data.get("weight"),
    "bmi": data.get("bmi"),
    "temperature": data.get("temperature"),
    "spo2": data.get("spo2") or data.get("spo2_percent"),
    "pulse_rate": data.get("pulse_rate") or data.get("pulse") or data.get("bpm"),
    "respiratory_rate": data.get("respiratory_rate") or data.get("rr"),
}
```

F.5 QR Report Generation Excerpt

```
def generate_session_report(session_data):
    """Generate PDF report and QR code after measurement completion."""
    report_path = report_service.create_pdf_report(session_data)
    qr_path = qr_service.create_qr_for_report(report_path)
    database_service.save_session(session_data, report_path, qr_path)
    return report_path, qr_path

# The generated PDF report contains measured values, health categories,
# diagnosis summary, recommendations and emergency guidance.
```

F.6 Voice-Assisted Workflow Trigger Excerpt

```
def play_voice_prompt(stage_name: str):
    """Play the audio prompt linked to the current measurement stage."""
    prompt_map = {
        "welcome": "001_welcome.mp3",
        "temperature": "002_adjust_temperature_screw.mp3",
        "spo2": "003_place_left_index_finger.mp3",
        "pulse": "004_place_right_index_finger.mp3",
        "weight": "005_stand_still_for_weight.mp3",
        "height": "006_stand_straight_for_height.mp3",
        "result": "007_view_result_and_scan_qr.mp3",
    }
    audio_file = prompt_map.get(stage_name)
    if audio_file:
        dfplayer_service.play(audio_file)
```

F.7 Appendix Note

Only the selected excerpts have been included in the hard copy of the report in order to keep it brief. The entire folder containing the source code must accompany the report when submitted via GitHub, Google Drive, or zipped source code.